

Laboratory Filtration Products



Filtration and ultrafiltration are essential process steps in nearly all environmental, chemistry and bioscientific laboratory applications.

Sartorius supplies a wide range of individual filter papers, microporous membranes, filtration devices, ultrafiltration units and protein purification devices to suit these applications. This catalog provides a condensed overview of the Sartorius Lab Filtration product range. Please contact us directly for specialty catalogs – available for in-depth technical information.

Table of Contents

Ultrafiltration

- 10 Protein Concentration
- 46 DNA Concentration
- 52 Protein Purification
- 62 Virus Purification and Concentration

Filtration Devices

- 64 Minisart® Syringe Filters
- 80 Claristep® Syringeless Filtration
- 86 Sartolab® P20 and Sartolab® P20 Plus
- 89 Sartolab® RF|BT
- 92 NEW! Sartolab® RF 50
- 94 NEW! Sartolab® Multistation
- 96 Sartoclear Dynamics® Lab P15
- 98 NEW! Sartoclear Dynamics® Lab V
- 102 Chemical Compatibility

Basic Filtration

- 108 Filter Papers
- 112 Glass and Quartz Microfiber Filters
- 116 Membrane Filters
- 124 Blotting | Chromatography Papers & Blotting Membranes
- 126 Filtration Equipment
- 136 Combisart® Manifolds









Ultrafiltration

Table of Contents

- 8 Membrane Selection Guide
- Protein Concentration
 - 10 Vivaspin[®] 500
 - 12 Vivaspin® 2
 - 15 Centrisart® I
 - 17 Vivaspin[®] Turbo 4
 - 19 Vivaspin® 6
 - 21 Vivaspin® 15R
- 23 Vivaspin® Turbo 15
- 25 Vivaspin® 20
- 28 Vivaclear Centrifugal Filters
- 29 Vivacell 100
- 33 Vivaflow® 50
- 36 Vivaflow® 50R
- 38 Vivaflow® 200
- 41 Vivapore® Solvent Absorption Concentrators
- 43 Ultrafiltration Membrane Filters
- DNA Concentration
 - 47 Vivacon® 500
 - 50 Vivacon® 2
- Protein Purification
 - 53 Vivapure® Ion Exchange Protein Purification Products
- Virus Purification and Concentration
 - 55 Vivapure® Virus Purification and Concentration Kits
 - 56 Adenovirus Purification with AdenoPACK Kits
 - 58 Lentivirus Purification with LentiSELECT Kit



Ultrafiltration is a convective process using anisotropic semi-permeable membranes to separate macromolecular species and solvents – primarily on the basis of size. By allowing solvents and salts to pass the ultrafiltration membrane while macromolecules are hindered to pass it, they are being concentrated. However, ultrafiltration can also be applied in solvent exchange applications. Multiple concentration and refilling steps will progressively lead to a buffer exchange, replacing lengthy techniques like dialysis. Although ultrafiltration is not a standard method for separating and fractionating macromolecules, it can be used as such if the macromolecules differ at least 10 times in size. Ultrafiltration is a gentle, non-denaturing method that is more efficient and flexible than other processes.

Ultrafiltration Methods

Sartorius offers you a comprehensive range of ultrafiltration process methods for the concentration of your biological samples.

- Centrifugal Concentration (100 μl to 100 mL Starting Volume)
- Pressure Ultrafiltration(5 mL to 100 mL Starting Volume)
- Crossflow (Tangential Flow)(100 mL to 5l Starting Volumes)
- Static Absorption(3 mL to 20 mL Starting Volume)

Further information about the operational methods can be found on page 8.

Typical Applications for Ultrafiltration

- Concentration | desalting of proteins, enzymes, DNA, monoclonal antibodies, immunoglobulins, viruses and nanoparticles
- Bence Jones Protein concentration from urine samples prior to capillary electrophoresis
- Forensic DNA sample concentration prior to sequencing reaction
- Peptide fractionation in FASP (filter-aided sample preparation)
- Free drug | hormone assays
- Removal of primers from PCR amplified DNA
- Removal of labeled amino acids and nucleotides
- Deproteinization of samples
- General purpose laboratory concentration and desalting of proteins, enzymes, DNA, biomolecules, viruses, antibodies and immunoglobulins

Membrane Performance Characteristics

Sartorius offers an extended range of membranes to cover the great majority of ultrafiltration requirements.

- Polyethersulfone (PES)
- Cellulose Triacetate (CTA)
- Hydrosart[®]

Further information about the properties of the different membrane types can be found on page 9.

Process Optimization

When the highest recoveries are crucial, particularly with solute quantities in the microgram range, Sartorius recommends considering the following tips for optimal ultrafiltration results:

- Select the lowest MWCO membrane that suits your application. For the highest recovery, choose a membrane MWCO which is at least half of the molecular weight of the solute to be retained.
- Avoid over-concentration. The smaller the final concentrate volume, the more difficult it is to achieve complete recovery. If feasible, rinse the device with one or more drops of buffer after the first concentration cycle and then recover it again.
- Pretreat the device overnight in distilled water with a passivation solution such as 5% SDS, Tween 20 or Triton X. Rinse thoroughly before use.



Membrane Selection Guide

Ultrafiltration Methods

Sartorius offers you a comprehensive range of ultrafiltration process methods for the concentration of your biological samples. The guide below will help you select the most suitable device according to sample volume and available equipment, as well as your desired filtration speed and process control.



Centrifugal Concentration (100 µl to 100 mL Starting Volume)

Driven by the centrifugal force, ultrafiltration can be used to purify and concentrate proteins, macromolecules and nanoparticles. This gentle process is quick to set up and offers fast filtration speeds for most solutions. Sartorius offers seven

Vivaspin® devices and the Vivacell® 100 and Centrisart® I for protein concentration, as well as Vivacon® devices for DNA and peptide concentration | fractionation.



Pressure Ultrafiltration (5 mL to 100 mL Starting Volume)

Pressurized air or inert gas provide the filtration vector for pressure filtration. To speed up the filtration process, the pressurized Vivacell® 100 or the Vivaspin® 20 can be placed on an orbital laboratory shaker. Agitation helps to prevent membrane blocking and ensures high filtration speed. Vivaspin® 20, and Vivacell® 100 can be run with gas pressure. The Vivaspin® 20 can also be run by


pressure-fugation, a unique Sartorius method combining gas and pressure with centrifugation.

This is the fastest of all methods, providing process times that are typically 30 to 50 percent faster than centrifugation.



Crossflow (Tangential Flow) (100 mL to 5 L Starting Volumes)

The sample is pumped across an ultrafiltration membrane and then returned to the original reservoir by building up pressure at the outlet of the device. The solution is progressively concentrated as solvent and micromolecules pass through the membrane into a separate filtrate vessel. Reusable Vivaflow® 50R and Vivaflow® 200, as well as disposable Vivaflow® 50, are offered for your dedicated laboratory crossflow filtration.

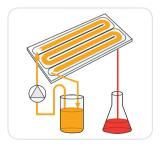

Static Absorption (3 mL to 20 mL Starting Volume)

This technique uses an absorbent cellulose pad mounted behind the ultrafiltration membrane to draw solvents and micromolecules through the membrane. The retained macromolecules thus concentrate at the bottom of the sample container.

No additional equipment is needed. These devices are ideal for clinical applications like urine concentration prior to further analysis. Both Vivapore® 5 and Vivapore® 20 offer this procedure.



Membrane Performance Characteristics

Sartorius offers an extended range of membranes to cover the great majority of ultrafiltration requirements. To select the most appropriate membranes for your application please refer to the following guide. Please note however that membrane behavior and ultimate performance, largely depends on the specific characteristics of the solution being processed. Where available alternatives membranes should be tested to optimize the process performance.

Polyethersulfone (PES)

Concentration | desalting of column eluates, cell culture supernatants, etc. This is a general purpose membrane that provides excellent performance with most solutions when retentate recovery is of primary importance. Polyethersulfone membranes are usually preferred for their low fouling characteristics, exceptional flux and broad pH range.

Cellulose Triacetate (CTA)

Free | bound drug studies and whenever the filtrate is being analyzed High hydrophilicity and very low non-specific binding characterize this membrane. Cast without any membrane support that could trap or bind passing micro solutes, these membranes are preferred for sample cleaning and protein removal and when high recovery of the filtrate solution is of primary importance.

Hydrosart[®]

Concentration | desalting of column eluates, Hydrosart® membrane evaluation for upscaling Hydrosart® demonstrates the same properties as regenerated cellulose, but with the added benefit of enhanced performance characteristics and extremely low protein binding, making it the membrane of choice for applications such as concentration and desalting of immunoglobulin fractions.

Membrane Selection Guide

The molecular weight cut-off (MWCO) is the molecular weight of molecules (e.g. globular proteins) which are retained by the membrane to an extent of 90%. Therefore, to ensure the highest recovery, select a membrane with a MWCO that represents half the molecular weight of the solute to be retained.

Sartorius Vivaspin® ultrafiltration units are designed to concentrate protein solutions. Therefore, the membranes in the devices are tested for the retention of proteins and not the passage of proteins into the filtrate vessel. PES and Hydrosart® membranes have support structures, which might lead to some loss of protein after the passage through the membrane.

Recommended MWCO (Da):

Application	<5,000	10,000	30,000	50,000	100,000	>300.000
Bacteria					•	•
DNA fragments		•	•	•	•	
Enzymes	•	•				
Growth factors	•	•				
mAB			•	•	•	
Nucleic acids	•	•	•	•	•	
Oligonucleotides	•					
Peptides	•					
Virus			•	•	•	
Yeast					•	•

Vivaspin® 500



mL centrifuge tubes.

100 μ L to 500 μ L Samples

Vivaspin® 500 μL centrifugal filter units

offer a simple, one-step procedure for


sample preparation. They can effectively

be used in fixed-angle rotors accepting 2.2

Specifications



Please use a fixed angle rotor for 2 mL reaction vials.

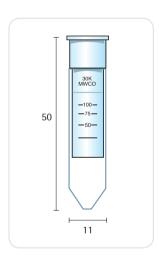

Vivaspin [®] 500		
Concentrator capacity	Swing-bucket rotor Fixed-angle rotor	do not use 500 µL
Dimensions	Total length Width Active membrane area Hold-up volume, membrane and support Dead-stop volume	50 mm 11 mm 0.5 cm ² < 5 μL 5 μL
Materials of construction	Body Filtrate vessel Concentrator cap Membrane	Polycarbonate Polypropylene Polycarbonate Polyethersulfone

The patented vertical membrane design

concentrations – even with particle-laden

and thin channel filtration chamber

(US 5,647,990) minimizes membrane

fouling and provides high-speed

solutions.

Typical Performance Characteristics

	Time to concentrate up to 30× [min.] at 20°C and solute recovery % Fixed angle		
Rotor			
Centrifugal force	12,000 g		
Start volume	500 μΙ		
	Min.	Rec.	
Aprotinin 0.25 mg/mL (6,500 MW)			
3,000 MWCO PES	30	96%	
BSA 1.0 mg/mL (66,000 MW)			
5,000 MWCO PES	15	96%	
10,000 MWCO PES	5	96%	
30,000 MWCO PES	5	96%	
IgG 0.25 mg/mL (160,000 MW)			
30,000 MWCO PES	10	96%	
50,000 MWCO PES	10	96%	
100,000 MWCO PES	10	96%	

11

Ordering Information

Vivaspin® 500 Polyethersulfone	Qty./Pkg.	Prod. No.
3,000 MWCO	25	VS0191
3,000 MWC0	100	VS0192
5,000 MWCO	25	VS0111
5,000 MWCO	100	VS0112
10,000 MWCO	25	VS0101
10,000 MWCO	100	VS0102
30,000 MWCO	25	VS0121
30,000 MWCO	100	VS0122
50,000 MWCO	25	VS0131
50,000 MWCO	100	VS0132
100,000 MWCO	25	VS0141
100,000 MWCO	100	VS0142
300,000 MWCO	25	VS0151
300,000 MWCO	100	VS0152
1,000,000 MWCO	25	VS0161
1,000,000 MWCO	100	VS0162
0.2 μm	25	VS0171
0.2 μm	100	VS0172
Starter pack	25	VS01S1
(5 of each 5 kDa, 10 kDa, 30 kDa, 50 kDa, 100	kDa1	

(5 of each 5 kDa, 10 kDa, 30 kDa, 50 kDa, 100 kDa)


Visit us at www.sartorius.com/Vivaspin500 to get additional info. Find instructions on how to use Vivaspin® 500 for

- Desalting and buffer exchange
- Preparation of biological nanoparticles and medical nanocarriers
- Concentration and purification of viruses
- Urine protein concentration
- Concentration to a predefined volume
- Concentration of diluted samples with increased recovery

Vivaspin® 2

Choice of Membranes



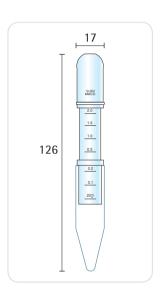

0.4 mL to 2 mL Samples

The Vivaspin® 2 bridges the gap between the 500 μ L and 4 mL centrifugal concentrators. This device combines the speed of the classic Vivaspin® products with low internal surface and membrane area for superior recoveries from very dilute solutions.

Available with a choice of PES, Cellulose Triacetate and Hydrosart® membranes, Vivaspin® 2 offers the highest flexibility for process optimization.

Also unique to the Vivaspin® 2 is the choice of directly pipetting the concentrate from the dead-stop pocket built into the bottom of the concentrator or alternatively reverse spinning into the concentrate recovery cap. Both methods result in nearly total concentrate recoveries.

Specifications

Vivaspin® 2

Concentrator capacity	Swing-bucket rotor	3 mL
	Fixed-angle rotor	2 mL
Dimensions	Total length	126 mm
	Width	17 mm
	Active membrane area	1.2 cm ²
	Hold-up volume of membrane	< 10 μL
	Dead-stop volume	8 μL
Materials of construction	Body	Polycarbonate
	Filtrate vessel	Polycarbonate
	Concentrator cap	Polycarbonate
	Membrane	PES, CTA, HY

Performance Characteristics

	Time to concentrate up to 30× [min.] at 20°C and solute recovery % Fixed angle		
Rotor			
Centrifugal force	5,000 g		
Start volume	2 mL		
	Min.	Rec.	
Aprotinin 0.25 mg/mL (6,500 MW)			
3,000 MWCO PES	50	96%	
BSA 1.0 mg/mL (66,000 MW)			
5,000 MWCO PES	12	98%	
5,000 MWCO CTA	50		
5,000 MWCO Hydrosart®	22	98%	
10,000 MWCO PES	8	98%	
10,000 MWCO CTA	10	96%	
10,000 MWCO Hydrosart®	12	98%	
20,000 MWCO CTA	5	97%	
30,000 MWCO PES	8	97%	
30,000 MWCO Hydrosart®	5	95%	

13

Typical Performance Characteristics

	Time to concentrate up to 30× [min.] at 20°C and solute recovery %		
Rotor	Fixed angle	Fixed angle	
Centrifugal force	5,000 g		
Start volume	2 mL		
	Min.	Rec.	
IgG 0.25 mg/mL (160,000 MW)			
20,000 MWCO CTA	6	97%	
30,000 MWCO PES	10	96%	
50,000 MWCO PES	10	96%	
100,000 MWCO PES	8	95%	

Ordering Information

Vivaspin® 2 Polyethersulfone	Qty./Pkg.	Prod. No.
3,000 MWCO	25	VS0291
3,000 MWC0	100	VS0292
5,000 MWCO	25	VS0211
5,000 MWCO	100	VS0212
10,000 MWCO	25	VS0201
10,000 MWCO	100	VS0202
30,000 MWCO	25	VS0221
30,000 MWCO	100	VS0222
50,000 MWCO	25	VS0231
50,000 MWC0	100	VS0232
100,000 MWCO	25	VS0241
100,000 MWC0	100	VS0242
300,000 MWCO	25	VS0251
300,000 MWC0	100	VS0252
1,000,000 MWCO	25	VS0261
1,000,000 MWCO	100	VS0262
0.2 μm	25	VS0271
0.2 μm	100	VS0272
Starter pack (5 of each 5 kDa, 10 kDa, 30 kDa, 50 kDa, 100 kDa)	25	VS02S1



PES, CTA, or Hydrosart® membranes; Filtrate container fits standard 15 mL tube carriers



Direct pipette recovery or choice of reverse spinning concentrate into sample cap

Vivaspin® 2 Cellulose Triacetate	Qty./Pkg.	Prod. No.
5,000 MWCO	25	VS02U1
5,000 MWCO	100	VS02U2
10,000 MWCO	25	VS02V1
10,000 MWCO	100	VS02V2
20,000 MWCO	25	VS02X1
20,000 MWC0	100	VS02X2

Vivaspin [®] 2 Hydrosart [®]	Qty./Pkg.	Prod. No.
2,000 MWCO	25	VS02H91
2,000 MWCO	100	VS02H92
5,000 MWCO	25	VS02H11
5,000 MWCO	100	VS02H12
10,000 MWCO	25	VS02H01
10,000 MWCO	100	VS02H02
30,000 MWCO	25	VS02H21
30,000 MWC0	100	VS02H22

Ordering Tips


- Choose a membrane pore size at least 50% smaller than the size of the molecule to be retained.
- It is usually best to select polyethersulfone membranes to achieve the fastest concentrations.
- Usually choose cellulose triacetate for protein removal or ultrafiltrate recovery.
- Usually choose Hydrosart® membranes for the highest recovery with Ig fractions.

Visit us at www.sartorius.com/Vivaspin2 to get additional info. Find instructions on how to use Vivaspin® 2 for

- Desalting and buffer exchange
- Preparation of biological nanoparticles and medical nanocarriers
- Concentration and purification of viruses
- Urine protein concentration
- Concentration of diluted samples with increased recovery
- Sample preparation for radio immunoassay

Centrisart® I



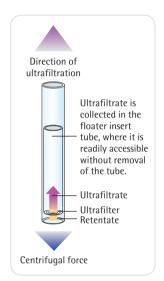

Centrisart[®] I* is a ready-to-use unit for small-volume, centrifugal ultrafiltration to separate proteins from low molecular weight substances in biological samples.

Centrisart® I features a unique design that enables ultrafiltration in the direction opposite to centrifugal force. This is so effective in preventing premature blockage of the filter that even whole blood samples can be deproteinized.

The ultrafiltrate is collected in the floater insert tube, where it is readily accessible without removing the tube.

Centrisart® I is ideal for the following applications:

- Drug binding studies
- Determination of metabolites in serum
- Protein removal from blood samples
- Cleaning of liposomes
- Virus removal

Specifications

Centrisart® I

Centribure 1		
Concentrator capacity	Swing-bucket rotor	2.5 mL
	Fixed-angle rotor	2.5 mL
Dimensions	Total length	93 mm
	Width	14 mm
	Active membrane area	0.79 cm^2
	Hold-up volume of membrane	< 5 μL
	Dead-stop volume	100 μL
Materials of construction	Centrifuge tube	Polystyrene
	Floater tube	Cellulose propionate
	Cap	Polyethylene
	Membrane	CTA, PES

Typical Performance Characteristics

	Time to filter 50% of sample volume	Time to filter 90% of sample volume	Passage of sample species volume
BSA 1.0 mg/mL (66,000 MW)			
5,000 MWCO	300 min	N A	0%
10,000 MWCO	35 min	80 min	2%
20,000 MWC0	9 min	20 min	2%
IgG 0.25 mg/mL (160,000 MW)			
100,000 MWCO	13 min	35 min	3%
Blue Dextran 0.1 mg/mL (2,000,000 MW)			
300,000 MWCO	9 min	25 min	28%

2.5 mL samples were loaded into each device. The devices were centrifuged at 2,000 g until the required filtrate volumes had been reached.

^{*} Centrisart is a registered trademark in the U.S. and the European Union

Ordering Information

	Qty./Pkg.	Prod. No.
5,000 MWCO CTA	12	13229E
10,000 MWCO CTA	12	13239E
20,000 MWCO CTA	12	13249E
100,000 MWCO PES	12	13269E
300,000 MWCO PES	12	13279E
Starter pack	12	13209E
(3 units each of 5 kDa, 10 kDa, 20 kDa, 100 kDa)		

References

P. Nebinger and M. Koel: Determination of acyclovir by ultrafiltration and high-performance liquid chromatography.

J. Chromatography 619, 342-344 (1993)

F. da Fonseca-Wollheim, K.-G. Heinze, K. Lomsky and H. Schreiner: Serum ultrafiltration for the elimination of endogenous interfering substances in creatinine determination. J. Clin. Chem. Clin. Biochem. 26, 523–525 (1988) R. H. Christenson, S. D. Studenberg, S. Beck-Davis and F. A. Sedor: Digoxin-like immunoreactivity eliminated from serum by centrifugal ultrafiltration before fluorescence polarization immunoassay of digoxin. Clinical Chemistry 33, 606-608 (1987)

Vivaspin® Turbo 4



2 mL to 4 mL Samples


Vivaspin® Turbo 4 is the newest member of the ultrafiltration family and allows the fastest sample concentration with the highest recoveries.

This device can handle up to 4 mL sample volumes in swing-bucket rotors and in fixed-angle rotors that accept 15 mL centrifuge tubes.

The Vivaspin® Turbo 4 optimized design and sleek internal profile ensure maximum process speeds all the way down to the last few microliters, resulting in more than 100-fold concentration.

17

The UV joining technology ensures smooth joint transition between the membrane and the plastic housing – allowing removal of the entire sample concentrated in the unique, pipette-friendly dead-stop pocket.

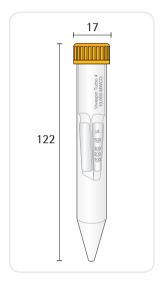

Specifications

Vivaspin® Turbo 4

Concentrator capacity	Swing-bucket rotor	4 mL
	Fixed-angle rotor	4 mL
Dimensions	Total length	122.5 mm
	Width	17 mm
	Active membrane area	3.2 cm ²
	Hold-up volume of membrane	< 10 µL
	Dead-stop volume	40 μL
	swing-bucket rotor	
	Dead-stop volume fixed-angle rotor (25°)	60 μL
Materials of construction	Body	Styrene butadiene copolymere
	Filtrate vessel	Polypropylene
	Concentrator cap	Polypropylene
	Membrane	Polyethersulfone
		•

Typical Performance Characteristics

	Time to concentrate up to 30× [min.] at 20°C and solute recovery %			
Rotor	Swing out	t	Fixed ang	gle (25°)
Centrifuge speed	4,000 xg		7,500 xg	
Starting vol.	4 mL		4 mL	
	Min.	Rec.	Min.	Rec.
Cytochrome c (12,400 MW) 3,000 MWCO PES 5,000 MWCO PES	60 40	98% 95%	80 50	96% 94%
Lysozyme (14,300 MW) 3,000 MWCO PES 5,000 MWCO PES	65 50	95% 94%	70 60	93% 92%
α-Chymotrypsin (25,000 MW) 10,000 MWCO PES	10	95%	8	95%
BSA (66,000 MW) 10,000 MWCO PES 30,000 MWCO PES	10 8	98% 96%	7 6	97% 97%

Performance Characteristics

	Time to concentrate up to 30× [min.] at 20°C and solute recovery %			
	Min.	Rec.	Min.	Rec.
IgG (160,000 MW)				
30,000 MWCO PES	18	94%	13	92%
50,000 MWCO PES	16	93%	12	90%
100,000 MWCO PES*	17	94%	13	92%

^{* 3,000} xg swing-out | 5,000 xg fixed angle

Ordering Information

Vivaspin® Turbo 4 Polyethersulfone	Qty./Pkg.	Prod. No.
3,000 MWCO	25	VS04T91
3,000 MWCO	100	VS04T92
5,000 MWCO	25	VS04T11
5,000 MWCO	100	VS04T12
10,000 MWCO	25	VS04T01
10,000 MWCO	100	VS04T02
30,000 MWCO	25	VS04T21
30,000 MWCO	100	VS04T22
50,000 MWCO	25	VS04T31
50,000 MWCO	100	VS04T32
100,000 MWCO	25	VS04T41
100,000 MWCO	100	VS04T42

Visit us at www.sartorius.com/VivaspinTurbo4 to get additional info. Find instructions on how to use $Vivaspin^{\circ}$ Turbo 4 for

- Desalting and buffer exchange
- Preparation of biological nanoparticles and medical nanocarriers
- Concentration and purification of viruses
- Urine protein concentration
- Separation of proteins and metabolites for disease detection

19

Vivaspin® 6




2 mL to 6 mL Samples

Vivaspin® 6 mL concentrators have been developed to offer increased volume flexibility and performance.

Vivaspin® 6 can process an impressive 6 mL in either swing-bucket or fixed-angle rotors accepting standard 15 mL conical bottom test tubes.

The Vivaspin® 6 features twin vertical membranes for unparalleled filtration speeds and 100× plus concentrations. The remaining volume is easy to read off the printed scale on the side of the concentrator and the modified dead-stop pocket further simplifies direct pipette recovery of the final concentrate.

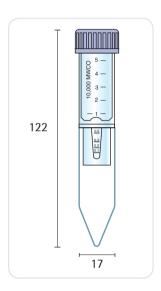

Specifications

Vivaspin® 6

Concentrator capacity	Swing-bucket rotor	6 mL
	Fixed-angle rotor	6 mL
Dimensions	Total length	122 mm
	Width	17 mm
	Active membrane area	2.5 cm ²
	Hold-up volume of membrane	< 10 μL
	Dead-stop volume	30 μL
Materials of construction	Body	Polycarbonate
	Filtrate vessel	Polycarbonate
	Concentrator cap	Polypropylene
	Membrane	Polyethersulfone

Typical Performance Characteristics

		centrate up to lute recovery		:
Rotor	Swing bucket		25° Fixed ang	le
Centrifugal force	3,000 g		7,500 g	
Start volume	6 mL		6 mL	
	Min.	Rec.	Min.	Rec.
Cytochrome c 0.25 mg/mL (12,400 MW)				
5,000 MWCO PES	_	_	90	97%
BSA 1.0 mg/mL (66,000 MW) 5,000 MWCO PES	20	98%	12	98%
10,000 MWCO PES	13	98%	10	98%
30,000 MWCO PES	12	98%	9	97%
IgG 0.25 mg/mL (160,000 MW)				
30,000 MWCO PES	18	96%	15	95%
50,000 MWCO PES	17	96%	14	95%
100,000 MWCO PES	15	91%	12	91%
Latex beads 0.004% in DMEM + 10% FCS (0.055 μm)				
300,000 MWCO PES	_	-	25	99%
Latex beads 0.004% in DMEM + 10% FCS (0.24 μm)				
1,000,000 MWCO PES	_	_	4	99%
Yeast 1.0 mg/mL (S. Cerevisiae) 0.2 µm PES	4	97%	3	97%

Ordering Information

Vivaspin® 6 Polyethersulfone	Qty./Pkg.	Prod. No.
3,000 MWCO	25	VS0691
3,000 MWCO	100	VS0692
5,000 MWCO	25	VS0611
5,000 MWCO	100	VS0612
10,000 MWCO	25	VS0601
10,000 MWCO	100	VS0602
30,000 MWCO	25	VS0621
30,000 MWCO	100	VS0622
50,000 MWCO	25	VS0631
50,000 MWCO	100	VS0632
100,000 MWCO	25	VS0641
100,000 MWCO	100	VS0642
300,000 MWCO	25	VS0651
300,000 MWCO	100	VS0652
1,000,000 MWCO	25	VS0661
1,000,000 MWCO	100	VS0662
0.2 μm	25	VS0671
0.2 μm	100	VS0672
Starter pack (5 of each 5 kDa, 10 kDa, 30 kDa, 50 kDa, 100 kDa)	25	VS06S1

Visit us at www.sartorius.com/Vivaspin6 to get additional info. Find instructions on how to use Vivaspin° 6 for

- Desalting and buffer exchange
- Preparation of biological nanoparticles and medical nanocarriers
- Concentration and purification of viruses
- Urine protein concentration
- Concentration of diluted samples with increased recovery

Vivaspin® 15R




4 to 15 mL Samples

Vivaspin® 15R is designed for the volume segment of 2 to 15 mL and features a modified regenerated cellulose membrane; Hydrosart®. This membrane is ideal where extremely high recovery with very low adsorption is needed. Examples of these applications include desalting and concentration of Ig fractions.

Advantages

- Ultimate recovery with low adsorption (95-98%)
- Exceptionally fast concentration time (30 × in 15 min.)
- Convenient application protocol with easy handling
- Easy scale-up to Vivaflow 200 with Hydrosart[®] membrane for volumes up to 5 liters
- Very low hold-up volume (< 20 μL)

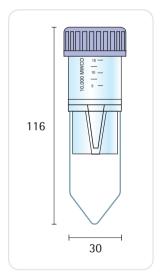

Specifications

Typical Performance Characteristics

Vivaspin® 15R

5,000 MWCO

(45,000 MW) 10,000 MWCO

30,000 MWCO

10,000 MWCO

Ovalbumin 1.0 mg/mL*

vivaspiii 15it		
Concentrator capacity	Swing-bucket rotor	15 mL
	Fixed-angle rotor	12.5 mL
Dimensions	Total length	116 mm
	Width	30 mm
	Active membrane area	3.9 cm ²
	Hold-up volume of membrane	< 20 μL
	Dead-stop volume	30 μL
Materials of construction	Body	Polycarbonate
	Filtrate vessel	Polypropylene
	Concentrator cap	Polypropylene
	Membrane	Hydrosart®

Time to concentrate up to $30 \times$ [min.] at

20°C and solute recovery %

98%

98%

98%

94%

45

18

14

12

98%

98%

98%

94%



Spin

Rotor Swing bucket 25° Fixed angle Centrifugal force 3,000 g 6,000 g Start volume 15 mL 12.5 mL Min. Min. Rec. Rec. Aprotinin 0.1 mg/mL* (6,500 MW) 5,000 MWCO 95% 45 95% 47 Cytochrome c 0.25 mg/mL* (12,400 MW) 5,000 MWCO 45 96% 45 96% 10,000 MWCO 25 94% 18 94% α-Chymotrypsin 0.25 mg/mL* (25,000 MW)

50

25

20

15



Performance Characteristics

	Time to concentrate up to 30× [min.] at 20°C and solute recovery %			.] at
	Min.	Rec.	Min.	Rec.
BSA 1.0 mg/mL* (66,000 MW) 30,000 MWCO	18	98%	15	98%
lgG 0.1 mg/mL* in DMEM (160,000 MW) 30,000 MWCO	30	98%	25	96%

^{*} proteins other than IgG made up in 50 mM potassium phosphate, 150 mM sodium chloride, pH 7.4

Ordering Information


Vivaspin® 15R Hydrosart®	Qty./Pkg.	Prod. No.
2,000 MWCO	12	VS15RH91
2,000 MWCO	48	VS15RH92
5,000 MWCO	12	VS15RH11
5,000 MWCO	48	VS15RH12
10,000 MWC0	12	VS15RH01
10,000 MWCO	48	VS15RH02
30,000 MWCO	12	VS15RH21
30,000 MWC0	48	VS15RH22

Visit us at www.sartorius.com/Vivaspin15R to get additional info. Find instructions on how to use $Vivaspin^{\circ}$ 15R for

- Desalting and buffer exchange
- Preparation of biological nanoparticles and medical nanocarriers
- Concentration and purification of viruses
- Urine protein concentration
- Concentration of diluted samples with increased recovery

Vivaspin® Turbo 15



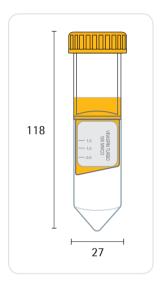

4 to 15 mL Samples

Vivaspin® Turbo 15 enables the fastest sample concentration with the highest recoveries. This device can handle a sample volume of up to 15 mL in swing-bucket rotors and of up to 11 mL in fixed-angle rotors that accept 50 mL centrifuge tubes.

The optimized design and sleek internal profile of the Vivaspin® Turbo 15 ensure maximum process speeds all the way down to the last few microliters, which results in more than 100-fold concentration.

The UV joining technology ensures smooth joint transition between the membrane and the plastic housing – allowing removal of the entire sample concentrated in the unique, pipette-friendly dead-stop pocket.

23

Specifications

Vivaspin® Turbo 15

vivaspiii Turoo 15		
Concentrator capacity	Swing-bucket rotor	15 mL
	Fixed-angle rotor (25°)	9 mL
Dimensions	Total length (concentrator insert)	77 mm
	Total length (in tube with cap)	118 mm
	Diameter (concentrator insert)	27 mm
	Active membrane area	7.2 cm ²
	Hold-up volume of membrane	<10 μΙ
	Dead-stop volume for swing-bucket rotor	100 μΙ
	Dead-stop volume for fixed-angle rotor	60 μΙ
Materials	Body	Styrene butadiene
		copolymer
	Filtrate vessel	Polypropylene
	Concentrator cap	Polypropylene
	Membrane	Polyethersulfone (PES)

Typical Performance Characteristics

		concentrate u I solute recove	_	.] at
Rotor	Swing bu	cket	Fixed ang	Jle (25°)
Centrifugal speed	4,000 g		4,000 g	
Start volume	15 mL		9 mL	
	Min.	Rec.	Min.	Rec.
Cytochrome c* (12,400 MW)				
3,000 MWCO PES	61	98%	86	97%
5,000 MWCO PES	30	98%	50	98%
Lysozyme* (14,300 MW)				
3,000 MWCO PES	56	98%	87	97%
5,000 MWCO PES	33	96%	50	96%
α-Chymotrypsin** (25,000 MW)			
10,000 MWCO PES	10	95%	10	95%
BSA** (66,000 MW)				
10,000 MWCO PES	10	99%	10	99%
30,000 MWCO PES	8	98%	10	98%

Typical Performance Characteristics

	Time to concentrate up to 20× [min.] at 20°C and solute recovery %			
	Min.	Rec.	Min.	Rec.
IgG** (160,000 MW)				
30,000 MWCO PES	23	95%	17	95%
50,000 MWCO PES	20	94%	15	94%

^{* 0.25} mg/mL

Ordering Information

Vivaspin® Turbo 15 Polyethersulfone	Qty./Pkg.	Prod. No.
3,000 MWCO	12	VS15T91
3,000 MWCO	48	VS15T92
5,000 MWCO	12	VS15T11
5,000 MWCO	48	VS15T12
10,000 MWCO	12	VS15T01
10,000 MWCO	48	VS15T02
30,000 MWCO	12	VS15T21
30,000 MWCO	48	VS15T22
50,000 MWCO	12	VS15T31
50,000 MWCO	48	VS15T32
100,000 MWCO	12	VS15T41
100,000 MWCO	48	VS15T42

Visit us at www.sartorius.com/VivaspinTurbo15 to get additional info. Find instructions on how to use Vivaspin® Turbo 15 for

- Desalting and buffer exchange
- Preparation of biological nanoparticles and medical nanocarriers
- Concentration and purification of viruses
- Urine protein concentration
- Concentration of diluted samples with increased recovery
- Concentration to a predefined volume
- Depyrogenation of the ultrafiltration devices
- Concentration of mammalian cell culture supernatants

^{** 1} mg/mL

Vivaspin® 20





Air pressure controller, VCA002

5 to 20 mL Samples

Vivaspin® 20 mL centrifugal concentrators have been developed to offer increased volume flexibility and performance.

Vivaspin® 20 handles up to 20 mL in swing-bucket centrifuges and 14 mL in 25° fixed-angle rotors that accept 50 mL centrifuge tubes.

Featuring twin vertical membranes for unparalleled filtration speeds, the Vivaspin® 20 can achieve 100x plus concentrations.


The remaining volume is easy to read off the printed scale on the side of the concentrator and the modified dead-stop pocket further simplifies direct pipette recovery of the final concentrate.

More Process Flexibility

Vivaspin® 20 is available with unique accessories and operating methods that are designed to provide more process flexibility and further time savings.

Gas Pressure Filtration

When an appropriate centrifuge is unavailable or for single sample processing, Vivaspin® 20 can be filled with up to 15 mL and then pressurized for bench-top concentration. For even faster processing, gas pressure can be combined with centrifugal force. "Pressure-fugation" is particularly suitable for difficult or viscous samples, such as serum, or for use of a low process temperature, which reduces filtration speed, and generally when minimum process time is essential.

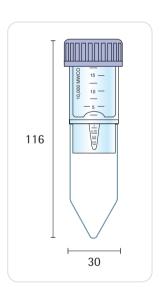

Specifications

Vivaspin® 20

Concentrator conceity	Curing buoket notes	20 mL
Concentrator capacity	Swing-bucket rotor	
	Fixed-angle rotor	14 mL
	With pressure head	15 mL
Dimensions	Total length	116 mm
	-	125 mm with pressure head
	Width	30 mm
	Active membrane area	6.0 cm ²
	Hold-up volume of membrane	< 20 μL
	Dead-stop volume	50 μL
Materials of construction	Body	Polycarbonate
	Filtrate vessel	Polycarbonate
	Concentrator cap	Polypropylene
	Pressure head	Acetal Aluminum
	Membrane	Polyethersulfone

Typical Performance Characteristics

	Time to concentrate up to 30× [min.] at 20°C and solute recovery %					
Mode	Centrifuge	Centrifuge	Bench top	Press-fuge		
Rotor	Swing bucket	25° Fixed angle	Pressure	Swing bucket		
Centrifugal speed pressure	3,000 g	6,000 g	4 bar	3,000 g + 4 bar		
Start volume	20 mL	14 mL	10 mL	10 mL		

Typical Performance Characteristics

Typical i ci ioi ilianee charae	ctciisti	CS						
		Time to concentrate up to 30× [min.] at 20°C and solute recovery %						
	Min.	Rec.	Min.	Rec.	Min.	Rec.	Min.	Rec.
Cytochrome c 0.25 mg/mL (12,400 MW) 3,000 MWCO PES	110	97%	180	96%	60	96%	_	_
BSA 1.0 mg/mL (66,000 MW)	110	37 70	100	3070		3070		
5,000 MWCO PES 10,000 MWCO PES 30,000 MWCO PES	23 16 13	99% 98% 98%	29 17 15	99% 98% 98%	50 32 32	98% 97% 97%	14 8 8	98% 97% 97%
IgG 0.25 mg/mL (160,000 MW)								
30,000 MWCO PES	27	97%	20	95%	46	94%	13	97%
50,000 MWCO PES	27	96%	22	95%	46	93%	13	96%
100,000 MWCO PES	25	91%	20	90%	42	88%	12	94%
Latex beads 0.004% in DMEN + 10% FCS (0.055 μm) 300,000 MWCO PES	20	99%	35	99%	10	99%	_	_
Latex beads 0.004% in DMEN + 10% FCS (0.24 μm) 1,000,000 MWCO PES	1 4	99%	12	99%	4	99%	_	_
Yeast 1.0 mg/mL (S. Cerevisiae) 0.2 µm PES	15	95%	5	95%	20	95%	2	95%


Ordering Information

Vivaspin® 20 Polyethersulfone	Qty./Pkg.	Prod. No.
3,000 MWCO	12	VS2091
3,000 MWCO	48	VS2092
5,000 MWCO	12	VS2011
5,000 MWCO	48	VS2012
10,000 MWCO	12	VS2001
10,000 MWCO	48	VS2002
30,000 MWCO	12	VS2021
30,000 MWCO	48	VS2022
50,000 MWCO	12	VS2031
50,000 MWCO	48	VS2032
100,000 MWCO	12	VS2041
100,000 MWCO	48	VS2042
300,000 MWCO	12	VS2051
300,000 MWCO	48	VS2052
1,000,000 MWCO	12	VS2061
1,000,000 MWCO	48	VS2062
0.2 μm	12	VS2071
0.2 μm	48	VS2072
Starter pack (2 each of 5 kDa, 10 kDa, 30 kDa, 50 kDa,	12	VS20S1

(2 each of 5 kDa, 10 kDa, 30 kDa, 50 kDa, 100 kDa, 0.2 μm)

27

Vivaspin® 20 Accessories	Qty./Pkg.	Prod. No.
Air pressure controller (APC)	1	VCA002
Charge valve for pressure head	1	VCA005
Diafiltration cups	12	VSA005
Female connector	1	VCA010
Male connector	1	VCA011
4 mm OD pneumatic tube (3 m)	1	VCA012
Vivaspin® 20 pressure head	1	VCA200

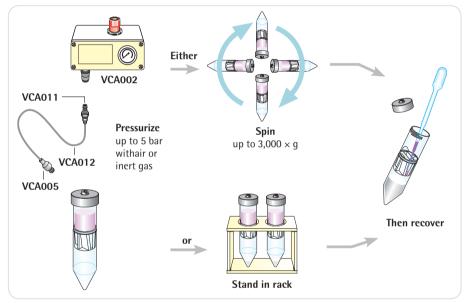

Using the Vivaspin® 20 pressure cap

Visit us at www.sartorius.com/Vivaspin20 to get additional info. Find instructions on how to use Vivaspin® 20 for

- Desalting and buffer exchange
- Preparation of biological nanoparticles and medical nanocarriers
- Concentration and purification of viruses
- Urine protein concentration
- Concentration of diluted samples with increased recovery
- The workflow in protein research laboratories

Vivaclear Centrifugal Filters



Vivaclear centrifugal filters are disposable microfiltration devices for the fast and reliable clarification|filtration of biological samples in the range 100 μL to 500 μL . They can be used in fixed-angle rotors accepting 2.2 mL centrifuge tubes.

Product Features

- High-flux polyethersulfone membrane
- 0.8 μm pore size
- Low hold-up volume (<5 μL)
- Fast and reproducible performance

Applications

- Clarification of samples before loading in Vivapure® protein purification spin columns
- Removal of particles and precipitates
- Filtration of plasma and serum
- Filtration of cells or cell debris

Specifications

Vivaclear Centrifugal Filters

Rotor	40–45° fixed-angle rotor, 500 μL	40–45° fixed-angle rotor, 500 μL				
Pore size	0.8 μm	0.8 μm				
Dimensions	Total length	43 mm				
	Filtrate collection tube diameter	11 mm				
	Active membrane area	0.34 cm^2				
	Hold-up volume,					
	membrane plus support	< 5 μL				
	Maximum RCF	$2,000 \times g$				
Materials of construction	Body	Polypropylene				
	Membrane	Polyethersulfone				
	Filtrate collection tube	Polypropylene				

Ordering Information

	Qty./Pkg.	Prod. No.
Vivaclear Mini 0.8 μm PES	100	VK01P042

Vivacell 100





Air pressure controller, VCA002

Vivacell 100 is a unique and innovative concentrator for volumes from 20 mL to 100 mL, and utilizes pressure, centrifuge or pressure-shake to rapidly concentrate even samples with very high-particle loading.

Vivacell 100 is designed for centrifugal concentration of samples up to 100 mL, which makes it the largest centrifugal unit available. At the same time, its new design allows for maximum centrifugal force of $2,000 \times g$ to be used for even faster concentration.

The patented vertical membrane design ensures the highest performance and unmatched flexibility.

Vivacell 100 Utilizes:


- Pressure
- Centrifuge
- Pressure-shake

Vivacell 100, when used as a centrifugal device, fits only into rotors that accept 250 mL bottles.

Vivacell 100 units can also be used for single or extremely sensitive samples in the pressurized mode only and left on a bench or placed on a swing-bucket laboratory shaker for faster concentration. It can also be kept in a pressurized mode in the refrigerator. Handling is facilitated by the use of quick connectors. In whichever mode Vivacell 100 is used, the vertical membrane design inhibits membrane fouling, while the built-in dead stop impedes concentration to dryness and loss of sample.



Patented vertical membrane design

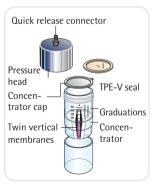

Filtrate container fits standard 250 mL rotors

Specifications

Vivacell 100

Concentrator capacity	Swing-bucket rotor With pressure head	90 mL 98 mL
Dimensions	Total length	123 mm centrifugal 197 mm with pressure head
	Width	62 mm
	Active membrane area	23.5 cm ²
	Hold-up volume of membrane	< 250 μL
	Dead-stop volume	350 μL
Operating requirements	Rotor type	Swing-bucket
	Rotor cavity	To fit 250 mL (62 mm)
		centrifuge bottles
		(maximum cavity depth
		105 mm)
	Maximum speed	2,000 g
	Maximum pressure	5 bar (75 psi)
Materials of construction	Body	Polycarbonate
	Filtrate vessel	Polycarbonate
	Concentrator cap	Polypropylene
	Pressure head seal	TPE-V
	Pressure head	Acetal
	Membrane	Polyethersulfone

Typical Performance Characteristics

-, prear : errormance enarace					
	Time to concentrate up to 30× [min.] at 20°C				
90 mL Start volume	In centrifuge 2,000 g swing-out rotor		As pressure cell 4 bar (60 psi) pressure		
		No agitation	Orbital shake	%	
BSA 1.0 mg/mL (66,000 MW)					
5,000 MWCO PES	22	75	25	96%	
10,000 MWCO PES	16	60	20	96%	
30,000 MWCO PES	16	60	20	94%	
IgG 0.25 mg/mL (160,000 MW)					
50,000 MWCO PES	20	70	30	94%	
100,000 MWCO PES	20	85	30	90%	
Latex beads 0.004% in DMEM + 10% FCS (0.055 µm)					
300,000 MWCO PES	35	_	120	99%	
Latex beads 0.004% in DMEM + 10% FCS (0.24 µm) 1.000.000 MWCO* PES	4	5	4	99%	
1,000,000 1111100 1 15	•	· ·	•	00 10	

^{* 2,000} g in centrifuge, 2 bar (29 psi) pressure

Ordering Information

Vivacell 100 Polyethersulfone with Polypropylene Concentrator Cap	Oty./Pkg.	Prod. No.
5,000 MWCO	2	VC1011
5,000 MWCO	10	VC1012
10,000 MWCO	2	VC1001
10,000 MWCO	10	VC1002
30,000 MWCO	2	VC1021
30,000 MWCO	10	VC1022
50,000 MWCO	2	VC1031
50,000 MWCO	10	VC1032
100,000 MWCO	2	VC1041
100,000 MWCO	10	VC1042
300,000 MWCO	2	VC1051
300,000 MWCO	10	VC1052
1,000,000 MWCO	2	VC1061
1,000,000 MWCO	10	VC1062
0.2 μm	2	VC1071
0.2 μm	10	VC1072

Vivacell 100 Accessories	Oty./Pkg.	Prod. No.
Air pressure controller (APC), complete with pressure gauge, regulator, over-pressure safety valve, female connector, 1 m extension line (4 mm pressure tubing) with male and female connectors and 1 m of 6 mm inlet tubing	1	VCA002
Female connector	1	VCA010
Male connector	1	VCA011
4 mm pressure tubing (3 m)	1	VCA012
Replacement seal made from TPE-V	10	VCA014
Vivacell 100 pressure head with replacement seals (5)	1	VCA800



Centrifuge use with polypropylene concentrator cap in swing out rotor

- Process convenience
- Low shear, non-foaming
- Less visual control



Pressure

use with pressure head VCA800

- Simplicity and the highest process control


 - Ideal for refrigerated use
- Slower concentrations



Pressure-Shake use with pressure

head VCA800

- Speed and process control
- Ideal for single samples

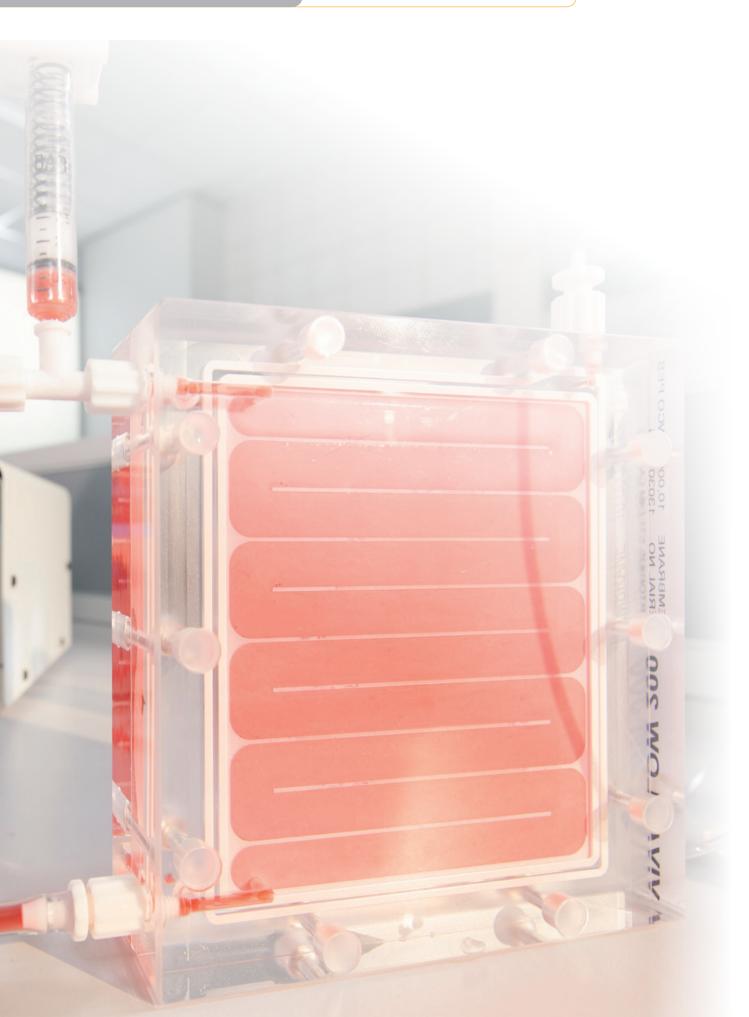


Vivaflow[®] 50






Flip-flow recirculation path

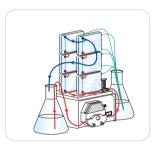

Multiple modules

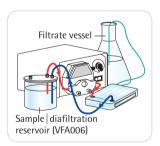

Single module

100 mL to 3 L Samples

The unique, patented Vivaflow** 50 system provides ease of use, performance, flexibility and economy that are unrivaled by any laboratory or pilot-scale filtration system on the market.

Unique Features

- The thin-channel, flip-flow recirculation path provides high crossflow velocities with minimum pump requirements
- No need for pressure holders
- Crystal clear for simple checking of remaining hold-up volume and membrane status
- Unique interlocking modules with series connectors for easy scale up
- Disposable

Unique Performance

- A single 50 cm² module will typically reduce 500 mL to less than 15 mL in under 50 minutes
- Less than 10 mL minimum system recirculation for the highest concentrations
- Less than 500 μL non-recoverable hold-up volume
- Nearly total recoveries achievable with a single 10 mL rinse

Unique "flip-flow" thin channel flow path results in high turbulence and linear velocity for exceptional flux even at high concentrations

Specifications

Vivaflow[®] 50

Vivaflow [®] 50		
Dimensions	Overall L×H×W	107 mm×84 mm×25 mm
	Channel W×H	15 mm × 0.3 mm
	Active membrane area	50 cm ²
	Hold-up volume (module)	1.5 mL
	Minimum recirculation volume	< 10 mL
	Non-recoverable hold-up	< 0.5 mL
Operating conditions	Pump flow	200 mL/min to 400 mL/min
	Maximum pressure	3 bar (45 psi)
	Maximum temperature	60°C
Materials of construction	Main housing	Polycarbonate
	Flow channel	TPX (PMP)
	Membrane support	TPX (PMP)
	Seals and O-rings	Silicone
	Pressure indicator	Polypropylene, SS spring
	Flow restrictor	Polypropylene
	Fittings	Nylon
	Tubing	PVC (medical grade)

Performance Characteristics

	Time to conce	entrate up to 2 essure, 20°C	20× [min.] at	
	Single device	Three devices	Solute recover	y %
	250 mL start volume	1 L start volume	Direct	10 mL rinse
BSA 1.0 mg/mL (66,000 MW)				
5,000 MWCO PES	34	49	96%	> 99%
10,000 MWCO PES	22	32	94%	> 99%
10,000 MWCO RC	38	55	96%	> 99%
30,000 MWCO PES	22	32	92%	99%
50,000 MWCO PES	20	29	92%	98%
γ Globulins 1.0 mg/mL (160,000 MW)				
100,000 MWCO PES	43	62	92%	98%
100,000 MWCO RC	40	58	92%	98%
Yeast 1.0 mg/mL (S. Cerevisiae)				
0.2 μm PES	33	47	92%	98%

Ordering Information

Vivaflow [®] 50*	Qty./Pkg.	Prod. No.
3,000 MWCO PES	2	VF05P9
5,000 MWCO PES	2	VF05P1
10,000 MWCO PES	2	VF05P0
30,000 MWCO PES	2	VF05P2
50,000 MWCO PES	2	VF05P3
100,000 MWCO PES	2	VF05P4
1,000,000 MWCO PES	2	VF05P6
0.2 μm PES	2	VF05P7
100,000 MWCO RC	2	VF05C4

^{*)} Vivaflow® 50 modules include filtrate tube, size 16 peristaltic tubing, flow restrictors and fittings.

Vivaflow[®] 50 Complete Set of Accessories

Pump (230 V), Easy Load pump head (size 16), tubing, 500 mL sample diafiltration reservoir, module stand, pressure indicator, T-connectors, series interconnectors	1	VFS502
Pump (115 V), Easy Load pump head (size 16), tubing,	1	VFS504
500 mL sample diafiltration reservoir, module stand,		
pressure indicator, T-connectors, series interconnectors		


VFA004
VFA009
VFA030
VFA031
VFA032
VFA034
VFA035
VFP001
VFP002
VFA006
VFA012
VFA016
VFA020

Visit us at www.sartorius.com/Vivaflow50 to get additional info. Here you can find instructions on how to use Vivaflow® 50 for


- Measurement of soluble trace metals in seawater
- Preparation of biological nanoparticles and medical nanocarriers
- Concentration and purification of viruses

Vivaflow[®] 50R



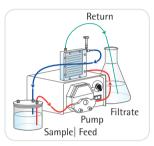

Vivaflow® 50R - Single module

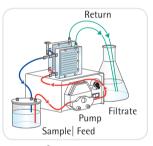

Vivaflow[®] 50R – Two modules

100 mL to 1 L Samples

Concentrate 100 mL to under 20 mL in just a few minutes or concentrate one liter 50 times in less than 60 minutes. Alternatively, speed up your process by using two Vivaflow® 50R units in parallel and concentrate 1 liters in under 30 min.

Vivaflow® 50R is a plug-and-play laboratory crossflow cassette for concentrating up to 1 L aqueous samples. The active membrane area per device is 50 cm².

One unit comes with all the necessary accessories for running the device with a laboratory pump and a size 16 pump head. For speeding up concentration, two cassettes can be run simultaneously.

- Fast and easy protein sample concentration
- Reusable
- Concentrates volumes from 0.1 L to 1 L
- Optimal for concentration of culture supernatants and viruses
- The most compact crossflow cassette with a premium Hydrosart® membrane

Specifications

Vivaflow® 50R

Dimensions	Overall L×H×W Channel W×H Active membrane area Hold-up volume (module) Min. recirculation volume	100 mm × 100 mm × 24 mm 7.5 mm × 0.4 mm 50 cm ² 1.7 mL 10 mL
Operating conditions	Non-recoverable hold-up Pump flow	< 0.5 mL 200 mL/min to 400 mL/min
operating conditions	Maximum pressure Maximum temperature	4 bar (60 psi) 60°C
Materials of construction	Main housing Flow channel Membrane support Seals and O-rings Pressure indicator Flow restrictor Fittings Tubing	Acrylic Acrylic Polypropylene Silicone Polypropylene, SS spring Polypropylene Nylon PVC (medical grade)

Performance Characteristics

	Time to concentrate up to $20 \times$ [min.] at 3.0 bar inlet 2.5 bar outlet pressure, 20° C			
	Start volume	Average flux	Recovery %	
	250 mL	mL/min	Direct	25 mL rinse
Lysozyme 0.25 mg/mL (14,000 MW)				
5,000 MWCO Hydrosart [®]	70	3.4	96%	98%
10,000 MWCO Hydrosart®	23	10.3	94%	96%
BSA 1.0 mg/mL (66,000 MW)				
10,000 MWCO Hydrosart®	24	9.9	98%	> 99%
30,000 MWCO Hydrosart®	15	15.8	97%	> 99%
γ Globulins 1.0 mg/mL (150,000 MW)	46	E 2	0.70/	> 000/
100,000 MWCO Hydrosart®	46	5.2	97%	> 99%

37

Performance Characteristics

	Time to concentrate up to 20× [min.] at 3.0 bar inlet 2.5 bar outlet pressure, 20°C			
	Start volume	Average flux	Recovery %	
	250 mL	mL/min	Direct	25 mL rinse
Start volume 1 L (one Vivaflow® 50R at 3 bar) 10,000 MWCO Hydrosart®				
BSA 1.0 mg/mL	95	10.0	98%	> 99%
Start volume 1 L (two Vivaflow® 50R in paralle at 3 bar) 10,000 MWCO Hydrosart®	I			
BSA 1.0 mg/mL	48	19.8	98%	> 99%

Ordering Information

Vivaflow® 50R*	Qty./Pkg.	Prod. No.
5,000 MWCO Hydrosart®	1	VF05H1
10,000 MWCO Hydrosart [®]	1	VF05H0
30,000 MWCO Hydrosart [®]	1	VF05H2
100,000 MWCO Hydrosart®	1	VF05H4

^{*)} Vivaflow® 50R modules include pressure indicator, flow restrictor and size 16 pvc peristaltic tubing and fittings.

Vivaflow[®] 50R Complete Set of Accessories

Pump (230 V), Easy Load pump head (size 16), tubing, 500 mL sample diafiltration reservoir	1	VFS202
Pump (115 V), Easy Load pump head (size 16), tubing,	1	VFS204
500 mL sample diafiltration reservoir		

Tubing and Fittings

Size 16 pvc pump tubing and Luer fittings	1	VFA004
(3 m, 3.2 × 1.6 mm)		
T-connectors for running 2 units	2	VFA030
Flow restrictor set (2 × 0.4 mm, 0.6 mm, 0.8 mm)	6	VFA009
Female luer fittings – size 16	10	VFA032
Flow restrictors 0.6 mm	6	VFA035
Female luer fittings – size 15	10	VFA036

Accessories


Masterflex economy drive variable speed peristaltic pump (230 V)	1	VFP001
Masterflex economy drive variable speed peristaltic pump (115 V)	1	VFP002
500 mL sample and or diafiltration reservoir	1	VFA006
Masterflex Fasy Load nump head – size 16	1	VFA012

Visit us at www.sartorius.com/Vivaflow50R to get additional info. Here you can find instructions on how to use Vivaflow® 50R for

- Preparation of biological nanoparticles and medical nanocarriers
- Concentration and purification of viruses

Vivaflow® 200



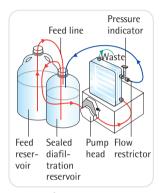

Vivaflow® 200 setup for diafiltration

500 mL to 5 L Samples

Concentrate 250 mL to under 20 mL in just a few minutes or concentrate one liter 50 times in less than 30 minutes.

Alternatively, use two Vivaflow® 200 units in parallel and concentrate 5 liters in under 75 minutes.

Nearly total sample recoveries can be expected with most solutions.

The standard economy package comes complete with tubing, pressure indicator, flow restrictor and high-pressure pump tubing. All you need is a peristaltic pump capable of handling 6.4 mm OD (size 16) tubing. Should your pump head require larger tubing, just use the interconnecter provided to attach your own peristaltic tubing to this standard product.

Two modules in parallel will concentrate 5 liters in under 75 minutes

Specifications

Vivaflow 200

Dimensions	Overall L×H×W	126 mm × 138 mm × 38 mm
J	Channel W×H	10 mm×0.4 mm
	Active membrane area	200 cm ²
	Hold-up volume (module)	5.3 mL
	Min. recirculation volume	< 20 mL
	Non-recoverable hold-up	< 2 mL
Materials of construction	Main housing	Acrylic
	Flow channel	Acrylic
	Membrane support	Polypropylene
	Seals and O-rings	Silicone
	Pressure indicator	Polypropylene,
		SS spring
	Flow restrictor	Polypropylene
	Fittings	Nylon
	Tubing	PVC (medical grade)
Operating conditions	Pump flow	200 mL/min to 400 mL/min
	Maximum pressure	4 bar (60 psi)
	Maximum temperature	60°C

39

Performance Characteristics

		entrate up to 2	20× [min.] at	
	3 bar inlet pr	3 bar inlet pressure, 20°C		
	Start volume	Average flux	Recovery %	
	1 L	mL/min	Direct	25 mL rinse
Lysozyme 0.25 mg/mL				
(14,000 MW)				
2,000 MWCO Hydrosart®	160	6	97%	> 99%
3,000 MWCO PES	180	5	97%	> 99%
BSA 1.0 mg/mL (66,000 MW)				> 99%
5,000 MWCO PES	29	33	98%	> 99%
5,000 MWCO Hydrosart®	70	14	98%	> 99%
10,000 MWCO PES	23	41	96%	> 99%
10,000 MWCO Hydrosart®	35	27	98%	
30,000 MWCO PES	25	38	96%	99%
30,000 MWCO Hydrosart®	20	48	96%	> 99%
50,000 MWCO PES	22	43	96%	98%
γ Globulins 1.0 mg/mL				
(average 160,000 MW)				
100,000 MWCO PES	54	18	96%	99%
Yeast 1.0 mg/mL				
(S. Cerevisiae)				
0.2 μm PES	11	86	92%	98%
Dilute solute concentration,			32 70	3070
start volume 1 L at 3 bar,				
10,000 MWCO PES				
BSA 0.001 mg/mL	18	52	90%	98%
BSA 0.001 mg/mL	20	47	92%	98%
BSA 0.01 mg/mL	21	45	94%	99%
	21	45	94%	99%
Start volume 5 L				
(two Vivaflow® 200				
in parallel at 3 bar)				
10,000 MWCO PES				
BSA 1.0 mg/mL (66,000 MW)	67	70	97%	> 99%

Ordering Information

Vivaflow® 200*	Qty./Pkg.	Prod. No.
3,000 MWCO PES	1	VF20P9
5,000 MWCO PES	1	VF20P1
10,000 MWCO PES	1	VF20P0
30,000 MWCO PES	1	VF20P2
50,000 MWCO PES	1	VF20P3
100,000 MWCO PES	1	VF20P4
0.2 μm PES	1	VF20P7
2,000 MWCO Hydrosart®	1	VF20H9
5,000 MWCO Hydrosart®	1	VF20H1
10,000 MWCO Hydrosart®	1	VF20H0
30,000 MWCO Hydrosart®	1	VF20H2
100,000 MWCO Hydrosart®	1	VF20H4

^{*)} Vivaflow® 200 modules include pressure indicator, flow restrictor and size 16 PVC peristaltic tubing and fittings.

Vivaflow® 200 Complete Set of Accessories

Pump (230 V), Easy Load pump head (size 16), tubing,	1	VFS202
500 mL sample diafiltration reservoir		VI 3202
Pump (115 V), Easy Load pump head (size 16), tubing,	1	VFS204
500 mL sample diafiltration reservoir		

Accessories

Masterflex economy drive variable speed peristaltic pump (240 V)	VFP001
Masterflex economy drive variable speed peristaltic pump (115 V)	VFP002
500 mL sample and or diafiltration reservoir	VFA006
Masterflex Easy Load pump head – size 16	VFA012
Masterflex Easy Load pump head – size 15	VFA013

Tubing and Fittings

Size 15 PVC pump tubing and Luer fittings (3 m, 4.8 mm × 2.6 mm)	VFA003
Size 16 PVC pump tubing and Luer fittings (3 m, 3.2 mm × 1.6 mm)	VFA004
Y-connector (size 15 to 2×size 16)	VFA005
Flow restrictor set (2 × 0.4 mm, 0.6 mm, 0.8 mm)	VFA009
Female luer fittings, size 16 (10 units)	VFA032
Flow restrictors 0.6 mm (6 units)	VFA035
Female luer fittings, size 15 (10 units)	VFA036

Visit us at www.sartorius.com/Vivaflow200 to get additional info. Find instructions on how to use Vivaflow® 200 for

- The measurement of soluble trace metals in seawater
- The workflow in protein research laboratories
- Preparation of biological nanoparticles and medical nanocarriers
- Concentration and purification of viruses
- Concentrating hybridoma supernatants prior to to affinity chromatography

Vivapore® Solvent Absorption Concentrators





Vivapore® 5



Vivapore® 10 | 20

3 mL to 20 mL Samples

With no need for additional equipment, pressure or vacuum, solvent absorption is the most economic and user-friendly concentration technique available to the clinician and research scientist.

Just fill the unit with the solution to be concentrated, wait for the desired concentration level to be achieved and then pipette the concentrated sample from the bottom of the reservoir.

Vivapore®* is ideal for general-purpose laboratory concentration and purification prior to further analysis. It is particularly suited for labile solutions that can denature with alternative shear- or pressure-inducing methods or that require processing in a cold room environment.

Vivapore® concentrators extend the solvent absorption technique to a totally new level of performance, application potential and ease of use.

Specifications

Vivapore® 5	Vivapore® 10 20
PES	PES
7,500	7,500
20 cm ²	28 cm ²
SAN	SAN
1 mL to 5 mL	2 mL to 10 mL 20 mL*
50 μL	50 μL
42	46
	7,500 20 cm ² SAN 1 mL to 5 mL

Typical Performance Characteristics

	Time to concentrate up to 10× [min.]			Concentrate recovery %		very %
Product	VP5	VP10 20	VP10 20*	VP5	VP10 20	VP10 20*
Start volume	5 mL	10 mL	20 mL	5 mL	10 mL	20 mL
Cytochrome c (12,600 MW)	0.25 mg/mL	0.25 mg/mL	0.1 mg/mL	0.25 mg/mL	0.25 mg/mL	0.1 mg/mL
7,500 MWCO PES	35	75	150	90%	90%	92%
BSA (66,000 MW) 7,500 MWCO PES	30	55	115	92%	92%	92%
lgG (160,000 MW) 7,500 MWCO PES	40	70	160	75%	77%	78%

^{*} with additional reservoir

Performance Characteristics

· ciioiiiiaiice diiaia		•				
		to concentrate 50× [min.]		Concentrate recovery %		ery %
Product	VP5	VP10 20	VP10 20*	VP5	VP10 20	VP10 20*
Start volume	5 mL	10 mL	20 mL	5 mL	10 mL	20 mL
Cytochrome c (12,600 MW) 7,500 MWCO PES	65	70	160	91%	88%	90%
BSA (66,000 MW) 7,500 MWCO PES	45	50	105	90%	90%	92%
IgG (160,000 MW) 7,500 MWCO PES	50	65	140	53%	65%	74%


^{*} with additional reservoir

Ordering Information

Vivapore® 5*	Qty./Pkg.	Prod. No.
7,500 MWCO PES	4	VP0503
7,500 MWCO PES	30	VP0501
* includes stand and recovery pipettes		
Requires Stand		
7,500 MWCO PES	100	VP0502
Vivapore® 10 20*		
7,500 MWCO PES	4	VP2003
7,500 MWCO PES	30	VP2001
* includes stand and recovery pipettes		
Requires Stand		
7,500 MWCO PES	100	VP2002
Accessories		
Disposable stands for 4 units	6	VPA002
Plastic recovery pipettes (Vivapore® 10 20)	100	VPA005
10 mL expansion reservoir (Vivapore® 10 20)	10	VPA006
Plastic recovery pipettes (Vivapore® 5)	100	VPA007

Ultrafiltration Membrane Filters

PES 146, CTA 145 and Hydrosart® 144

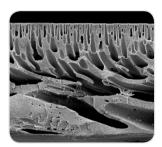

Polyethersulfone (PES)

This is a general purpose membrane that provides excellent performance with most solutions when retentate recovery is of primary importance. Polyethersulfone membranes exhibit no hydrophobic or hydrophilic interactions and are usually preferred for their low fouling characteristics, exceptional flux and broad pH range.

Cellulose Triacetate (CTA)

High hydrophilicity and exceptionally low non-specific binding characterize this membrane. Cast without any membrane support that could trap or bind passing microsolutes, these membranes are to be ideal for sample cleaning and protein removal and when high recovery of the filtrate solution is of primary importance.

Hydrosart®

These membranes are also highly hydrophilic and are often preferred for their high protein recovery when processing very dilute solutions. Resistance to autoclaving, ease of cleaning and extended chemical resistance also characterize this type of membrane.

Specifications

Specifications for Polyethersulfone, Type 146

•		
Thickness	120 μm	
pH range	1–14	
Water flux	MWCO 10,000	0.2 mL/min/cm ²
Protein retention	Cytochrome C	95%

Specifications for Cellulose Triacetate, Type 145

Thickness	120 μm	
pH range	4–8	
Water flux	MWCO 10,000	0.11 mL/min/cm ²
Protein retention	Cytochrome C	90%

Specifications for Hydrosart[®], Type 144

. ,		
Thickness	180 μm	
pH range	1–13	
Water flux	MWCO 10,000	0.08 mL/min/cm ²
Protein retention	Cytochrome C	99%

Ordering Information

Polyethersulfone Membrane Filters, Type 146

Diameter in mm	MWCO	Qty./Pkg.	Prod. No.
47	1,000 daltons	10	1460947D
63	1,000 daltons	10	1460963D
76	1,000 daltons	10	1460976D
25	5,000 daltons	10	1462925D
47	5,000 daltons	10	1462947D
63	5,000 daltons	10	1462963D
76	5,000 daltons	10	1462976D
25	10,000 daltons	10	1463925D
63	10,000 daltons	10	1463963D
76	10,000 daltons	10	1463976D
150	10,000 daltons	10	14639-150D
25	30,000 daltons	10	1465925D
63	30,000 daltons	10	1465963D
76	30,000 daltons	10	1465976D
25	50,000 daltons	10	1465025D
47	50,000 daltons	10	1465047D
76	50,000 daltons	10	1465076D
25	300,000 daltons	10	1467925D
47	300,000 daltons	10	1467947D
76	300,000 daltons	10	1467976D

Cellulose Triacetate Membrane Filters, Type 145

Diameter in mm	MWCO	Qty./Pkg.	Prod. No.
25	5,000 daltons	10	1452925D
47	5,000 daltons	10	1452947D
25	10,000 daltons	10	1453925D
47	10,000 daltons	10	1453947D
50	10,000 daltons	10	1453950D
25	20,000 daltons	10	1454925D
43	20,000 daltons	10	1454943D
47	20,000 daltons	10	1454947D
47	20,000 daltons	100	1454947N
63	20,000 daltons	10	1454963D

45

Hydrosart® Membrane Filters, Type 144

Diameter in mm	MWCO	Qty./Pkg.	Prod. No.
25	5,000 daltons	10	1442925D
44	5,000 daltons	10	1442944D
63	5,000 daltons	10	1442963D
76	5,000 daltons	10	1442976D
25	10,000 daltons	10	1443925D
47	10,000 daltons	10	1443947D
63	10,000 daltons	10	1443963D
76	10,000 daltons	10	1443976D
25	30,000 daltons	10	1445925D
47	30,000 daltons	10	1445947D
63	30,000 daltons	10	1445963D
76	30.000 daltons	10	1445976D



Vivacon® 500

For DNA Sample Desalting and Concentration



Reproducible DNA and Protein Sample Desalting and Concentration

Vivacon® 500 centrifugal concentrators offer the optimal solution for DNA and protein concentration and buffer exchange applications. For optimal performance with highly dilute samples, e.g. forensic samples, Vivacon® 500 incorporates the patented regenerated cellulose membrane Hydrosart®.

High recoveries and excellent reproducibilities are combined with convenience offered by the molecular weight cutoff printed on the individual Vivacon® 500 units.

As Vivacon® 500 can be respun after sample processing, this ensures complete concentrate recovery, which is especially important when working with low-sample concentrations.

Vivacon® 500-PCR Grade

To use DNA amplification technologies, any traces of DNA originating from the equipment need to be eliminated.

Vivacon® 500-PCR Grade units are treated with ethylene oxide (ETO) in a validated process to deactivate all traces of DNA that might interfere with subsequent amplification procedures.

References: K. Shaw et al., Int. J. Legal Med. (2008) 122: 29–33

Specifications

Vivacon® 500

VIVACUII 300			
Concentrator capacity	Fixed-angle rotor	0.5 mL	
Dimensions	Total length (concentration)	45 mm	
	Total length (backspin)	47.5 mm	
	Width	12.4 mm	
	Active membrane area	0.32 cm ²	
	Hold-up volume of membrane		
	and support	< 5 μL	
	Dead-stop volume (40° rotor)	5 μL	
Materials of construction	Body	Polycarbonate	
	Filtrate vessel	Polypropylene	
	Membrane	Hydrosart [®] ,	
		Cellulose Acetate	

Conversion Table for Hydrosart® MWCO to Nucleotide Cutoff

Membrane	MWCO	Double-Stranded Nucleotide Cutoff (bp)
Hydrosart®	2 kDa	> 10
Hydrosart®	10 kDa	> 30
Hydrosart®	30 kDa	> 50
Hydrosart®	50 kDa	> 300
Hydrosart®	100 kDa	> 600
Cellulose Acetate	125 kDa	> 650

Performance Characteristics for DNA

Start volume 0.5 mL, sample concentration 50 mg/mL.

	Sample size (bp)	Time to concentrate up to 30× [min.] at 20°C	Concentrate recovery %	g-force (xg)
2,000 MWC0	10	60 min.	93%	7,500
10,000 MWCO	30	25 min.	94%	7,500
30,000 MWCO	50	18 min.	88%	5,000
50,000 MWCO	300	18 min.	91%	5,000
100,000 MWCO	600	10 min.	87%	3,000
125,000 MWCO	650	12 min.	85%	2,000
125,000 MWCO	900	9 min.	94%	3,000

Performance Characteristics for Proteins

Start volume 0.5 mL, sample and concentration of proteins as specified in table.

	Sample	Time to concentrate up to 30× [min.] at 20°C	Concentrate recovery %	g-force (xg)
2,000 MWC0	0.25 mg/mL cytochrome c	30 min.	95%	14,000
10,000 MWCO	0.25 mg/mL cytochrome c	15 min.	92%	14,000
30,000 MWCO	1.0 mg/mL BSA	10 min.	95%	14,000
50,000 MWCO	1.0 mg/mL BSA	10 min.	92%	14,000
100,000 MWCO	1.0 mg/mL bovine lgG	11 min.	90%	8,000
125,000 MWC0	1.0 mg/mL bovine lgG	10 min.	81%	8,000

Ordering Information

Vivacon® 500	Oty./Pkg.	Prod. No.
2,000 MWCO	25	VN01H91
2,000 MWCO	100	VN01H92
10,000 MWCO	25	VN01H01
10,000 MWCO	100	VN01H02
30,000 MWCO	25	VN01H21
30,000 MWCO	100	VN01H22
50,000 MWCO	25	VN01H31
50,000 MWCO	100	VN01H32
100,000 MWCO	25	VN01H41
100,000 MWCO	100	VN01H42
125,000 MWCO	25	VN01H81
125,000 MWCO	100	VN01H82

49

Vivacon® 500	Qty./Pkg.	Prod. No.
Sample kit L (4 units each of 2 kDa, 10 kDa, 30 kDa)	12	VN01HL12
Sample kit H (4 units each of 30 kDa, 50 kDa, 100 kDa)	12	VN01HH12
Vivacon® 500-PCR Grade	Qty./Pkg.	Prod. No.
30,000 MWCO	25	VN01H21ETO
30,000 MWCO	100	VN01H22ETO
50,000 MWCO	25	VN01H31ETO
50,000 MWCO	100	VN01H32ETO
100,000 MWCO	25	VN01H41ETO
100,000 MWCO	100	VN01H42ETO
125,000 MWCO	25	VN01H81ETO
125,000 MWC0	100	VN01H82ETO
Accessories	Qty./Pkg.	Prod. No.
Collection tube	100	VNCT01

Visit us at www.sartorius.com/Vivacon500 to get additional info. Find instructions on how to use $Vivacon^{\circ}$ 500 for

- Primer removal after a PCR reaction
- Filter aided sample preparation (FASP) for proteomic analysis by mass spectrometry

Vivacon[®] 2

For DNA Sample Desalting and Concentration



Reproducible DNA Sample Desalting and Concentration

Vivacon® 2 centrifugal concentrators offer the optimal solution for DNA and protein concentration and buffer exchange applications. For optimal performance with highly dilute samples, e.g. forensic samples, Vivacon® 2 incorporates the patented regenerated cellulose membrane Hydrosart®.

High recoveries and excellent reproducibilities are combined with the convenience provided by the volume

graduation and molecular weight cutoff printed on the individual Vivacon® 2 units.

As Vivacon® 2 can be respun after sample processing, this ensures complete concentrate recovery, which is especially important when working with low-sample concentrations.

Vivacon® 2-PCR Grade

Vivacon® 2-PCR Grade units are treated with ethylene oxide (ETO) in a validated process to deactivate all traces of DNA that might interfere with subsequent amplification procedures.

Specifications

Concentrator capacity	Fixed-angle rotor	2 mL
Dimensions	Total length (concentration)	125 mm
	Total length (backspin)	115 mm
	Width	16 mm
	Active membrane area	0.95 cm ²
	Hold-up volume membrane and support	10 μL
	Dead-stop volume (25° rotor)	55 μL
Materials of construction	Body	Polycarbonate
	Filtrate vessel	Polypropylene
	Backspin vial	Polypropylene
	Concentrator cap	Polypropylene
	Membrane	Hydrosart [®] ,
		Cellulose Acetat

Conversion Table for Hydrosart® MWCO to Nucleotide Cutoff

Membrane	MWCO	Double-Stranded Nucleotide Cutoff (bp)
Hydrosart®	2 kDa	> 10
Hydrosart®	10 kDa	> 30
Hydrosart®	30 kDa	> 50
Hydrosart®	50 kDa	> 300
Hydrosart®	100 kDa	> 600
Cellulose Acetate	125 kDa	> 650

Performance Characteristics for DNA

Start volume 2 mL, sample concentration 50 mg/mL.

	Sample size (bp)	Time to concentrate up to 30× [min.] at 20°C	Concentrate recovery %	g-force (xg)
2,000 MWC0	10	120 min.	92%	7,500
10,000 MWC0	30	60 min.	94%	5,000
30,000 MWC0	50	60 min.	95%	2,500
50,000 MWC0	300	45 min.	96%	2,500
100,000 MWCO	600	30 min.	93%	2,500
125,000 MWCO	650	30 min.	88%	2,500
125,000 MWCO	900	30 min.	89%	2,500

51

Performance Characteristics for Proteins

Start volume 2 mL, sample and concentration of proteins as specified in table.

	Sample	Time to concentrate up to 30× [min.] at 20°C	Concentrate recovery %	g-force (xg)
2,000 MWC0	0.25 mg/mL cytochrome c	120 min.	95%	7,500
10,000 MWCO	0.25 mg/mL cytochrome c	90 min.	96%	5,000
30,000 MWC0	1.0 mg/mL BSA	40 min.	96%	5,000
50,000 MWC0	1.0 mg/mL BSA	30 min.	94%	5,000
100,000 MWC0	1.0 mg/mL bovine lgG	30 min.	92%	5,000
125,000 MWC0	1.0 mg/mL bovine lgG	27 min.	81%	5,000

Ordering Information


Vivacon® 2	Qty./Pkg.	Prod. No.
2,000 MWCO	25	VN02H91
2,000 MWCO	100	VN02H92
10,000 MWCO	25	VN02H01
10,000 MWCO	100	VN02H02
30,000 MWC0	25	VN02H21
30,000 MWC0	100	VN02H22
50,000 MWC0	25	VN02H31
50,000 MWC0	100	VN02H32
100,000 MWCO	25	VN02H41
100,000 MWCO	100	VN02H42
125,000 MWCO	25	VN02H81
125,000 MWCO	100	VN02H82

Vivacon® 2-PCR Grade	Qty./Pkg.	Prod. No.
30,000 MWCO	25	VN02H21ETO
30,000 MWCO	100	VN02H22ETO
50,000 MWCO	25	VN02H31ETO
50,000 MWCO	100	VN02H32ETO
100,000 MWCO	25	VN02H41ETO
100,000 MWCO	100	VN02H42ETO
125,000 MWC0	25	VN02H81ETO
125,000 MWCO	100	VN02H82ETO



Vivapure[®]

Ion Exchange Protein Purification Products

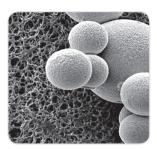

Chromatography gel beads (right) are shown on top of a membrane adsorber in this SEM. The membrane adsorber pores are more than 50 times larger than bead pores.

Fast and Easy-to-Use Spin Columns

Vivapure® Ion Exchange (IEX) spin columns are centrifugal devices that incorporate Sartobind® membrane adsorber technology as their chromatography matrix. Vivapure® IEX spin columns make protein purification as easy as filtration. The devices are ready to use and eliminate the risk of running dry. For many protein purification applications, they replace time-consuming and tedious column chromatography.

The rapid Vivapure® IEX 1-2-3 bind-wash-elute protocol especially lends itself to screening applications, where many different samples are processed in parallel.

The Sartobind® Membrane Adsorber Matrix

Sartobind® IEX membrane adsorbers are based on stabilized regenerated cellulose and exhibit a microporous structure with a pore size of $> 3 \mu m$, which is orders of magnitude larger than conventional chromatographic gel materials. This allows molecules to be transported to the ligands immobilized on the membrane adsorber by convective flow, resulting in exceptionally high flow rates.

By contrast, gel chromatography is slowed down due to diffusion limitations, as the molecules need to enter the small bead pores in order to be bound by the ligands. The porous membrane adsorber enables fast, reproducible and scalable protein purification.

Fast and Simple-to-Use Spin Columns

- Devices are ready to use
- They make protein purification as simple as filtration

Reproducible Results

- No column packing necessary as the devices are ready to use
- Membrane adsorber spin columns cannot crack or run dry

Centrifugal Devices


- Offer the possibility of working in parallel

Low Bed Volume

 Small membrane adsorber bed volumes allow working with lower buffer amounts, resulting in concentrated elution fractions

Scalable Product Range

 Process-scale modules are available with the same Sartobind® IEX membrane adsorber matrix

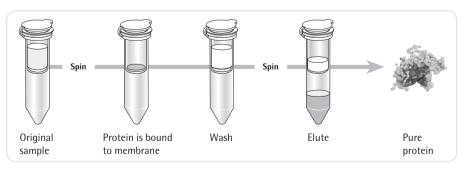

Fast and easy protein purification with Vivapure® spin columns



Vivapure® Mini-400 | 500 μL Binding capacities: 1 mg to 4 mg



Vivapure® Maxi-19 | 20 mL Binding capacities: 15 mg to 80 mg

Specifications

Available Formats

Vivapure® IEX Products	Application
Vivapure® Mini spin columns	Sample fractionationPurification condition scoutingSmall-scale purification
Vivapure® Maxi spin columns	Large-scale sample fractionationOne-step protein purification concentrationPolishing of His-tagged protein

Membrane Availability

Functional Groups	Ion Exchanger Type	
Sulfonic acid (S)	Strong acidic cation exchanger:	R-CH ₂ -SO ₃ -Na+
Quaternary ammonium (Q)	Strong basic anion exchanger:	R-CH ₂ -N ⁺ -(CH ₃) ₃ Cl ⁻
Diethylamine (D)	Weak basic anion exchanger:	R-CH ₂ -NH ⁺ -(CH ₂ H ₅) ₂

Performance Characteristics

Vivapure® Spin Columns	Protein Binding Capacity* (mg)	Max. Volume per Centrifuge Run Using a Swing-Out Rotor (ml)	Max. Volume per Centrifuge Using a Fixed-Angle Rotor Run (ml)
Vivapure [®] Mini H	4	0.4	
Vivapure® Maxi H	60 to 80	19	10.5

^{*} Actual yields depend on specific protein sample and selected pH and salt conditions. Yields established using 1 mg/mL BSA in 25 mM Tris/HCL pH 8.0 with Vivapure® Q € D spin columns and 1 mg/mL cytochrome c in 25 mM sodium acetate buffer pH 5.5 with Vivapure® S spin columns.

Ordering Information

Description	Spin Columns	Centrifuge Tubes	Prod. No.						
Vivapure® Mini Ion Exchange Spin Columns (up to 0.5 mL)									
Vivapure® Mini S&Q H starter kit	16	32	VS-IX01SQ16						
Vivapure® D Mini H	24	48	VS-IX01DH24						
Vivapure® Q Mini H	24	48	VS-IX01QH24						
Vivapure® S Mini H	24	48	VS-IX01SH24						
Vivapure® Maxi Ion Exchange Spi	n Columns (up	to 20 mL)							
Vivapure® D Maxi H	8	16	VS-IX20DH08						
Vivapure® Q Maxi H	8	16	VS-IX20QH08						
Vivapure [®] S Maxi H	8	16	VS-IX20SH08						

Virus Purification and Concentration

Vivapure® Virus Purification and Concentration Kits



Recombinant virus vectors are the preferred method for a wide range of gene delivery applications. Especially adenovirus type 5 and VSV-G pseudotyped lentivirus are two frequently utilized viral vectors for in vitro and in vivo applications.

Recombinant Adenovirus Vectors

Recombinant adenovirus vectors are versatile tools in research and therapeutic applications for gene transfer and protein expression in cell lines that have low transfection efficiency with liposomes. After entering cells, the virus remains epichromosomal – i.e., does not integrate into the host chromosome, leaving the host genome unaffected. The delivery of RNAi into cells is becoming a major application for adenovirus vectors.

Lentivirus Vectors

Lentivirus vectors are frequently used in gene transfer studies, due to their ability of gene transfer and integration into dividing and non-dividing cells. The pseudotyped envelope with vesicular stomatitis virus envelope G (VSV-G) protein broadens their target cell range. Lentiviral vectors have been shown to deliver genes into cell types (e.g. neurons, lymphocytes and macrophages) which other retrovirus vectors could not be used for. The lentivirus vector is increasingly used to integrate siRNA efficiently in a wide variety of cell lines and primary cells, both in vitro and in vivo.

Rapid Virus Purification by Membrane Chromatography

The Sartobind[®] ion exchange membrane adsorber technology used in AdenoPACK and LentiSELECT is unique in its capability to efficiently and rapidly capture and recover large virus particles. Compared with chromatography media, membrane adsorbers provide large 3,000 nm pores, allowing unrestricted access and recovery of virus from the charged adsorber surface. Convective flow through the syringe filter devices provides high-speed separations not possible with traditional chromatography, cesium chloride density gradients and ultracentrifugation methods. Sartorius membrane adsorbers with porous matrices, high capacities, low differential pressures, high flow rates and low unspecific adsorption show excellent performance in small-scale virus purification. In addition, these syringe filter devices are scalable and comply with cGMP requirements for large-volume, high-performance separation, reducing final process time ten-fold.

Adenovirus Purification with Vivapure® AdenoPACK Kits

AdenoPACK 20 100 500

The AdenoPACK adenovirus purification and concentration kits offer researchers who need to recover up to 3×10^{13} purified recombinant adenovirus particles for in vitro transfection a fast, safe and easy-to-use solution. The kits include all reagents and devices necessary for clarification, purification and concentration of adenovirus type 5 from HEK293 cell cultures – all within just two hours. These straightforward kits replace time-consuming and labor-intensive 48-hour CsCl density gradients.

AdenoPACK kits are offered as AdenoPACK 20, AdenoPACK 100 and AdenoPACK 500 for the purification and concentration of adenovirus type 5 from 20 mL to 500 mL cell cultures, resulting in 1 × 10¹¹ to 3 × 10¹³ purified viral particles. For each sample volume, the most convenient handling method is provided for ultimate convenience.

To this end, preparations using AdenoPACK 20 are supplied in a spin column format for centrifuges. AdenoPACK 100 is a manually operated kit in a syringe filter format* and AdenoPACK 500 is a pump-driven kit.

* Vivapure® AdenoPACK 100 can be alternatively be operated with a laboratory pump or an infusion pump, for which protocols are provided on our web page at www.sartorius-stedim.com. Additionally, the tubes and adaptors needed for these operating modes can be ordered.

AdenoPACK Advantages Fast and Easy Virus Purification

- Purification completed in just 2 hours
- Convenient, over 10 x faster alternative to CsCl density gradient

Quantitative Yields

 In contrast to CsCl density gradient, the complete cell culture is used for virus purification and not only the viral pellet

Flexible Product Range

 Applicable from initial construct screening to large-scale virus production

Complete Kit

 Including filtration devices, AdenoPACK units for virus purification, Vivaspin® and all buffers

Low Endotoxin Levels

 High cell viability and infection rates due to endotoxin levels of < 0.025 EU/mL

Specifications

Adenovirus Purification Kit Specifications

Product	AdenoPACK 20	AdenoPACK 100	AdenoPACK 500
Sample size	20 mL cell culture	20 mL to 200 mL of cell culture	500 mL of cell culture
Number of purifications	6×20 mL	2×20 mL to 60 mL 1×200 mL	1 × 500 mL
Virus particles (VP) per mL	Typically up to 1×10^{11} - 10^{12}	Typically up to 1×10^{13}	Typically up to 3×10^{13}
VP/IU	50 to 100	20 to 50	20 to 50
Processing time	Typically one hour	Typically two hours	
Endotoxin level	<0.025 EU/mL	<0.025 EU/mL	<0.025 EU/mL



Ordering Information

Vivapure® AdenoPACK 20

Vivapure® AdenoPACK 20	VS-AVPQ020	
Vivapure® AdenoPACK 20 RT*	VS-AVPQ022	

^{*} AdenoPACK 20 RT does not contain Benzonase®***



Vivapure® AdenoPACK 100

Vivapure® AdenoPACK 100	VS-AVPQ101
Vivapure® AdenoPACK 100 RT*	VS-AVPQ102

AdenoPACK 100 Accessories

Pump tubing set for	VS-AVPA001	
Vivapure® AdenoPACK 100		

^{*} AdenoPACK 100 RT does not contain Benzonase®**



Vivapure® AdenoPACK 500

Vivapure® AdenoPACK 500	VS-AVPQ501
Vivapure® AdenoPACK 500 RT*	VS-AVPQ502

^{*} AdenoPACK RT-kits do not contain Benzonase®***

^{**} Benzonase® is a registered trademark of Merck

Lentivirus Purification with the Vivapure® LentiSELECT Kit

LentiSELECT 40 500 1000

The LentiSELECT kits for lentivirus purification and concentration offer researchers who need to recover up to 5×10^9 infective lentivirus particles per mL for in vitro transfection or animal studies a fast and easy-to-use solution.

These straightforward kits replace time-consuming ultracentrifugation protocols, which typically take approximately one day for large sample volumes. Vivapure® LentiSELECT thus reduces purification time to just a few hours.

LentiSELECT kits are offered as LentiSELECT 40, LentiSELECT 500 and LentiSELECT 1000 for the purification and concentration of VSV-G pseudotyped lentivirus from 40 mL to 1,000 mL cell cultures, resulting in 8×10^8 to 1×10^{10} purified infective particles. The most convenient handling method is provided for each sample volume. To this end, 40 mL sample volumes are processed manually with LentiSELECT 40, while LentiSELECT 500 and 1000 are pump-driven kits.

LentiSELECT Advantages

Fast and Easy Virus Purification

- Purification completed in less than one to six hours, depending on sample volume
- Kit is as easy to use as filtration

No Need for Expensive Instruments

 Lentivirus purification with LentiSELECT is independent of equipment, such as ultracentrifuges

High Virus Purity

 Achieve pure virus based on a chromatographic method for your experiments instead of a crude and variable cell culture supernatant pellet

Optimal for Multiple Virus Construct Screening

With LentiSELECT 40, four purification runs can be conducted in parallel with one kit

Complete Kits

 Including LentiSELECT units for virus purification, Vivaspin® units for concentration | buffer exchange and all buffers and syringes necessary

Low Endotoxin Levels

 High cell viability and infection rates due to endotoxin levels of < 0.025 EU/mL

Specifications

Product	LentiSELECT 40	LentiSELECT 500	LentiSELECT 1000
Sample size	40 mL cell culture	500 mL of cell culture	1,000 mL of cell culture
Number of purifications	4 × 40 mL	1×500 mL	1×1,000 mL
Virus particles (VP) per mL	Typically up to 3×10^9	Typically up to 2-5×10 ⁹	Typically up to 4-6×10 ¹³
VP/IU	5 to 15	5 to 15	20 to 50
Processing time	Typically up to 45 min	Typically up to 3 hours	Typically up to 6 hours
Endotoxin level	<0.025 EU/mL	<0.025 EU/mL	<0.025 EU/mL



Ordering Information

Vivapure® LentiSELECT 40

Vivapure® LentiSELECT 40 VS-LVPQ040



Vivapure® LentiSELECT 500

Vivapure® LentiSELECT 500 VS-LVPQ500



Vivapure® LentiSELECT 1000

Vivapure® LentiSELECT 1000 VS-LVPQ1000





Filtration Devices

Table of Contents

- 64 Minisart® Syringe Filters
- 80 Claristep® Syringeless Filtration
- 86 Sartolab® P20 and Sartolab® P20 Plus
- 89 Sartolab® RF|BT
- 92 NEW! Sartolab® RF 50
- 94 NEW! Sartolab® Multistation
- 96 Sartoclear Dynamics® Lab P15
- 98 NEW! Sartoclear Dynamics® Lab V
- 102 Chemical Compatibility



Syringe filters are used for many routine preparation steps in laboratories all over the world. They are convenient, ready-to-use disposables for sterile filtration of liquids and removal of particles from solutions and gases. Depending on the reagents filtered, syringe filters have to fulfill certain requirements to best serve customer's application. Sartorius offers Minisart® syringe filters and filters optimized for a wide range of relatively large volumes. The filters are clean and safe as they are virtually free of leachables and extractables and reliably remove particles and microorganisms with no leakage. If you need to rely on the quality of your filtrate whether it needs to be sterile prior to use or particulatefree before analysis – field-proven, high-quality Sartorius filter syringes are the No. 1 choice for reliable, convenient preparation steps.

Our Product Range

Minisart® syringe filters feature a housing made of medical-grade acrylic (MBS), which makes them the perfect choice for sterile filtration and clarification of additives, buffers, reagents, liquids and gases. Their effective filtration area of 6.2 cm² for liquids is the largest among premium syringe filters. The combination of the MBS housing material and surfactant-free cellulose acetate, PES or PTFE membranes, provides the lowest non-specific binding. As a result, Minisart® delivers excellent performance in filtering liquids or gases for sterilization or particle removal prior to use.

Minisart® syringe filters with a polypropylene housing are optimized for filtration prior to analysis and withstand even harsh solvents and chemicals. Due to the typical range of volumes from less than 1 mL to 100 mL, these filters are available in three different diameters with an effective filtration area of 0.07 cm², 1.7 cm² and 4.8 cm² for the optimum filter capacity and the lowest hold-up volume.

Sartorius has developed a new, easy-to-use and straightforward filtration setup. The manually operated Claristep® Filtration System consisting of a station and filter units offers a novel way for clarifying your samples prior to analysis.

Claristep® Filter units are processed without syringe and are made of the purest materials. Another major benefit is that the contact time of the samples with the filters and the caps is extremely short, ensuring optimal, contamination-free results. The Claristep® Station consists of a base, a lid and an exchangeable tray for easy and accurate positioning of sample vials and Claristep® Filter units.

Claristep® syringeless filter units with RC membranes are optimized for solvents and aqueous solutions. They provide maximum chemical compatibility and exceptionally low non-specific binding of analytes.

Sartolab® vacuum filtration devices with 0.1 μ m and 0.22 μ m PES membranes for convenient filtration of 150 mL up to 1 L are ready to use and sterile. Sartolab® RF is a complete system that includes a receiver flask. Sartolab® BT Sterile is a bottle top filter without a receiver flasks. This enables customers to use a receiver bottle of their choice and to even expand filtration capacity, depending on the particle load of the filtered liquid by filling more than one receiver flask. Sartolab® 150V is a disposable vacuum filter with a pleated 0.22 μ m PES membrane, which is suitable for up to 15 L of liquid.

Sartolab® P20 pressure filtration devices with a 0.2 μ m SFCA or 0.2 μ m and 0.45 μ m PES membrane are available with or without a glass fiber prefilter, depending on your needs. Sartolab P20 is designed for up to 3 L volumes and can also be used in-line. The polycarbonate housing and membrane components are ideal for filtering liquids. The glass fiber prefilter types are ideal for filtering environmental samples that have a high particle load prior to analyzing such samples.

Typical Applications for Filtration Devices

- Sterile filtration of liquids and gases with virtually no effect on the ingredients
- Particle removal from liquids and gases prior to downstream processes
- Venting of vials, bottles, containers, bags and bioreactors and fermenters
- Removal of precipitates and coagulates from solutions prior to use

Minisart® Syringe Filters

Removal of Particles and Microorganisms from Liquids and Gases

Sartorius offers Minisart® syringe filters for a wide range of applications. The filters are clean and safe as they are virtually free of leachables and extractables. In addition, they reliably remove particles and microorganisms – without any leakage. Minisart® with a polypropylene housing is optimized for filtration prior to analysis and withstands even harsh solvents and chemicals. Minisart® with a housing made of medical-grade acrylic (MBS) is the perfect choice for sterile filtration and clarification of additives, buffers, reagents, liquids and gases.

A Full Range of Filters Dedicated for Various Filtration Applications



Sample Preparation HPLC | UHPLC | Analytics

Elimination of particles from your samples prior to HPLC or other chromatographic analysis is essential in order to maintain the integrity of your chromatography column and to maximize its operating lifetime. Minisart® syringe filters optimized for sample preparation consist of a

polypropylene housing and membrane components featuring maximum chemical compatibility and minimum extractables to ensure excellent results. Due to the typical range of volumes from less than 1 mL to 100 mL, these filters are available in three different diameters with an effective filtration area of 0.07 cm², 1.7 cm² and 4.8 cm². For a selection quide, please see page 65.



Filtration of Aqueous Liquids Clarification | Sterile Filtration

For clarification and sterilization of liquids, filtration is the optimal method. It removes all microorganisms and particles reliably, without any effects on the ingredients due to adsorption or decomposition. For optimal results, Minisart® syringe filters with an MBS housing provide a choice of

membranes with pore sizes ranging from 0.1 μ m to 5 μ m for high flow rates and the lowest adsorption characteristics. The effective filtration area of 6.2 cm² for the fastest filtration is the largest among premium syringe filters available, and the MBS housing is color-coded for easy pore size identification. For a list of the types offered, please see page 69.



Sterile Filtration and Sterile Venting Applications

Minisart® syringe filters are ideal for clarification of liquids laden with particles. For sterilization and removal of particles from air and other gases, syringe filters are optimal for

sterile venting of containers, bioreactors,

fermenters and tubing systems in medical devices. Minisart® syringe filters are available with a wide choice of membranes, connectors and

housing materials. For a list of the types offered, please see page 72.

Minisart® Selection Guide

Please refer to Minisart® RC, NY or SRP for the highest chemical compatibility on page 67. Please refer to Minisart® NML or Minisart® High Flow on page 70. Several additional Minisart® units for venting and special purposes are shown on page 74.

	Aqu	eous	Aqueous Solvents				
Composition	All Aqueous Solutions Buffers, Protein Analysis	All Aqueous Solutions Tissue Culture Media	Aqueous Solvent Mixtures Solvents	Solvent Mixtures Solvents	Solvents Gases Acids Bases		
Sample C	SFCA Surfactant-Free Cellulose Acetate	PES Polyethersulfone	RC Regenerated Cellulose	NY Polyamide, Nylon	PTFE Polytetrafluoroethylene		
		Hydrophobic					





Sample Preparation for Analytical Procedures

Reliable Removal of Particles and Microorganisms from Liquids and Gases


Particle removal by filtration before analysis substantially increases the lifetime of your columns. Minisart® RC is optimized for aqueous liquids and solvents and is compatible with DMSO, other amides, ketones, esters and ethers. Minisart® NY is exceptionally pure compared with other common polyamide (nylon) filters and competitor products. For this product raw materials are used which do not interfere with standard analytical methods.

Our coating-free hydrophobic PTFE membrane used in Minisart® SRP is suitable for venting applications and leachables-free clarification of exceptionally harsh chemicals.

Minisart® Features

- Low adsorption of analytes
- Maximum chemical compatibility
- Minimum extractables or leachables
- Bidirectional use possible
- 100% optical integrity test

HPLC Certification

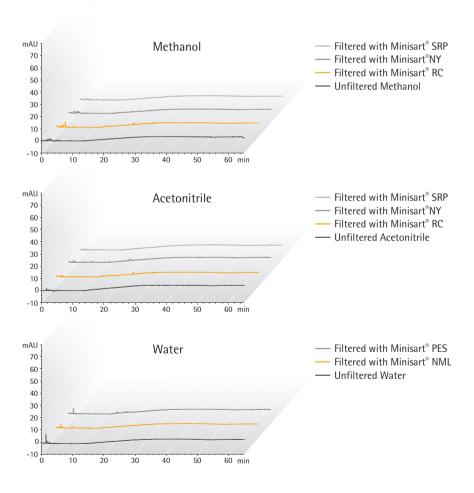

HPLC Procedure

Column C18: 250×4.6 mm, Flow Rate: 1 mL/min, Wavelength: 220 nm **HPLC** Injection Volume: 20 μ L, Analysis Time: 65 min, Temperature: 40°C, Mobile Phases: A) Acetonitrile | B) Water, Gradient: Hold 60% A for 10 min, 60% to 95% A in 20 min, 95% to 100% A in 35 min

Sample Preparation for Chromatography

Ordering Information

Dia. in mm | EFA¹ Membrane Housing Pore Size Connector Outlet Color | Printing Sterile* Oty./Pkg. Order No.

25 mm	RC	PP	0.2 μm	Male Luer Slip	White, Printed	Yes	50	17764ACK
25 mm	RC	PP	0.2 μm	Male Luer Slip	White, Printed	No	50	17764K
25 mm	RC	PP	0.2 μm	Male Luer Slip	White, Printed	No	200	17764S
25 mm	RC	PP	0.2 μm	Male Luer Slip	White, Printed	No	500	177640
25 mm	RC	PP	0.45 μm	Male Luer Slip	White, Printed	No	50	17765K
25 mm	RC	PP	0.45 μm	Male Luer Slip	White, Printed	No	200	17765S
25 mm	RC	PP	0.45 μm	Male Luer Slip	White, Printed	No	500	17765Q
15 mm	RC	PP	0.2 μm	Male Luer Slip	White, Printed	Yes	50	17761ACK
15 mm	RC	PP	0.2 μm	Male Luer Slip	White, Printed	No	50	17761K
15 mm	RC	PP	0.2 μm	Male Luer Slip	White, Printed	No	500	17761Q
15 mm	RC	PP	0.45 μm	Male Luer Slip	White, Printed	No	50	17762K
15 mm	RC	PP	0.45 μm	Male Luer Slip	White, Printed	No	500	177620
4 mm	RC	PP	0.2 μm	Male Luer Slip	Blue Tray	No	50	17821K
4 mm	RC	PP	0.2 μm	Male Luer Slip	Blue Tray	No	500	178210
4 mm	RC	PP	0.45 μm	Male Luer Slip	Yellow Tray	No	50	17822K
4 mm	RC	PP	0.45 μm	Male Luer Slip	Yellow Tray	No	500	17822Q
Minisart® SI	RP (Hydrophol	hic PTFF)						
25 mm	PTFE	PP	0.2 μm	Male Luer Slip	White, Printed	Yes	50	S7575FXOSK
25 mm	PTFE	PP	0.2 μm	Male Luer Slip	White, Printed	No	50	17575K
	PTFE PTFE	PP PP	0.2 μm 0.2 μm	Male Luer Slip Male Luer Slip	White, Printed White, Printed	No No	50 200	17575K 17575S
25 mm			<u> </u>	· · · · · · · · · · · · · · · · · · ·				
25 mm 25 mm	PTFE	PP	0.2 μm	Male Luer Slip	White, Printed	No	200	17575S
25 mm 25 mm 25 mm	PTFE PTFE	PP PP	0.2 μm 0.2 μm	Male Luer Slip Male Luer Slip	White, Printed White, Printed	No No	200 500	17575S 17575Q
25 mm 25 mm 25 mm 25 mm 25 mm 25 mm	PTFE PTFE PTFE	PP PP PP	0.2 μm 0.2 μm 0.45 μm	Male Luer Slip Male Luer Slip Male Luer Slip	White, Printed White, Printed White, Printed	No No	200 500 50	17575S 17575Q 17576K
25 mm 25 mm 25 mm 25 mm	PTFE PTFE PTFE PTFE	PP PP PP	0.2 μm 0.2 μm 0.45 μm 0.45 μm	Male Luer Slip Male Luer Slip Male Luer Slip Male Luer Slip	White, Printed White, Printed White, Printed White, Printed	No No No	200 500 50 200	17575S 17575
25 mm 25 mm 25 mm 25 mm 25 mm 15 mm	PTFE PTFE PTFE PTFE	PP PP PP PP	0.2 μm 0.2 μm 0.45 μm 0.45 μm 0.45 μm	Male Luer Slip	White, Printed White, Printed White, Printed White, Printed White, Printed	No No No No	200 500 50 200 500	17575S 17575
25 mm 25 mm 25 mm 25 mm 25 mm 15 mm	PTFE PTFE PTFE PTFE PTFE	PP PP PP PP PP	0.2 μm 0.2 μm 0.45 μm 0.45 μm 0.45 μm 0.2 μm	Male Luer Slip Male Spike	White, Printed White, Printed White, Printed White, Printed White, Printed White, Printed	No No No No No	200 500 50 200 500 50	17575S 17575
25 mm 25 mm 25 mm 25 mm 25 mm 15 mm 15 mm	PTFE PTFE PTFE PTFE PTFE PTFE	PP PP PP PP PP PP	0.2 μm 0.2 μm 0.45 μm 0.45 μm 0.45 μm 0.2 μm	Male Luer Slip Male Spike Male Spike	White, Printed	No No No No No No	200 500 50 200 500 50 50	17575S 17575
25 mm 25 mm 25 mm 25 mm 25 mm 15 mm 15 mm 15 mm	PTFE PTFE PTFE PTFE PTFE PTFE PTFE	PP PP PP PP PP PP PP	0.2 μm 0.2 μm 0.45 μm 0.45 μm 0.45 μm 0.2 μm 0.2 μm 0.2 μm	Male Luer Slip Male Spike Male Spike Male Spike	White, Printed	No No No No No No No Vo	200 500 50 200 500 50 50 50	17575S 17575
25 mm 25 mm 25 mm 25 mm 25 mm 15 mm 15 mm 15 mm 15 mm	PTFE PTFE PTFE PTFE PTFE PTFE PTFE PTFE	PP PP PP PP PP PP PP PP	0.2 μm 0.2 μm 0.45 μm 0.45 μm 0.45 μm 0.2 μm 0.2 μm 0.2 μm	Male Luer Slip Male Spike Male Spike Male Luer Slip Male Luer Slip Male Luer Slip	White, Printed	No No No No No No No No No	200 500 50 200 500 50 50 50 50	17575S 17575
25 mm 25 mm 25 mm 25 mm 25 mm 15 mm 15 mm 15 mm 15 mm 15 mm	PTFE PTFE PTFE PTFE PTFE PTFE PTFE PTFE	PP PP PP PP PP PP PP PP PP	0.2 μm 0.2 μm 0.45 μm 0.45 μm 0.45 μm 0.2 μm 0.2 μm 0.2 μm 0.2 μm 0.2 μm	Male Luer Slip Male Spike Male Spike Male Spike Male Luer Slip Male Luer Slip Male Luer Slip	White, Printed	No No No No No No No No No No No	200 500 50 200 500 50 50 50 50 50	17575S 17575
25 mm 25 mm 25 mm 25 mm 25 mm 15 mm 15 mm 15 mm 15 mm 15 mm 15 mm	PTFE PTFE PTFE PTFE PTFE PTFE PTFE PTFE	PP	0.2 μm 0.2 μm 0.45 μm 0.45 μm 0.45 μm 0.2 μm 0.2 μm 0.2 μm 0.2 μm 0.2 μm 0.45 μm	Male Luer Slip Male Spike Male Spike Male Spike Male Luer Slip Male Luer Slip Male Luer Slip Male Luer Slip Male Spike	White, Printed	No Yes No No No	200 500 50 200 500 50 50 50 50 50 50	17575S 17575
25 mm 25 mm 25 mm 25 mm 25 mm 25 mm 15 mm	PTFE PTFE PTFE PTFE PTFE PTFE PTFE PTFE	PP	0.2 μm 0.2 μm 0.45 μm 0.45 μm 0.45 μm 0.2 μm 0.2 μm 0.2 μm 0.2 μm 0.2 μm 0.45 μm	Male Luer Slip Male Spike Male Spike Male Spike Male Luer Slip Male Luer Slip Male Luer Slip Male Luer Slip Male Spike Male Spike	White, Printed	No N	200 500 50 200 500 50 50 50 50 50 50 50	17575S 17575
25 mm 25 mm 25 mm 25 mm 25 mm 15 mm	PTFE PTFE PTFE PTFE PTFE PTFE PTFE PTFE	PP	0.2 μm 0.2 μm 0.45 μm 0.45 μm 0.45 μm 0.2 μm 0.2 μm 0.2 μm 0.2 μm 0.2 μm 0.45 μm 0.45 μm	Male Luer Slip Male Spike Male Spike Male Luer Slip Male Spike Male Spike Male Spike	White, Printed	No N	200 500 50 200 500 50 50 50 50 50 50 50 50	17575S 17575S 17576
25 mm 25 mm 25 mm 25 mm 25 mm	PTFE PTFE PTFE PTFE PTFE PTFE PTFE PTFE	PP	0.2 μm 0.2 μm 0.45 μm 0.45 μm 0.45 μm 0.2 μm 0.2 μm 0.2 μm 0.2 μm 0.2 μm 0.45 μm 0.45 μm 0.45 μm	Male Luer Slip Male Spike Male Spike Male Luer Slip Male Spike Male Spike Male Spike Male Spike Male Spike	White, Printed	No N	200 500 50 200 500 50 50 50 50 50 50 50 50	17575S 17575S 17576

Dia. in mm | EFA¹ Membrane Housing Pore Size Connector Outlet Color | Printing Sterile* Oty./Pkg. Order No.

Minisart[®] NY (Nylon) and NY25 Plus (Glass Fiber 0.7 μm² + Nylon)

25 mm	Nylon	PP	0.2 μm	Male Luer Slip	White, Printed	Yes	50	17845ACK
25 mm	Nylon	PP	0.2 μm	Male Luer Slip	White, Printed	No	500	17845Q
25 mm	Nylon	PP	0.45 μm	Male Luer Slip	White, Printed	Yes	50	17846ACK
25 mm	Nylon	PP	0.45 μm	Male Luer Slip	White, Printed	No	500	17846Q
15 mm	Nylon	PP	0.2 μm	Male Luer Slip	White, Printed	No	50	1776BK
15 mm	Nylon	PP	0.2 μm	Male Luer Slip	White, Printed	No	500	1776BQ
15 mm	Nylon	PP	0.45 μm	Male Luer Slip	White, Printed	No	50	1776CK
15 mm	Nylon	PP	0.45 μm	Male Luer Slip	White, Printed	No	500	1776CQ
25 mm	GF+Nylon	PP	0.2 μm	Male Luer Slip	White, Printed	No	50	1784BK
25 mm	GF+Nylon	PP	0.2 μm	Male Luer Slip	White, Printed	No	500	1784BQ
25 mm	GF+Nylon	PP	0.45 μm	Male Luer Slip	White, Printed	No	50	1784CK
25 mm	GF+Nylon	PP	0.45 μm	Male Luer Slip	White, Printed	No	500	1784CQ

^{*} Sterile Minisart® syringe filters are individually packaged. If not stated otherwise, Minisart® units have been sterilized by ethylene oxide. Non-presterilized Minisart® units: RC, PTFE and nylon can be sterilized by autoclaving at 121°C for 30 min/or by using ethylene oxide (EO).

For technical product specifications, please see page 76.

¹ Diameter of EFA – Effective Filtration Area

 $^{^{2}}$ 0.7 µm = GF particle retention ≠ pore size!



Minisart® High Flow with PES



Minisart® NML with SFCA



Minisart® HY with PTFE

Clarification and Sterilization by Filtration

Filtration Is the Optimal Method for Clarification and Sterilization of Liquids and Gases


Sterilization by filtration is the fastest method for removal of bacterial cells from liquids, while minimizing the effects on ingredients. Minisart® NML with surfactant-free cellulose acetate (SFCA) is the best choice for all aqueous solutions with a pH of 4 to 8. It combines fast flow rates with an exceptionally pure and leachables-free membrane and is available in many different pore sizes - also for the removal of larger particles. Minisart® High Flow with polyethersulfone (PES) is optimal for delivering the highest flow rates and for a broad pH compatibility range from 1 to 13. Due to the asymmetric membrane structure, the PES surface virtually behaves like a prefilter. Both Minisart® types -NML and High Flow - can be sterilized by

ethylene oxide (EO) or gamma irradiation and are suitable for gluing and assembling on tubing. Minisart HY with hydrophobic PTFE membrane is the perfect choice for oily solutions or gases.

Minisart® Features

- The largest effective filtration area (EFA) of 6.2 cm² for NML and High Flow
- The lowest adsorption
- Superior flow rate
- High total throughput
- Low hold-up volume
- Minimum extractables
- PVC-free
- Gamma-irradiated or EO-sterilized
- Bidirectional use possible
- 100% optical integrity test


Water Flow Rates at 1 bar (mL/min); 15.4 psi; 0.45 µm Hydrophilic Membranes

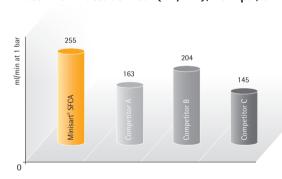

0.45 µm SFCA membrane compared with CA (Cellulose Acetate) or similar hydrophilic membranes of the same pore size



0.45 µm PES membrane compared with other PES membranes of the same pore size

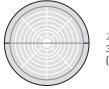

28 mm EFA 33 mm housing diameter (for NML and High Flow)

Preparation of Aqueous Liquids

Ordering Information

Dia. in mm EFA ¹	Membrane	Housing	Pore Size	Connector Outlet	Color Printing	Sterile*	Qty./ Pkg.	Order No.
Minisart [®]	High Flow (PES – Poly	ethersulfo	one)				
28 mm	PES	MBS	0.1 μm	Male Luer Lock	Dark Red	Yes	50	16553K
28 mm	PES	MBS	0.22 μm	Male Luer Lock	Royal Blue	Yes#	50	16532GUK
28 mm	PES	MBS	0.22 μm	Male Luer Lock	Royal Blue	Yes	50	16532K
28 mm	PES	MBS	0.22 μm	Male Luer Slip	Royal Blue	Yes	50	16541K
28 mm	PES	MBS	0.22 μm	Male Luer Lock	Royal Blue	No	500	16532Q
28 mm	PES	MBS	0.22 μm	Male Luer Slip	Royal Blue	No	500	16541Q
28 mm	PES	MBS	0.45 μm	Male Luer Lock	Amber	Yes	50	16537K
28 mm	PES	MBS	0.45 μm	Male Luer Lock	Amber	No	500	16537Q
28 mm	PES	MBS	0.45 μm	Male Luer Slip	Amber	Yes#	50	16533GUK
28 mm	PES	MBS	0.45 μm	Male Luer Slip	Amber	Yes	50	16533K
28 mm	PES	MBS	0.45 μm	Male Luer Slip	Amber	No	500	16533Q
28 mm	SFCA	MBS	0.2 μm	Male Luer Lock	Blue	Yes Yes#	50	S6534FMOSK
28 mm	SFCA	MBS	0.2 μm	Male Luer Lock	Blue	Yes#	50	S6534FMGU
28 mm	SFCA	MBS	0.2 μm	Male Luer Lock	Blue	No	500	S6534FXQ
28 mm	SFCA	MBS	0.2 μm	Male Luer Slip	Blue	Yes	50	S7597FXOSK
28 mm	SFCA	MBS	0.2 μm	Male Luer Slip	Blue	No	500	S7597FXQ
28 mm	SFCA	MBS	0.45 μm	Male Luer Lock	Yellow	Yes	50	S6555FMOSk
28 mm	SFCA	MBS	0.45 μm	Male Luer Lock	Yellow	Yes#	50	S6555FMGU
28 mm	SFCA	MBS	0.45 μm	Male Luer Lock	Yellow	No	500	S6555FMQ
28 mm	SFCA	MBS	0.45 μm	Male Luer Slip	Yellow	Yes	50	S7598FXOSK
28 mm	SFCA	MBS	0.45 μm	Male Luer Slip	Yellow	No	500	S7598FXQ
28 mm	SFCA	MBS	0.65 μm	Male Luer Slip	Pink	Yes	50	16569K
28 mm	SFCA	MBS	0.8 µm	Male Luer Lock	Green	Yes	50	16592K
28 mm	SFCA	MBS	0.8 µm	Male Luer Lock	Green	Yes#	50	16592GUK
28 mm	SFCA	MBS	0.8 µm	Male Luer Lock	Green	No	500	16592Q
28 mm	SFCA	MBS	1.2 μm	Male Luer Lock	Red	Yes	50	17593K
28 mm	SFCA	MBS	1.2 μm	Male Luer Lock	Red	No	500	17593Q
28 mm	SFCA	MBS	5 μm	Male Luer Lock	Brown	Yes	50	S7594FMOSK
28 mm	SFCA	MBS	5 μm	Male Luer Lock	Brown	No	500	17594Q

Dia. in mm EFA ¹	Membrane	Housing	Pore Size	Connector Outlet	Color Printing	Sterile*	Qty./ Pkg.	Order No.		
Minisart [®]	NML Plus (Glass Fiber	0.7 μm ² +	SFCA)						
28 mm	GF+SFCA	MBS	0.2 μm	Male Luer Lock	Blue	Yes	50	17823K		
28 mm	GF+SFCA	MBS	0.2 μm	Male Luer Lock	Blue	No	500	17823Q		
28 mm	GF+SFCA	MBS	0.45 μm	Male Luer Lock	Yellow	Yes	50	17829K		
28 mm	GF+SFCA	MBS	0.45 μm	Male Luer Lock	Yellow	No	500	17829Q		
28 mm	GF+SFCA	MBS	1.2 μm	Male Luer Lock	Red	No	500	17825Q		
28 mm	GF	MBS	$0.7~\mu m^2$	Male Luer Lock	White	No	50	17824K		
28 mm	GF	MBS	0.7 μm ²	Male Luer Lock	White	No	500	17824Q		
Minisart®	Minisart® HY (hydrophobic PTFE), for Venting and Gas Filtration									
26 mm	PTFE	MBS	0.2 μm	Male Luer Lock	Clear	Yes	50	S6596FMOSK		
26 mm	PTFE	MBS	0.2 μm	Male Luer Lock	Clear	No	500	S6596FMQ		

^{*} Sterilized Minisart® units are individually packaged. If not stated otherwise, Minisart® are sterilized by ethylene oxide. #-Mark indicates sterilization by gamma irradiation.

For technical product specifications, please see page 78.



Non-presterilized Minisart® units: High Flow, NML, NML Plus and HY can be sterilized by ethylene oxide; High Flow, NML and NML Plus can also be sterilized by gamma irradiation

¹ Diameter of EFA – Effective Filtration Area

 $^{^{2}}$ 0.7 µm = GF particle retention ≠ pore size!



Minisart® available in many configurations

Special Applications – Sterile Filtration and Sterile Venting

Choose from a Broad Range of Pore Sizes, Materials and Formats


Bacterial cell removal or particulate removal from liquids can easily be performed with Minisart® syringe filters. Minisart® has a minimal effect on the ingredients of a filtered solution. Minisart® NML and Ophthalsart with surfactant-free cellulose acetate (SFCA) and Minisart® HY and SRP with hydrophobic PTFE are frequently used for sterile filtration of aqueous and oily solutions and other liquids. Minisart® NML with a 5 μm pore size removes particulates or coagulates, offering the highest total throughput under sterile conditions without any early clogging. Hydrophobic PTFE filters are suitable for venting purposes and are additionally available in special

formats with activated carbon.

Minisart® Features

- 100% optical integrity test
- Low adsorption
- Minimum extractables
- Particulate-free
- PVC-free
- Gamma-irradiated or EO-sterilized
- Bidirectional use possible
- Many configurations available

Pressure-Hold Test

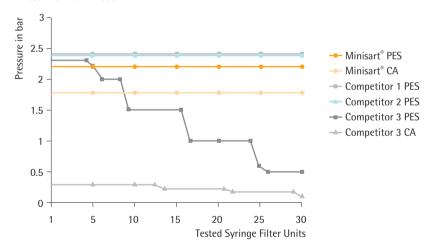

Method: Pressure-hold tests were performed by connecting pre-wetted syringe filter units to a 10-fold pressure device with a pressure gauge. The pressure applied equaled 4/5 of the bubble point. Units failed the test if they released air bubbles before reaching 4/5 of the minimum bubble point. These units were remeasured to determine at which pressure air bubbles appear.

Result: Testing the pressure-hold capability of syringe filters revealed that many competitor 3's filters are dysfunctional and not intact. Filtration performed with such filters will result in unsterile filtrates or filtrates without an appropriately reduced level of particles.



Biocompatibility Certificate Minisart® HY



Biocompatibility Certificate Minisart® NML



ISO 13485:2016

Minisart® Syringe Filters – Special

Ordering Information

Dia. in mm EFA ¹	Membrane	Housing	Pore Size	Connector Outlet	Color Printing	Sterile*	Oty./ Pkg.	Order No.
Minisart®	NML (SFCA	– Cellulos	e Acetate)	Aqueous Filtration	1			
28 mm	SFCA	MBS	0.2 μm	Male Luer Lock	Blue	Yes	50	16534K*
28 mm	SFCA	MBS	0.2 μm	Male Luer Lock	Blue	Yes#	50	16534GUK*
28 mm	SFCA	MBS	0.2 μm	Male Luer Lock	Blue	No	500	16534Q*
28 mm	SFCA	MBS	0.2 μm	Male Luer Slip	Blue	Yes	50	17597K *
28 mm	SFCA	MBS	0.2 μm	Male Luer Slip	Blue	No	500	17597Q*
28 mm	SFCA	MBS	0.45 μm	Male Luer Lock	Yellow	Yes	50	16555K*
28 mm	SFCA	MBS	0.45 μm	Male Luer Lock	Yellow	Yes#	50	16555GUK*
28 mm	SFCA	MBS	0.45 μm	Male Luer Lock	Yellow	No	500	16555Q*
28 mm	SFCA	MBS	0.45 μm	Male Luer Slip	Yellow	Yes	50	17598K*
28 mm	SFCA	MBS	0.45 μm	Male Luer Slip	Yellow	No	500	17598Q*
28 mm	SFCA	MBS	5 μm	Male Luer Lock	Brown	Yes	50	17594K*
Minisart [®]	Ophthalsar	t (SFCA –	Cellulose A	cetate) Aqueous F	iltration			
28 mm	Ophthalsart	•	0.2 μm	Male Luer Slip	Pink	Yes	50	17528K*
Minisart [®]	High Flow (PES – Poly	vethersulfo	ne) Aqueous Filtra	ntion			
28 mm	PES	MBS	0.1 μm	Male Luer Lock	Dark Red	Yes	50	16553K
	(
	PES (Polyeti				14/1 14			47700 401/
15 mm	PES	PP	0.22 μm	Male Luer Slip	White	Yes	50	1776DACK
15 mm	PES	PP	0.22 μm	Male Luer Slip	White	No	50	1776DQ
Minisart®	PES (hydro	phobic PES	S), for Ven	ting and Gas Filt	ration			
25 mm	PES-	PP	0.2 μm	Male Luer Slip	White	No	50	1757HK
25 mm	PES-	PP	0.2 μm	Male Luer Slip	White	No	500	1757HQ
25 mm	PES-	PP	0.2 μm	Hose Barb ^b	White	No	50	1757GK
25 mm	PES-	PP	0.2 μm	Hose Barb ^b	White	No	500	1757GQ
Minisort®	Air (Hydrop	hobio DTEI	E) Vantina					
15 mm	PTFE	MBS	0.2 µm	Male Luer Slip	Yellow	No	500	1751AQ
15 mm	PTFE	MBS	0.2 μm	Male Luer Slip +		Yes#	500	16596HNK
15 111111	FIFE	IVIDS	υ.2 μπ	Needle	Tellow	105#	50	10390ПІЛК
Minisart®	HY (hydron	nobic PTFF	for Vent	ing and Gas Filtrat	ion			
26 mm	PTFE	MBS	0.2 μm	Male Luer Lock	Clear	Yes	50	16596HYK*
26 mm	PTFE	MBS	0.2 μm	Male Luer Lock	Clear	No	500	16596HYQ*
26 mm	PTFE	MBS	0.2 μm	Male Luer Lock ^a	Clear	No	500	16599HYQ*
26 mm	PTFE	MBS	1 μm	Male Luer Lock	Clear	No	500	1659AHYQ
20 111111	11111	כטואי	ιμιιι	wait Luci Luck	Cicai	INO	300	1000/11110

Minisart	* Acticosart	with Dome	e Reservoir	and Hydrophobic	PTFE for Ver	iting and	Ultracle	aning of Gases
26 mm	Active Carbon	MBS	0.45 μm	Male Luer Slip ^a	Blue	No	500	17840Q
Dia. in mm EFA	Membrane 1	Housing	Pore Size	Connector Outlet	Color Printing	Sterile**	Qty./ Pkg.	Order No.
Minisart	SRP (Hydro	phobic PTF	E) Venting	& Gas Filtration				
25 mm	PTFE	PP	0.2 μm	Male Luer Slip	White, Printed	Yes	50	17575ACK
25 mm	PTFE	PP	0.2 μm	Hose Barb	White	No	500	1757AQ

For technical product specifications, please see pages 76 and 78.

Do you need other quantities per package? Are you looking for special versions or do you have requirements? Do you need other inlet and or outlet connectors?

Please contact us to learn more about additional Minisart® configurations that are available.



^{*} Article numbers are not available in: Argentina, Bolivia, Chile, Paraguay, Peru, Uruguay, USA
** Sterilized Minisart® units are individually packaged. If not stated otherwise, Minisarts are sterilized by ethylene oxide. #-mark indicates sterilization by gamma irradiation

Non-presterilized Minisart* units: SFCA can be sterilized by ethylene oxide or gamma irradiation. PTFE can be sterilized by ethylene oxide.

^a Connector inlet: Male Luer slip (all other Minisart® types have female luer lock inlets)

b Hose barbs, inlet and outlet, 4 - 6 mm stepped

¹ Diameter of EFA – Effective Filtration Area

Minisart® with Polypropylene Housing

()	_	, C.	4.5	
	>ne	citic	OITC	nc
	-		.a uv	HJ

Cytotoxity (17575-ACK)

Minisart® RC SRP NY PES with	4 15 25 mm Dia	ımeter Membrar	e Filtration Area	a	
Housing material	Polypropylene (
Membranes		ed Cellulose NY obic PTFE = Polyte		e PES = Polyethers	sulfone
Glass fiber prefilter	NY Plus: Ultrapı	ıre quartz, 0.7 μm	particle retentio	n	
Application limits	Max. recommer	nded operating pr	essure 4.5 bar 65	i psi	
Housing burst pressure	>7 bar 102 psi				
Max. temperature	121°C, 30 min (autoclavable)			
Sterilization	Non-sterile Min	isart [®] can be auto	oclaved or sterilize	ed by ethylene oxid	le (EO)
Minisart® Membrane Types	RC 0.2 μm	RC 0.2 μm	RC 0.45 μm	SRP 0.2 μm	SRP 0.45 μm
Non-sterile packs: 50 (K), 200 (S), 500 (Q), 1000 (R) sterile packs: individual packaged, 50 (ACK)	K S Q R	ACK	K S Q R	K S Q ACK	K S Q
Bubble point (≥)	With water 3.0 bar 44 psi	With water 4.6 bar 67 psi	With water 2.0 bar 29 psi	With ethanol 1.1 bar 16 psi	With ethanol 0.9 bar 13 psi
Flow rate, 4 mm dia. = 0.07 cm ² fi - For water at 1 bar	lter area Hold-up 0.5	volume¹: 5 μL to –	0 10 μL (▶ mL/mir 1.5	n) _3	_3
- For methanol at 1 bar	1.5	_	3.0	2.0	4.5
- For air at 0.1 bar	_2	-	_2	30	60
Flow rate, 15 mm dia. = 1.7 cm ² fi	lter area Hold-up	o volume¹: 30 μL t	to 100 μL (▶ mL/r	nin)	
– For water at 1 bar	20	10	40	_3	_3
- For methanol at 1 bar	55	25	105	55	150
- For air at 0.1 bar	_2	_2	_2	800	1600
Flow rate, 25 mm dia. = 4.8 cm ² fi	Iter area Hold-up	o volume¹: 100 μl	. to 200 μL (▶ mL		
– For water at 1 bar	80	50	160	_3	_3
- For methanol at 1 bar	160	90	325	160	260
- For air at 0.1 bar	_2	_2	_2	1800	3000
Water penetration point	_	_	_	>4.0 bar 58 psi ³	> 3.0 bar 44 psi
Sterile filtration capability acc. to the bacteria challenge text	No ⁵	Yes	No	Yes	No
Pyrogen-free according				Yes	

No inhibition with MRC-5 (human lung cells) and L929

Minisart® Membrane Types	NY 0.2 μm	NY 0.45 μm	NY Plus 0.2 µm	NY Plus 0.45 µm	PES 0.2 μm	PES -0.2 μm
Non-sterile packs: 50 (K), 200 (S), 500 (Q), 1000 (R) sterile packs: individual packaged, 50 (ACK)	K Q R ACK	K Q R ACK	K Q	K Q	K Q ACK	K Q
Bubble point (≥)	With water	With water	With water	With water	With water	With ethanol
	3.0 bar 44 psi	2.0 bar 29 psi	3.0 bar 44 psi	2.0 bar 29 psi	3.2 bar 46 psi	0.95 bar 14 ps
Flow rate, 4 mm dia. = 0.07	cm² filter area	Hold-up volum	e¹: 5 μL to 10 μL (l	► mL/min)		
– For water at 1 bar	-	-	-	-	1.5	-
– For methanol at 1 bar	-	_	_	-	_4	_
- For air at 0.1 bar	_	_	_	-	_2	_
Flow rate, 15 mm dia. = 1.7 - For water at 1 bar	20	40	-		40	-
For water at 1 bar	20	40	-	-	40	-
- For methanol at 1 bar	40	110	-	_	_4	_
- For air at 0.1 bar	_2	_2	-	-	_2	-
Flow rate, 25 mm dia. = 4.8	cm² filter area	Hold-up volum	e¹: 100 μL to 200	μL (▶ mL/min)		
- For water at 1 bar	50	100	50	100	100	-
– For methanol at 1 bar	70	200	70	200	_4	_4
- For air at 0.1 bar	_2	_2	_2	_2	_2	1200
Water penetration point	_	_	_	_	_	>3 bar 44 psi ³
Sterile filtration capability acc. to the bacteria challenge test	Yes	No	Yes	No	Yes	Yes
Pyrogen-free according to the USP						-

Cytotoxity (17575-ACK) No inhibition with MRC-5 (human lung cells) and L929

² Hydrophilic membranes can filter dry air or gas but become impermeable to air or gas when wetted!

¹ Hold-up volume after air purge

³ Hydrophobic membranes cannot be wetted with aqueous solutions unless you overcome their water penetration point or pre-wet them using an organic solvent (e.g. ethanol)

⁴ PES is suitable for solutions only containing up to 30% MeOH

⁵ According to the bacterial challenge test (BCT) with 10⁷ Brevundimonas diminuta. Non-sterile RC Minisart® types are optimized for sample preparation and are not suitable for sterile filtration according to the bacteria challenge test. All other non-sterile Minisart® types with 0.2 mm pore size can be sterilized by autoclaving or EO before use for sterile filtration.

Minisart® with MBS Housing

Specifications

 $\label{eq:minisart} \textbf{Minisart}^* \ \textbf{High Flow} \ | \ \textbf{NML} \ | \$

Minisart[®] HY | Acticosart with 26 mm accessible membrane filtration area diameter, 100 μL to 150 μL hold-up volume¹ Minisart[®] Air with 15 mm accessible membrane filtration area diameter, 100 μm hold-up volume¹

Housing material	Methacrylate butadiene styrene (MBS)
Membranes	High Flow: PES = Polyethersulfone
	NML: (SF) CA = (Surfactant-free) Cellulose Acetate
	NML Plus: (SF) CA = (Surfactant-free) Cellulose Acetate
	HY Acticosart Air: Hydrophobic PTFE = Polytetrafluoroethylene
Glass fiber pre-filter	NML Plus: Binder-free GF, 0.7 μm particle retention
Application limits	High Flow: Max. recommended operating pressure 6.0 bar 87 psi
	NML, NML Plus, HY, Air: Max. recommended operating pressure 4.5 bar 65 psi
	Acticosart: Max. recommended operating pressure 1 bar 14.5 psi
Housing burst pressure	>7 bar 102 psi (not determined for Acticosart)
Max. temperature	60°C, not autoclavable
Sterilization	Non-sterile Minisart® High Flow, NML and NML Plus can be or sterilized by ethylene oxide
	(EO) or by gamma sterilization
	Non-sterile Minisart® HY, Acticosart, Air* can be sterilized by ethylene oxide (EO)

Minisart® Membrane Types	PES 0.1 μm	PES 0.2 μm	PES 0.45 μm	CA 0.2 μm	CA 0.45 μm	CA 0.65 μm	CA 0.8 μm	CA 1.2 μm	CA 5.0 μm
Non-sterile packages: 500 (Q, HYQ), 1000 (R), sterile packs: individually packaged, 50 (K, GUK, HYK, HNK)	K	K GUK Q	K GUK Q	K GUK Q	K GUK Q	K	K GUK Q	K Q	K Q
Bubble point (≥)	With water 5.0 bar 73 psi	With water 3.2 bar 46 psi	With water 2.0 bar 29 psi	With water 3.2 bar 46 psi	With water 2.0 bar 29 psi	With water 1.3 bar 19 psi	With water 0.8 bar 12 psi	With water 0.7 bar 10 psi	With water 0.4 bar 6 psi
Flow Rate for ^{2 3} (▶ mL/min)									
28 mm dia. for water at 1 bar	40	140	220	60	180	250	400	500	600
15 mm dia. for air at 0.1 bar	_	_	_	-	_	_	-	-	_
26 mm dia. for air at 0.1 bar	-	-	_	-	-	-	-	_	_
Water penetration point	_	_	_	_	_	_	_	_	_
Sterile filtration capability ⁴ acc. to the bacteria challenge test	Yes	Yes	No	Yes	No	No	No	No	No
Pyrogen-free according to the USP				Yes	Yes				Yes
Cytotoxity	No inhib	tion with I	MRC-5 (hu	man lung	cells) and	L929			

Minisart® Membrane Type	GF+CA 0.2 μm	GF + CA 0.45 μm	GF+CA 1.2 μm	GF 0.7 μm	PTFE 0.2 μm	PTFE 1.0 μm	Acticosart	PTFE (Air) 0.2 μm
Non-sterile packages: 500 (Q, HYQ), 1000 (R), sterile packs: individually packaged, 50 (K, GUK, HYK, HNK)	K Q	K Q	Q	K Q	HYK HYQ	НҮО	Q	Q HNK
Bubble point (≥)	With water 3.2 bar 46 psi	With water 2.0 bar 29 psi	With water 0.7 bar 10 psi	With water 0.5 bar 7 psi	With ethanol 1.4 bar 20 psi	With ethanol 0.5 bar 7 psi	With ethanol 0.9 bar 13 psi	With ethanol 1.0 bar 14 psi
Flow rate for ^{2 3} (► mL/min)								
28 mm dia. for water at 1 bar	60	160	350	450	-	-	-	-
15 mm dia. for air at 0.1 bar	_	_	_	_	_	_	_	800
26 mm dia. for air at 0.1 bar	-	-	-	-	2000	4000	2300	-
Water penetration point	-	_	_	_	> 4.0 bar 58 psi ³	> 1.5 bar 22 psi ³	N.a.	> 3.0 bar 44 psi ³
Sterile filtration capability ⁴ according to the bacteria challenge test	Yes	No	No	No	Yes	No	N.a.	Yes
Pyrogen-free according to the USP					Yes			

Cytotoxity

No inhibition with MRC-5 (human lung cells) and L929

¹ Hold-up volume after air purge

² Hydrophilic membranes can filter dry air or gas but become impermeable to air or gas when wetted!

³ Hydrophobic membranes cannot be wetted with aqueous solutions unless you overcome their water penetration point.
⁴ According to bacterial challenge test (BCT) with 10⁷ Brevundimonas diminuta. All non-sterile Minisart[®] types listed above can be sterilized according to the method recommended in this table.

^{*} Minisart® Air can be sterilized by gamma rays according to the following parameters: Range 25 – 40 kGy (validated with 50 kGy).

Claristep® Filtration System

The Power of Simplicity



The Claristep® Station consists of a base, a lid and an exchangeable tray for easy and accurate positioning of sample vials and Claristep® Filter units.

The patent-pending design features unique grooves in the station's lid and matching guide ridges on Claristep® Filter units to enable intuitively correct alignment and convenient handling of the system.

Preparing samples by clarification is an essential step prior to nearly all analytical techniques, such as high pressure liquid chromatography (HPLC). This filtration step to eliminate particles is crucial for maintaining the integrity of chromatography columns and for maximizing their operating life time.

In addition, as the sensitivity of automated analytical instruments continues to improve, they increasingly require less volume to operate in order to maximize throughput. Therefore, fast clarification of small volumes that does not add leachables or extractables to the original sample is indispensable for achieving the best analytical results.

To meet these requirements, Sartorius has developed a new, easy-to-use and straightforward filtration setup.

The manually operated Claristep® Filtration System consisting of a station and filter units offers a novel way for clarifying your samples prior to analysis.

- Up to 8 samples are processed simultaneously
- No syringe required
- No need for a vacuum source or a power supply
- For low sample volumes ranging from 60 μL to 600 μL
- Hold-up volume < 30 μL





The grooves automatically guide the filter unit caps into the correct positions for simultaneous and accurate cap closure.



Claristep® Filter units are made of the purest materials. Another major benefit is that the contact time of the samples with the filters and the caps is extremely short, ensuring optimal, contamination–free results. Filtered liquids are collected in any 12 × 32 mm outer diameter vials of your choice based on the analytical method to be performed.

Sample Preparation for Analytics

Use the Most Ergonomic Clarification Solution

Filter 8 samples simultaneously – without needing any power supply or a vacuum | pressure source. Simply place the filters on

your vials, gently close the station and press on the station lid to filter – that's it!



 Close the station lid. The grooves align the caps automatically, securely sealing every single Claristep® Filter unit for the most convenient processing.



2. Apply slight uniform pressure with your hand to start sample clarification. You will feel a certain resistance while liquid is pressed through each membrane.



3. Press down on the station lid so that the left and right corners touch the base plate. Hold the lid in place for 3 seconds to ensure all sample liquid is filtered through.



Claristep® Filter units press liquid through each membrane by an air pocket that forms over each filter unit when the station lid is closed. This air pocket is released when you stop holding down the lid – you will feel it in your fingertips!



Before clarification, the samples are pipetted in the filter reservoire.




After clarification, the filtrates are collected in sample vials.


Reliable Removal of Particles


Filter Samples Without Adding Extractables and Leachables

Claristep® Filter units with RC membranes are optimized for solvents and aqueous solutions. They provide maximum chemical compatibility and exceptionally low non-specific binding of analytes.

HPLC Certification

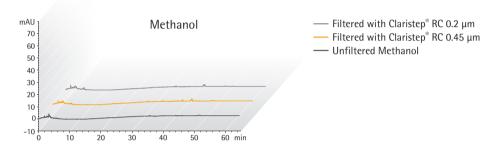

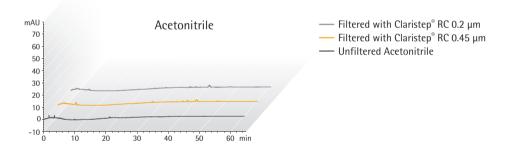

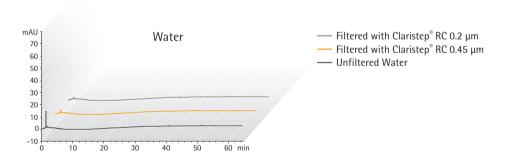

HPLC Procedure

Column: C18: $5 \mu m \times 250 \text{ mm} \times 4.0 \text{ mm}$, Flow Rate: 1 mL/min, Wavelength: 220 nm Injection Volume: 20 μ L, Analysis Time: 65 min, Temperature: $40 \,^{\circ}$ C, Mobile Phases: A) Acetonitrile | B) Water, Gradient: Hold 60% A for 10 min, 60% to 100% A in 20 min, $100 \,^{\circ}$ A for 30 min

Sample Preparation Techniques

Choose the Best Solution for Your Needs

Do you process dozens of samples each day? A syringeless solution will help you reduce time, effort and waste – and minimize hand stress. If you need to analyze only a few samples a day, you will benefit from our proven combination of a syringe and syringe filter. The choice is all yours!


Claristep® Filtration System



Syringe Filters

for processing one sample at a time



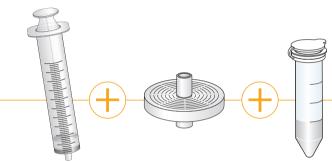

Analytical Sample Volumes Run Small

Get the Particle-free Volume You Really Need

If you need to fill only 12×32 mm vials, a syringeless solution will help you save time and reduce sample loss!


Syringeless Solution



Different Diameters

for different sample volumes



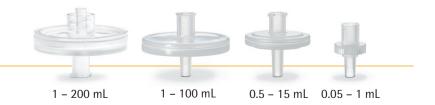



Claristep® Filters are availabe in a choice of two pore sizes

Ordering Information

Claristep® Filters

\emptyset mm EFD ¹	Membrane	Housing	Pore Size	Sterile	Qty Pk	Order No.
9.7 mm	RC	PP	0.2 μm	No	96	17C07FT96
9.7 mm	RC	PP	0.2 μm	No	480	17C07FT480
9.7 mm	RC	PP	0.45 μm	No	96	17C06FT96
9.7 mm	RC	PP	0.45 μm	No	480	17C06FT480

¹ Effective Filtration Diameter RC = Regenerated Cellulose

Claristep® System

Name	Qty Pk	Order No.
Claristep® Station complete	1	17CM8
Claristep® Single Tray	1	17CS1







The Tray can be removed and exchanged



 12×32 mm sample vials

Additional Components Needed


The free choice of 12×32 mm sample vials and lids is enabeling you to chose the right vial for your particular sample and application, e.g. for light sensitive substances you can use brown glass. For small sample volumes you can use vessels with inlays. You can use glass or plastic, screw caps and | or slid lids - whatever you prefer.



Sartolab® P20 and Sartolab® P20 Plus

Sterile Filtration of Sample Volumes of Up to 5 Liters



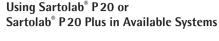

Systems with Luer Lock Connectors
Sartolab® units with a luer lock inlet fit
directly onto the corresponding connectors
of tubing from peristaltic pumps. This
configuration can also be used with luer

lock syringes.



Systems with Tubing

If the liquid to be filtered will be pumped through open tubing from a pressure tank or a peristaltic pump, then a Sartolab® unit with a hose nipple inlet is required. The stepwise increase in diameter of the nipple, from 6 mm to 12 mm, makes it suitable for a wide range of tubing.

Application

The Sartolab® P20 is a ready-to-use pressure filtration unit for sterile filtration of media and aqueous solutions in batches ranging from 100 mL to 5 L. For media that contain sera and difficult-to-filter solutions, a Sartolab® P20 Plus unit with an incorporated prefilter is also available.


Security


The effectiveness of every batch of 0.2 µm cellulose acetate membranes for sterile filtration is confirmed by bacteria challenge tests (HIMA) using Brevundimonas diminuta. Only biosafe

material is used in the filtration units; they have been proven non-toxic by passing the USP plastics test for toxicity. Tests with MRC-5 human lung cells on cellulose acetate membranes and glass fiber prefilters showed no cytotoxic effects. Finished units are also tested for their sterile filtration capability and for housing and membrane integrity.

Fast Filtration

The combination of a large filtration area (20 cm²) and the optimal design of the filter support guarantees high flow rates with high total throughputs. Automatic venting of any trapped air through the PTFE membrane-protected vent ports ensures that the entire filter surface is used for effective filtration.

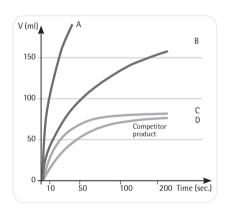

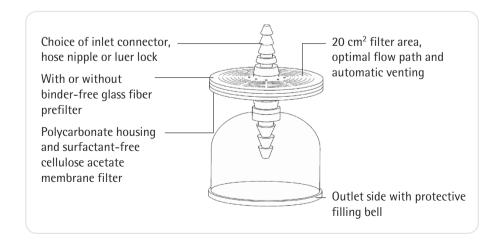

Specifications

Technical Specifications for Sartolab® P20 and Sartolab® P20 Plus Units with SFCA Membrane and PTFE Venting

Properties	Description	
	Sartolab® P20 18052 18053	Sartolab® P20 Plus 18056 18058
Filter material	SFCA, type 12587, 0.2 μm pore size and PTFE	SFCA, type 12587 0.2 µm pore size plus GF, 100% free of binding agents and PTFE
Housing material	Polycarbonate	Polycarbonate
Color code	Transparent	Transparent
Filter diameter	64 mm	64 mm
Connector inlet	Female luer lock or stepped hose nipple with 6 mm to 12 mm outer diameter	Female luer lock or stepped hose nipple with 6 mm to 12 mm outer diameter
Connector outlet	Hose nipple	Hose nipple
Filling bell	Yes	Yes
Filtration area	20 cm ²	20 cm ²
Hold-up volume before bubble point	1 mL	Approx 1.5 mL
Housing burst pressure	> 5 bar 72.5 psi	> 5 bar 72.5 psi
Bubble point	≥ 3.2 bar 46.4 psi	≥ 3.2 bar 46.4 psi
Max. recommended inlet pressure	3 bar 43.5 psi	3 bar 43.5 psi
Flow rate for water	≥ 250 mL/min at ∆p = 1 bar 14.5 psi	≥ 250 mL/min at ∆p = 1 bar 14.5 psi
Filtration range	100 mL to max. 5 L	100 mL to max. 10 L
pH-range	4–8	4–8
Non-specific protein adsorption	No loss of protein detectable (filtration of γ globulin, method acc. to Bradford)	$<$ 80 μg/cm 2 (filtration of γ globulin, method acc. to Bradford)
Sterilization	EO sterilization	EO sterilization
Biosafety	Class VI Plastics Test	Class VI Plastics Test
Operating instructions	Directions for use included in each box	Directions for use included in each box

Technical Specifications for Sartolab® P20 and Sartolab® P20 Plus Units with PES Membrane

	Sartolab® P20 18075	Sartolab® P20 Plus 18068	Sartolab® P 20 Plus 18076
Filter material	PES, type 15407 MI 0.2 μm pore size	PES, 0.45 μm pore size plus GF, 100% free of binding agents	PES, 0.45 μm pore size plus GF, 100% free of binding agents
Housing material	Polycarbonate	Polycarbonate	Polycarbonate
Color code	Transparent	Transparent	Transparent
Filter diameter	61 mm	61 mm	61 mm
Connector inlet	Female luer lock or stepped hose nipple	Female luer lock or stepped hose nipple	Stepped hose nipple with 6 mm to 12 mm outer
Connector outlet	with 6 mm to 12 mm outer diameter Hose nipple	with 6 mm to 12 mm outer diameter Hose nipple	diameter Hose nipple
Filling bell	No	Yes	No
Filtration area	20 cm ²	20 cm ²	20 cm ²
Hold-up volume before bubble point	1 mL	Approx 1.5 mL	1 mL
Housing burst pressure	> 5 bar 72.5 psi	> 5 bar 72.5 psi	> 5 bar 72.5 psi

Properties	Description		
	Sartolab® P20 Plus 18068	Sartolab® P20 18075	Sartolab® P 20 Plus 18076
Bubble point	≥ 3.2 bar 46.4 psi	≥ 3.2 bar 46.4 psi	≥ 3.2 bar 46.4 psi
Max. recommended inlet pressure	3 bar 43.5 psi	3 bar 43.5 psi	3 bar 43.5 psi
Flow rate for water	400 mL/min at Δp = 1 bar 14.5 psi	400 mL/min at Δp = 1 bar 14.5 psi	> 600 mL/min at Δp = 1 bar 14.5 psi
Filtration range	100 mL – max. 10 L	100 mL – max. 5 L	100 mL – max. 5 L
pH range	1 to 8	1 to 8	1 to 8
Non-specific protein adsorption	< 80 μ g/cm ² (filtration of γ globulin, method acc. to Bradford)	No loss of protein detectable (filtration of γ globulin, method acc. to Bradford)	No loss of protein detectable (filtration of y globulin, method acc. to Bradford
Sterilization	EO sterilization	EO sterilization	-
Biosafety	Class VI Plastics Test	Class VI Plastics Test	Class VI Plastics Test
Operating instructions	Directions for use included in each box	Directions for use included in each box	Directions for use included in each box

Ordering Information

Type	Membrane	Housing	Pore Size	Inlet	Outlet	Sterile	Qty./Pkg.	Order No.
Sartolab®	P20							
	SFCA	PC	0.2 μm	Hose Barb	Hose Barb	Yes	10	18052D
	SFCA	PC	0.2 μm	Luer Lock	Hose Barb	Yes	10	18053D
	PES	PC	0.2 μm	Luer Lock	Hose Barb	Yes	10	18075D ¹
Sartolab®	P20 Plus							
	SFCA + GF	PC	0.2 μm	Hose Barb	Hose Barb	Yes	10	18056D
	SFCA + GF	PC	0.2 μm	Luer Lock	Hose Barb	Yes	10	18058D
	PES + GF	PC	0.2 μm	Luer Lock	Hose Barb	Yes	10	18068D
	PES + GF	PC	0.45 μm	Hose Barb	Hose Barb	No	100	18076N ¹
	GF*	PC	n.a.	Luer Lock	Hose Barb	No	10	18072D ¹

 ${\sf SFCA-Cellulose\ Acetate,\ PES-Polyethersulfone,\ GF-Glass\ Fiber\ Prefilter,\ PC-Polycarbonate\ GF^*-Glass\ Fiber\ Filter\ only}$

¹ without filling bell

Sartolab® RF | BT

Disposable Filter Systems and Bottle-Top Filters





Sartolab® single-use sterile filter systems and bottle-top filters are designed for the vacuum filtration of tissue culture media and components, biological fluids, and other aqueous solutions.

The Sartolab® single-use 150 mL, 250 mL, 500 mL and 1,000 mL bottles are designed as storage containers for sterile media, buffers, or other aqueous solutions.

These products are for laboratory use only and not for human parenteral applications.

Materials

The filter funnels, dust covers and receiver bottles are manufactured from virgin, heavy metal-free polystyrene. The tubing adapters, filter adapters, and the plug seal caps are made of heavy metal-free polyethylene. Sartolab® filter systems are available with polyethersulfone and cellulose acetate membranes. All units are sterilized by gamma irradiation.

Performance

The filter units contain membranes integrally sealed to a support grid designed to maximize flow and reduce foaming and protein denaturation.

The membranes are compatible with most aqueous solutions and tested for use in cell culture applications.

Filter Systems

The filter adapter utilizes a gasket design to ensure a vacuum-tight seal on the receiver storage bottle. Each filter unit also contains a convenient tubing adapter that will fit most vacuum hoses.

The bottles are single-use containers. They cannot withstand autoclaving or use at temperatures greater than 70°C. The suitability of the bottles for storage of solutions below 0°C depends both on the solution and the storage conditions. Many aqueous solutions, including culture media, have been successfully frozen and stored at temperatures down to -20°C. However, a trial run under actual conditions is strongly recommended to test the suitability of the bottles for frozen storage.

Specifications

Pore Size	Membrane Material	Characteristics
0.1 μm	Polyethersulfone	Very low protein binding and low extractables, fast flow rate
0.2 μm	Cellulose Acetate	Optimized for aqueous liquids, very low nonspecific binding
0.45 μm	Polyethersulfone	Very low protein binding and low extractables, fast flow rate

Bottle-Top Filters

The filter adapter is available with a 45 mm thread finish, and is designed to ensure a vacuum-tight seal on customer-supplied bottles with the appropriate thread finish. Each filter unit also contains a tubing adapter that will fit most vacuum hoses.

Chemical Compatibility

The mechanical strength, color, appearance, and dimensional stability of filter systems, bottle-top filters and plastic bottles are affected to varying degrees by the chemicals with which they come in contact. Specific operating conditions, especially temperature, will also affect their chemical resistance. A table is provided to serve as a general guideline for the chemical resistance of Sartolab® single-use sterile filters and bottles.

Chemical Resistance of Sartolab® Filters

Chemical Class	Membrane (PES)	Membrane (CA)	Housing (PS)
Weak Acids	3	2	1
Strong Acids	3	2	2
Alcohols	2	2	2
Aldehydes	3	3	3
Aliphatic Amines	1	3	3
Aromatic Amines	3	3	3
Bases	3	3	1
Esters	3	3	3
Hydrocarbons	3	2	3
Ketones	3	3	3

Key: 1. recommended

- 2. may be suitable for some applications; a trial run is recommended
- 3. not recommended. PS, polystyrene; PES, polyethersulfone; CA, cellulose acetate.


Ordering Information

Volume	Membrane	Membrane Filter Area Oty./Pkg.		Order No.				
Sartolab® RF Filtration System Including Collection Bottle								
150 mL	0.22 μm PES	18 cm ²	12	180C1E				
250 mL	0.22 μm PES	24 cm ²	12	180C7E				
500 mL	0.22 μm PES	39 cm ²	12	180C2E				
1,000 mL	0.22 μm PES	62 cm ²	12	180C3E				
1,000 mL	0.1 μm PES	62 cm ²	12	180C8E				
250 mL	0.45 μm CA	24 cm ²	12	180A1E				
500 mL	0.45 μm CA	39 cm ²	12	180A2E				
1,000 mL	0,45 μm CA	62 cm ²	12	180A3E				

Sartolab® BT Bottle-Top Filters Without Integrated Collection Bottle

150 mL	0.22 μm PES	18 cm ²	48	180C4K
500 mL	0.22 μm PES	39 cm ²	12	180C5E
1,000 mL	0.22 μm PES	62 cm ²	12	180C6E
500 mL	0.45 μm CA	39 cm ²	12	180A4E

Water Throughput

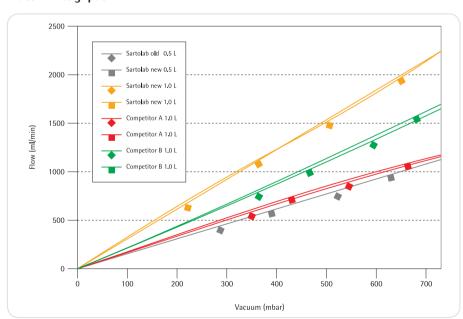

Sartolab® RF 50

Vacuum filtration unit for volumes of up to 50 mL



Sartolab® RF 50 vacuum filtration units are single-use units, designed for the filtration of sample volumes of up to 50 mL. They are available either with a 0.22 or a 0.45 µm polyethersulfone membrane which have been developed and manufactured by Sartorius at its own facilities.

Both membranes provide fastest flow rates due to their asymmetric structure and ensure lowest protein binding as well as low extractables. The 0.22 μ m version is ideal for sterile filtration of cell cultures, buffers and aqueous solutions; the 0.45 μ m for clarification. Both versions are individually-packed, sterile and ready to use.

Sartolab® RF 50 vacuum filtration unit is composed of a funnel, with dust cover, a 50 mL conical tube with graduation and writing field as well as a tubing connector for vacuum connection and a screw cap to store your filtrate in the conical tube. The design of the yellow adapter connecting the funnel to the conical tube ensures a vacuum-tight seal and enables the filtration unit to be used on the Sartolab® MultiStation for filtration of up to 6 samples in parallel with one vacuum source.


Sartolab® RF 50 can also be used alone when connecting the tubing connector delivered with each unit to your vacuum source.

The tubing connector and the screw cap for the conical tube are individually wrapped to maintain sterility until needed.

User Benefits

- Highest flow rates with an asymmetric membrane and a large surface area
- No loss of protein with a very low protein binding membrane
- Low dead volume thanks to an optimized membrane support
- Minimized risks of contamination with a complete ready-to-use unit (no further transfer of liquid for storage necessary)
- Designed as standalone system or for a parallel filtration of up to 6 samples with Sartolab® MultiStation

Filtration of 50 mL cell media + 10 % FBS in a 50 mL vacuum filtration unit (0.22 μ m)

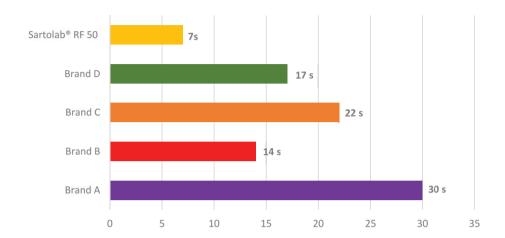

Materials

Funnel with dust cover	Styrene-acrylonitrile (SAN)
Membrane filter	0.22 μm polyethersulfone (order no. 180E01)
Funnel adapter	0.45 μm polyethersulfone (order no. 180F01)
Tubing connector for vacuum connection	Polyethylene
Conical tube	Polyethylene
Conical tube cap	Polypropylene

Specifications

Membrane diameter	58.5 mm
Effective filtration area	21 cm ²
Hold-up volume	1.2 mL
Filtration capacity	50 mL
Size of the 50 mL conical tube	Height: 115.5 mm
	External/internal diameter: 29.5 mm/27.48mm
Autoclavable	No
Sterilization method	E-Beam (beta) irradiation
Storage temperature of the conical tube	4°C to 30°C (short-term: -80°C to max. 100°C)
Packaging	Single-packaged, sterile
Operating pressure	-350 to -700 mbar

Ordering Information

Description	Quantity	Order Number	
Sartolab® RF 50, 0.22 μm, PES	24 units	180E012	
Sartolab® RF 50, 0.45 µm, PES	24 units	180F012	

Sartolab® Multistation

For hands-free parallel filtration of up to 6 samples



Sartolab® MultiStation is a stand specially designed to hold 1 to 6 vacuum filtration units, allowing simultaneous filtration of up to 6 samples.

The MultiStation is permanently connected to your vacuum source. Easily install your vacuum filtration units in the MultiStation for quick and easy filtration of samples without the need for installation of extra connectors and time-consuming stabilization.

Sartolab® MultiStation works with our Sartolab® RF 50 vacuum filtration units; the funnel adapter of these units is designed to fit perfectly in the bracket of the MultiStation.

With one click, connect the filtration unit to the device, assuring perfect filter stability. With a second click, engage the vacuum automatically and begin filtering.


Easily manipulate your samples with the rotating, multi-directional head, and easily keep track of your samples during filtration with the numbered brackets

User Benefits

- Single vacuum source enables simultaneous filtration of up to 6 samples
- Time-saving (no installation time for each filter unit before use)
- Hands-free filtration



The MultiStation advantage: Using the Sartolab® MultiStation to filter 6 x 50 mL samples saves handling time compared to standalone filtration

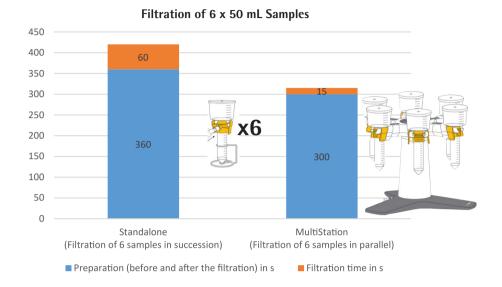

Specifications

Material (visible parts)	ABS/Aluminum/Stainless steel
Dimensions (Length x Width x Height)	307 x 348 x 281 mm
Weight	4.6 kg
Tubing connector	Designed for tubing with an inner diameter between 4 and 10 mm and with a wall thickness of minimum 3 mm

Ordering Information

Description	Quantity	Order Number	_
Sartolab® MultiStation	1 unit	SDLC01	

Sartoclear Dynamics® Lab P15

Clarification and Sterile Filtration of up to 15 mL Mammalian Cell Culture in One Step



Sartoclear Dynamics® Lab P15 is a new kit for single step harvesting of 15 mL animal cell cultures with even high cell densities. With this kit, the clarification and sterile filtration of mammalian cell culture is performed in a single pressure filtration step. Inspired by the plasma industry, Sartoclear Dynamics® is based on the principles of body feed filtration.

This ready to use kit combines a 20 mL syringe pre-filled with a 0.5 g filter aid and an integrated filter* for sterile filtration. The filter aid facilitates filtration through the sterile filter while allowing complete protein recovery.

A convenient filling tube can be connected to the syringe, for the easy recovery of samples from 50 mL Falcon tubes or ambr 15 bioreactors.

As a result, this method replaces centrifugation and subsequent sterile filtration steps, leading to clarified and sterilized cell culture harvest in minutes. Your cell culture harvest will be available for following sample concentration by ultrafiltration and downstream analytics in no time.

Sartoclear® Dynamics Lab P15 Features:

- Single step mammalian cell culture harvest
- Designed for cell densities up to 20×10^6 cells and even more
- Fast and effortless filtration
- Optimized for cell culture harvest from ambr[®] 15 bioreactors
- *The sterile filter included in the kit contains a 0.2 µm polyethersulfone membrane and a prefilter made of 100 % high-purity quartz: The choice of these materials, along with the larger surface area of the filter, enables higher flow rates.





1. Fill the syringe with cell culture broth.



2. Shake the syringe to mix the broth with the filter aid.



3. Connect the sterile filter to the syringe and filter.

Typical Results

Cell Type	Cell Density		Mab concentration before filtration	Mab concentration after filtration	Recovery Rate
CHO DG44	16×10^6 cells/mL	78%	6.02 g/L	5.77 g/L	96%
CHO DG44	38×10^5 cells/mL	48%	0.43 g/L	0.43 g/L	100%

Specifications

DE Syringe

Syringe material	Syringe barrel and plunger rod: polypropylene; stopper: latex-free elastomer
Syringe Cap	Polyamide
Filling tube material	Polypropylene
Filter aid	0.5 g highly pure diatomaceous earth (Celpure® C300 – pharmaceutical-grade*)

Sterile Filter

Sterile Filter	
Housing material	Polycarbonate
Prefilter material	100% high-purity quartz, binder-free
Filter material	0.2 μm polyethersulfone
Filter diameter	61 mm
Filtration area	20 cm ²
Connector inlet	Female luer lock
Connector outlet	Hose barb
Hold-up volume	Approx. 3 mL
Housing burst pressure	> 5 bar 72.5 psi
Packaging	Individually packed
Sterilization	EO sterilization

Ordering Information

Product Name	Number of units per box	Order No.
Sartoclear Dynamics [®] Lab P15	6 × 20 mL syringes, pre-filled with 0.5 g DE, including caps and filling tubes 6 × sterile filters	SDLP0015C

^{*} Celpure® is a trademark of Advanced Minerals

Sartoclear Dynamics® Lab V

Clarification and Sterile Filtration of 50 mL up to 1 L Mammalian Cell Culture in One Step



Sartoclear Dynamics® Lab V kits enable clarification and sterile filtration to be performed in a single step. These kits simplify the cell harvesting process by fully eliminating the centrifugation step otherwise needed for clarification. As a result, they enable cell cultures to be efficiently clarified and sterilized in minutes – quickly and easily.

Sartoclear Dynamics $^{\circ}$ Lab kits have been designed and optimized for harvesting mammalian cell cultures, such as CHO, HEK, hybridomas and many others, with cell densities of up to 20×10^6 cells/mL.

Each kit provides filter aid pouches for clarification and Sartolab® RF vacuum filtration units for sterile filtration.

The filter aid used in Sartoclear Dynamics® Lab products is made of highly-pure diatomaceous earth (DE) that is insoluble and inert. It is packed in ready-to-use pouches in pre-wetted condition to prevent the release of dust particles.

The DE pouches are gamma-irradiated to rule out any contamination.

Sartoclear® Dynamics Lab V Features:

- Single step mammalian cell culture harvest
- Designed for cell densities up to 20×10^6 cells
- Fast and effortless filtration

Quick and Easy Filtration



 Attach the Sartolab® RF unit to a vacuum source. Pour the filter aid into your cell culture.



2. Mix to homogenize. Pour the cell culture and filter aid mixture into the Sartolab® RF first ...



- **3.** ... then filter by simply starting the pump.
- **4.** After 5 15 min, the sample will be available for subsequent downstream processing steps, such as protein purification and concentration.
- **5.** Parallel Filtration of 50 mL samples with the Sartolab[®] MultiStation

Typical Results

Cell Type	Cell Density	Viabi- lity	mAb concentration before filtration	Mab concentration after filtration		Turbidity
СНО	14.46×10 ⁶ cells/mL	85.2%	5.2 g/L	5.15 g/L	99%	18 NTU
HEK	8×10 ⁶ cells/mL	70%	0.035 g/L	0.034 g/L	97%	8 NTU

Specifications

Clarification

Pouches of Filter Aid

Diatomaceous Earth (DE)	1 g, 5 g or 10 g highly pure diatomaceous earth, (Celpure® C300 – pharmaceutical grade)*, mixed with water in a ratio of 1 DE: 1.25 ultrapure water
Packaging Sterilization	Dust-free, gamma irradiated pouches

Filtration

Vacuum filtration units with receiver flasks (Sartolab® RF)

vacaam meracion ames wii	in receiver masks (surtoido m)
Funnel, dust cover, receiver bottles	Heavy metal-free polystyrene
Filter adapter, tubing adapter, cap	Heavy metal-free polyethylene
Filter material	0.22 μm polyethersulfone
Packaging Sterilization	Single-packaged, sterile

Filtration

Vacuum filtration units with conical tube (Sartolab® RF 50)

Funnel, dust cover	Styrene-acrylonitrile (SAN)
Funnel adapter, tubing connector, cap	Heavy metal-free polyethylene
Conical tube	Polypropylene
Filter material	0.22 μm polyethersulfone (order no. 180E01) 0.45 μm polyethersulfone (order no. 180F01)
Packaging	Single-packaged, sterile

^{*} Celpure® is a trademark of Advanced Minerals

Each Sartoclear Dynamics® Lab V kit is comprised of pouches of filter aid and Sartolab® RF vacuum filtration units that match your needs. Find the right kit in just two easy steps:

- 1. Determine the volume range of your sample to be filtered.
- 2. Then reference it to the cell density of your cell culture.

Volume*	Cell density**											
	< 5 million	5-10 million	10-20 million									
	cells/mL	cells/mL	cells/mL									
≤ 50 mL	SDLV-0050-01E0-2	SDLV-005	60-02E0-2									
> 50 –150 mL	SDLV-0150-02CE	SDLV-0150-05C2										
150 – 250 mL	SDLV-0250-05C2	SDLV-025	60-10C2									
250 – 500 mL	SDLV-0500-05C2	SDLV-0500-10C2	SDLV-0500-20CE									
500 –1000 mL	SDLV-1000-10C2	SDLV-1000-20CE	SDLV-1000-40CE									

^{*} Especially for cultures with high cell densities (>10 million cells/mL) > 50, we recommend that you use a Sartolab® RF unit nominally rated for the next higher volume of the solution to be filtered as it will accelerate your filtration run (larger membrane surface area | higher quantity of DE).

^{**} Tested with CHO cell lines with a cell viability of approx. 85%

Ordering Information

Description	Oty. of Units	Order No.
Sartoclear Dynamics® Lab V50 Kits – 0.22 μm PES		
Sartoclear Dynamics Lab V, 50 mL, 1 g		
Filtration of up to 50 mL with 1 g of DE per unit	24	SDLV-0050-01E0-2
Contents: 1 × 180E012		
(24 × Sartolab [®] RF 50, 0.22 μm PES)		
1 × SDLKG-01.02		
(24 × pouches of filter aid, 1 g)		
Sartoclear Dynamics Lab V, 50 mL, 2 g		
Filtration of up to 50 mL with 2 g of DE per unit	24	SDLV-0050-02E0-2
Contents: 1 × 180E012 (24 × Sartolab® RF 50, 0.22 μm PES)		
2 × SDLKG-01.02		
$(48 \times \text{pouches of filter aid, 1 g})$		
Sartoclear Dynamics® Lab V50 Kits – 0.45 μm PES		
Sartoclear Dynamics Lab V, 50 mL, 1 g		
Filtration of up to 50 mL with 1 g of DE per unit	24	SDLV-0050-01F0-2
Contents: 1 × 180F012		551, 5555 5.1.5 2
(24×Sartolab® RF 50, 0.45 μm PES)		
1×SDLKG-01.02		
(24×pouches of filter aid, 1 g)		
Sartoclear Dynamics Lab V, 50 mL, 2 g		
Filtration of up to 50 mL with 2 g of DE per unit	24	SDLV-0050-02F0-2
Contents: 1 × 180F01E		
(24 × Sartolab® RF 50, 0.45 μm PES)		
2 × SDLKG-01.02		
(48 \times pouches of filter aid, 1 g)		
Sartoclear Dynamics® Lab V150 Kits		
Sartoclear Dynamics® Lab V, 150 mL, 2 g		
Filtration of up to 150 mL with 2 g of DE per unit	12	SDLV-0150-02CE
Contents: 1 × 180C1E		
(12 × Sartolab® RF 150, 0.22 μm PES)		
1 × SDLKG-01.02		
(24 × pouches of filter aid, 1 g)		
Sartoclear Dynamics® Lab V, 150 mL, 5 g	24	CDIV 0150 050 2
Filtration of up to 150 mL with 5 g of DE per unit Contents: 2 × 180C1E	24	SDLV-0150-05C2
(24× Sartolab [®] RF 150, 0.22 μm PES)		
1×SDLKG-05.02		
(24×pouches of filter aid, 5 g)		
· 1		

Description

Qty.

Order No.

Description	of Units	Order No.
Sartoclear Dynamics® Lab V250 Kits		
Sartoclear Dynamics® Lab V, 250 mL, 5 g		
Filtration of up to 250 mL with 5 g of DE per unit	24	SDLV-0250-05C2
Contents: 2 × 180C7E		
(24×Sartolab® RF 250, 0.22 μm PES)		
1×SDLKG-05.02		
(24×pouches of filter aid, 5 g)		
Sartoclear Dynamics® Lab V, 250 mL, 10 g		
Filtration of up to 250 mL with 10 g of DE per unit	24	SDLV-0250-10C2
Contents: 2 × 180C7E (24×Sartolab® RF 250, 0.22 μm PES)		
(24× Sartolao - KF 250, 0.22 μm FES) 1×SDLKG-10.02		
(24×pouches of filter aid, 10 g)		
Sartoclear Dynamics® Lab V500 Kits		
Sartoclear Dynamics® Lab V, 500 mL, 5 q		
Filtration of up to 500 mL with 5 g DE per unit	24	SDLV-0500-05C2
Contents: 2 × 180C2E		, 000 2
(24 \times Sartolab® RF 500, 0.22 μ m PES)		
1 × SDLKG-05.02		
(24 × pouches of filter aid, 5 g)		
Sartoclear Dynamics® Lab V, 500 mL, 10 g		
Filtration of up to 500 mL with 10 g DE per unit	24	SDLV-0500-10C2
Contents: 2 × 180C2E		
(24 × Sartolab® RF 500, 0.22 μm PES) 1 × SDLKG-010.02		
$(24 \times \text{pouches of filter aid, } 10 \text{ g})$		
Sartoclear Dynamics® Lab V, 500 mL, 20 g		
Filtration of up to 500 mL with 20 g DE per unit	12	SDLV-0500-20CE
Contents: 1 × 180C2E		551, 0000 100 1
(12 \times Sartolab® RF 500, 0.22 μ m PES)		
1 × SDLKG-010.02		
$(24 \times pouches of filter aid, 10 g)$		
Sartoclear Dynamics® Lab V1000 Kits		
Sartoclear Dynamics® Lab V, 1,000 mL, 10 g		
Filtration of up to 1 L with 10 g of DE per unit Contents: 2 × 180C3E	24	SDLV-1000-10C2
(24 \times Sartolab $^{\circ}$ RF 1000, 0.22 μ m PES)		
1×SDLKG-10.02		
(24×pouches of filter aid, 10 g)		
Sartoclear Dynamics® Lab V, 1,000 mL, 20 g		
Filtration of up to 1 L with 20 g of DE per unit	12	SDLV-1000-20CE
Contents: 1 × 180C3E		
(12×Sartolab® RF 1000, 0.22 μm PES) 1×SDLKG-10.02		
$(24 \times \text{pouches of filter aid, } 10 \text{ g})$		
Sartoclear Dynamics® Lab V, 1,000 mL, 40 g		
Filtration of up to 1 L with 40 g of DE per unit	12	SDLV-1000-40CE
Contents: 1×180C3E	12	2014-1000-40CE
(12×Sartolab® RF 1000, 0.22 μm PES) 2×SDLKG-10.02		

Chemical Compatibility

	Mat	erial							Minisart® Types									
	PES membrane	SFCA membrane	PTFE membrane	RC membrane	Nylon membrane	GF depth filter	Housing MBS	Housing PP	Minisart [®] HighFlow	Minisart® NML Ophthalsart	Minisart® NML Plus	Minisart® NML GF	Minisart® HY Minisart® Air	Minisart® RC	Minisart® NY	Minisart® NY Plus	Minisart® SRP	Minisart® PES
Filter Membrane	PES	SFCA	PTFE	RC	PA				PES	SFCA			PTFE	RC	PA	PA	PTFE	PES
Pre-Filter						GF			-	-	GF	GF	-	_	-	GF	-	-
Housing Material							MBS	PP	MBS	MBS	MBS	MBS	MBS	PP	PP	PP	PP	PP
Sterilization																		
Ethylene oxide	++	++	++	++	++	++	++	++	++	++	++	++	++	++	++	++	++	++
Gamma irradiation	++	++	_1	++	-	++	++	-	++	++	++	++	_1	-	-	-	-	-
Autoclaving 121°C, 30 min	++	++	++	++	++	++	-	++	-		-			++	++	++	++	++
Solvents																		
Acetone	-		++	++	++	++		++	-	_		_		++	++	++	++	
Acetonitrile	-		++	++	++	++		++	-	-		_		++	++	++	++	
Gasoline	+	++	++	++	++	++	+	++	+	+	+	+	+	++	++	++	++	+
Benzene	+	+	++	++	++	++		++	-	-	_	-	-	++	++	++	++	+
Benzyl alcohol	+	+	++	++	++	++		+	-	-		_	-	++	++	++	++	+
n-Butyl acetate	-	_	++	++	++	++		++	-	_	_			++	++	++	++	
n-Butanol	++	++	++	++	++	++	+	++	+	+	+	+	+	++	++	++	++	++
Cellosolve	+	_	++	++	++	++	_	++	-	-	_	-	_	++	++	++	++	+
Chloroform	-	_	++	++	++	++	_	++	-	_	_	_		++	++	++	++	
Cyclohexane	-	_	++	++	++	++	+	+	-	_	_	+	+	+	+	+	+	
Cyclohexanone	-		++	++	++	++		+	-		_	_	_	+	+	+	+	
Diethylacetamide	-	_	++	++	++	++	_	++	-		_			++	++	++	++	-
Diethyl ether	-	+	++	++	++	++	_	++	-		-			++	++	++	++	-
Dimethyl formamide	-		++	+	+	++	_	++	-		_			+	+	+	++	
Dimethylsulfoxide	-		++	++	++	++		++	-		_			++	++	++	++	
Dioxane	-		++	++	++	++	_	++	-		_			++	++	++	++	
Ethanol, 98%	++	++	++	++	++	++		++	-		-			++	++	++	++	++
Ethyl acetate	-	-	++	++	++	++	-	+	-	-	-	-	-	+	+	+	+	-
Ethylene glycol	++	+	++	++	++	++	++	++	++	+	+	++	++	++	++	++	++	++
Formamide	++		++	+	++	++	++	++	++	-	-	++	++	+	++	++	++	++
Glycerin	++	++	++	++	++	++	++	++	++	++	++	++	++	++	++	++	++	++
n-Heptane	+	+	++	++	++	++	++	+	+	+	+	+	++	+	+	+	+	+
n-Hexane	+	+	++	++	++	++	++	+	+	+	+	+	++	+	+	+	+	+
Isobutanol	++	+	++	++	++	++	_	++	-	-	_	-	-	++	++	++	++	++
Isopropanol	++	++	++	++	++	++		++	-	_	_	-		++	++	++	++	++
Isopropyl acetate	-	_	++	++	++	++		++	-	_	_	-		++	++	++	++	
Methanol, 98%	+	_	++	++	++	++	_	++	-	-	_	-		++	++	++	++	+
Methyl acetate	-	-	++	++	++	++	_	+	-	-	-	-		+	+	+	+	
Methylene chloride	-	-	++	++	++	++	_	++	-	-	-	-	-	++	++	++	++	
Methyl ethyl ketone	-	+	++	++	++	++	_	+	-	-	-	-	-	+	+	+	+	-
Methyl isobutyl ketone	-	-	++	++	++	++	_	+	-	-	-	-	-	+	+	+	+	-
Monochlorobenzene	+	+	++	++	++	++	-	+	-	-	-	-	-	+	+	+	+	+
Nitrobenzene	-	-	++	++	+	++	-	+	-	-	-	-	-	+	+	+	+	-
n-Pentane	++	++	++	++	++	++	+	+	+	+	+	+	+	+	+	+	+	+
Perchloroethylene	-	-	++	++	++	++	-	+	-	-	_	-	-	+	+	+	+	-

	Mat	terial							Minisart® Types									
	PES membrane	SFCA membrane	PTFE membrane	RC membrane	Nylon membrane	GF depth filter	Housing MBS	Housing PP	Minisart ^{® HighFlow}	Minisart [®] NML Ophthalsart	Minisart® NML Plus	Minisart® NML GF	Minisart® HY Minisart® Air	Minisart® RC	Minisart® NY	Minisart® NY Plus	Minisart® SRP	Minisart® PES
Filter Membrane	PES	SFC/	A PTFE	RC	PA				PES	SFCA	SFCA		PTFE	RC	PA	PA	PTFE	PES
Prefilter						GF			-	-	GF	GF	-	-	-	GF	-	_
Housing Material							MBS	PP	MBS	S MBS	MBS	MBS	MBS	PP	PP	PP	PP	PP
Solvents (continued)																		
Pyridine	-	-	++	++	++	++	-	++	-	-	-	-	-	++	++	++	++	-
Carbon tetrachloride	-	-	++	++	++	++	-	+	-	-	-	-	-	+	+	+	+	-
Tetrahydrofuran	-	-	++	++	++	++	-	++	-	-	-	-	-	++	++	++	++	_
Toluene	-	+	++	++	++	++	-	+	-	-	-	-	-	+	+	+	+	_
Trichloroethane	-	-	++	++	+	++	-	+	-	-	-	-	-	+	+	+	+	_
Trichloroethylene	-	+	++	++	++	++	-	+	-	-	-	-	-	+	+	+	+	_
Xylene	-	+	++	++	++	++	-	+	-	-	-	-	-	+	+	+	+	_
Acids																		
Acetic acid, 25%	+	+	++	++	-	++	+	++	+	+	+	+	+	++	-	-	++	+
Acetic acid, 80%	-	-	++	+	-	++	-	+	-	-	-	-	-	+	-	-	+	_
Hydrofluoric acid, 25%	+	-	++	+	-	++	+	+	+	-	-	+	+	+	-	-	+	+
Hydrofluoric acid, 50%	+	-	++	+	-	++	-	+	-	-	-	-	-	+	-	-	+	+
Perchloric acid, 25%	-	-	++	-	-	++	-	+	-	-	-	-	-	-	-	-	+	_
Phosphoric acid, up to 10%	+	+	++	-	-	++	+	+	+	+	+	+	+	-	-	-	+	+
Phosphoric acid, 86%	+	+	++	_	_	++	_	+	_	_	_	_	_	_	_	_	+	+
Nitric acid, 30%	+	_	++	_	_	++	+	+	+	_	_	+	+	_	_	_	+	+
Nitric acid, conc.	_	_	++	_	_	++	_	_	_	_	_	_	_	_	_	_	_	_
Hydrochloric acid, 15%	++	+	++	_	_	++	+	+	+	+	+	+	+	_	_	_	+	+
Hydrochloric acid, 20%	++	_	++	_	_	++	+	+	+	_	_	+	+	_	_	_	+	+
Sulfuric acid, 25%	+	_	++	+	_	++	++	++	+	_	_	++	++	+	_	_	++	+
Sulfuric acid, 98%	_	_	++	_	_	++	_	_	_	_	_	_	_	_	_	_	_	_
Trichloroacetic acid, 25%	_	_	++	++	_	++	_	+	_	_	_	_	_	+	_	_	+	_
Bases																		
Ammonia, 1N	++	+	++	+	++	++	+	++	+	+	+	+	+	+	++	++	++	++
Ammonium hydroxide, 25%	+	+	++	+	++	+	_	+	_	_	_		_	+	+	+	+	+
Potassium hydroxide, 32%	++		++	_	+	+	_	++	_	_	_	_	_	<u> </u>	+	+	++	++
Sodium hydroxide, 1N	++	_	++	+	++	+	_	++	_	_	_	_	_	+	++	+	++	++
Sodium hydroxide, 32%	++	_	++		+		_	+	_	_	_	_	_	<u> </u>	+		+	+
Aqueous solutions					•			•									•	
Formaldehyde, 30%	+	++	++	+	++	++	+	+	+	+	+	+	+	+	+	+	+	+
Sodium hypochlorite, 5%	++	-	++			++	+	+	+		<u> </u>	+	+	<u> </u>	<u> </u>	<u> </u>	+	+
Hydrogen peroxide, 35%	++	_	++	_	_	++	+	++	+		_	+	+	_		_	++	++
pH range							•		<u>'</u>			•						
pH 1 to 14	_	_	++	_	_	++	-				Lec	end						
pH 1 to 13	++		++	_	_	++	_	++				<u> </u>	atibilit	у			++	_
pH 3 to 14		_		+	++	++	_	++	Limited compatibility +								_	
pH 3 to 12	++	_	++	++	++	++	+	++	Not compatible -									
pH 4 to 8	++	++	++	++	++	++	++	++			¹ G	amma i	rradiati	on fea	sible f	or Mir	nisart®	Air
ριι τιο ο	ГТ	гт	ГТ	ГТ	гт	ГТ	гт	7										

Contact time: 24 hours at 20°C. Chemical compatibilities can be influenced by various factors. Therefore, we recommend that you confirm compatibility with the liquid you want to filter by performing a trial filtration run before you start your actual filtration.



Basic Filtration

Table of Contents

- 112 Glass and Quartz Microfiber Filters
- 124 Blotting | Chromatography Papers & Blotting Membranes
- 126 Filtration Equipment





Filters are indispensable for your routine work in laboratory and industrial applications. Sartorius supplies you with a broad range of filters for a myriad of filtration tasks and supports you with all your filtration challenges.

Our Product Range Covers:

- Filter papers
- Glass and quartz microfiber filters
- Membrane filters
- Blotting & chromatography papers & membranes
- Filtration equipment

Quality Assurance and Quality Control

Sartorius pays particular attention to continuous in-process quality control. Regular checks and exact analyses of the raw materials and each finished product assure constant high quality and product uniformity. We meet the requirements set forth by the ISO 9001 quality management system and the ISO 14001 environmental management system.

How Do Filter Papers Work?

Filter papers are depth filters. Their efficiency is influenced by various parameters: the mechanical particulate retention, adsorption, pH, surface properties, thickness and strength of the filter paper as well as the shape, density and quantity of particles to be retained. The precipitates deposited on the filter form a "cake layer" which – depending on its density – increasingly affects the progress of an ongoing filtration and decisively affects the retention capability. Therefore, it is essential to select the perfect filter paper to ensure the best filtration results. This choice depends on the filtration method as well as on the amount and properties of the medium to be filtered, the size of the particulate solids to be removed and the required degree of clarification.

How Do Membrane Filters Work?

Membrane filters retain particles larger than their pore sizes. Smaller particles pass through the membrane or are captured in the membrane. Such filters are used for the filtration of smaller particles and for critical applications such as sterility testing. The choice of the right membrane type depends on the specifications of the solution to be filtered. The most important parameters for this are adsorption, chemical compatibility and the particle size to be retained.



Ash-free Filter Papers

For Quantitative and Gravimetric Analyses

These filter papers are used for quantitative and gravimetric analyses as well as for pressure or vacuum filtration. They are made out of 100% cotton linters with an α -cellulose content of > 98% and are acid-washed to make the papers ashless and achieve high purity.

Typical Values

Grade	Weight (g/m²)	Thickness (mm)	Particle Retention (μm)	Filtration (s)	Precipitates	Properties
388	84	0.21	12 to 15	10	Coarse crystalline	Wide-pore, loose structure, fast filtering
389 389	84	0.19	8 to 12	20	Medium-fine crystalline	Medium- to wide-pore, medium fast filtering
389 F	84	0.19	8 to 12	20	Medium-fine crystalline	Medium- to wide-pore, medium fast filtering
392 392	84	0.17	5 to 8	50	Fine crystalline	Medium dense, medium fast filtering
390 390	84	0.16	3 to 5	100	Fine crystalline	Narrow-pore, dense, slow filtering
391	84	0.15	2 to 3	180	Very fine crystalline	Fine-pore, dense, very slow filtering
393 393	100	0.17	1 to 2	300	Very fine crystalline	Very fine-pore, very dense, very slow filtering


Ordering Information

Filter Discs, 100 Pieces

Dia. in mm	Grade 388	Grade 389	Grade 389 F	Grade 390	Grade 391	Grade 392	Grade 393
55	FT-3-101-055	FT-3-102-055	FT-3-112-055	FT-3-103-055	FT-3-104-055	FT-3-105-055	FT-3-127-055
70	FT-3-101-070	FT-3-102-070		FT-3-103-070	FT-3-104-070	FT-3-105-070	FT-3-127-070
90	FT-3-101-090	FT-3-102-090	FT-3-112-090	FT-3-103-090	FT-3-104-090	FT-3-105-090	FT-3-127-090
110	FT-3-101-110	FT-3-102-110	FT-3-112-110	FT-3-103-110	FT-3-104-110	FT-3-105-110	FT-3-127-110
125	FT-3-101-125	FT-3-102-125	FT-3-112-125	FT-3-103-125	FT-3-104-125	FT-3-105-125	FT-3-127-125
150	FT-3-101-150	FT-3-102-150	FT-3-112-150	FT-3-103-150	FT-3-104-150	FT-3-105-150	FT-3-127-150
185	FT-3-101-185	FT-3-102-185	FT-3-112-185	FT-3-103-185	FT-3-104-185	FT-3-105-185	FT-3-127-185
240	FT-3-101-240	FT-3-102-240		FT-3-103-240	FT-3-104-240	FT-3-105-240	FT-3-127-240

Folded Filters, 100 Pieces

Dia. in mm	Grade 388	Grade 389	Grade 389 F	Grade 390	Grade 391	Grade 392
110	FT-4-101-110	FT-4-102-110		FT-4-103-110	FT-4-104-110	FT-4-105-110
125	FT-4-101-125	FT-4-102-125		FT-4-103-125	FT-4-104-125	FT-4-105-125
150	FT-4-101-150	FT-4-102-150		FT-4-103-150	FT-4-104-150	FT-4-105-150
185	FT-4-101-185	FT-4-102-185	FT-4-112-185	FT-4-103-185	FT-4-104-185	FT-4-105-185
240	FT-4-101-240	FT-4-102-240			FT-4-104-240	

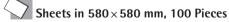

Grade 388	Grade 389	Grade 390	Grade 391	Grade 392	Grade 393
FT-2-101-580580	FT-2-102-580580	FT-2-103-580580	FT-2-104-580580	FT-2-105-580580	FT-2-127-580580

Filter Papers Basic Filtration

109

Wet-strengthened Filter Papers

For Qualitative Analyses

These qualitative filter papers are essentially used for analytical purposes and routine analyses, whenever no gravimetric analyses are required. They are wet-strengthened and can be used for pressure and vacuum filtration. They are made of refined pulp and linters with an >95% α -cellulose content and are very pure with an ash content $\le 0.1\%$.

Typical Values

Grade	Weight (g/m²)	Thickness (mm)	Particle Retention (μm)	Filtration (s)	Precipitates	Properties
1288	84	0.21	12 to 15	10	Coarse crystalline	Wide-pore, loose structure, fast filtering
1289	84	0.21	8 to 12	20	Medium-fine crystalline	Medium- to wide-pore, medium fast filtering
1292	84	0.17	5 to 8	50	Fine crystalline	Medium dense, medium fast filtering
1290	84	0.15	3 to 5	100	Fine crystalline	Narrow-pore, dense, slow filtering
1291	84	0.15	2 to 3	180	Very fine crystalline	Fine-pore, dense, very slow filtering
293	80	0.15	1 to 2	300	Very fine crystalline	Very fine-pore, very dense, very slow filtering

Ordering Information

Filter Discs, 100 Pieces

Dia. in mm	Grade 1288	Grade 1289	Grade 1290	Grade 1291	Grade 1292	Grade 293
55	FT-3-206-055	FT-3-207-055	FT-3-208-055	FT-3-209-055	FT-3-210-055	FT-3-211-055
70	FT-3-206-070	FT-3-207-070	FT-3-208-070	FT-3-209-070	FT-3-210-070	FT-3-211-070
90	FT-3-206-090	FT-3-207-090	FT-3-208-090	FT-3-209-090	FT-3-210-090	FT-3-211-090
110	FT-3-206-110	FT-3-207-110	FT-3-208-110	FT-3-209-110	FT-3-210-110	FT-3-211-110
125	FT-3-206-125	FT-3-207-125	FT-3-208-125	FT-3-209-125	FT-3-210-125	FT-3-211-125
150	FT-3-206-150	FT-3-207-150	FT-3-208-150	FT-3-209-150	FT-3-210-150	FT-3-211-150
185	FT-3-206-185	FT-3-207-185	FT-3-208-185	FT-3-209-185	FT-3-210-185	FT-3-211-185
240	FT-3-206-240	FT-3-207-240	FT-3-208-240	FT-3-209-240	FT-3-210-240	

Folded Filters, 100 Pieces

Dia. in mm	Grade 1288	Grade 1289	Grade 1290	Grade 1291	Grade 1292	Grade 293
110	FT-4-206-110	FT-4-207-110	FT-4-208-110	FT-4-209-110	FT-4-210-110	
125	FT-4-206-125	FT-4-207-125	FT-4-208-125	FT-4-209-125	FT-4-210-125	FT-4-211-125
150	FT-4-206-150	FT-4-207-150	FT-4-208-150	FT-4-209-150	FT-4-210-150	FT-4-211-150
185	FT-4-206-185	FT-4-207-185	FT-4-208-185	FT-4-209-185	FT-4-210-185	FT-4-210-185
240	FT-4-206-240	FT-4-207-240	FT-4-208-240	FT-4-209-240	FT-4-210-240	FT-4-211-240



Grade 1288	Grade 1289	Grade 1290	Grade 1291	Grade 1292	Grade 293
FT-2-206-580580	FT-2-207-580580	FT-2-208-580580	FT-2-209-580580	FT-2-210-580580	FT-2-211-580580

Other dimensions are available on request

High-Purity Filter Papers

For Qualitative Analyses

These paper grades are used for analytical purposes that require a low ash content. Grades 292 and 292a are especially suitable for soil analyses because they are low in nitrogen. For phosphate or sodium determination, we recommend grades 131 and 132.

Typical Values

Grade	Weight (g/m²)	Thickness (mm)	Particle Retention (μm)	Filtration (s)	Material
292	87	0.18	5 to 8	45	Cotton linters, low-nitrogen and nitrates, ash content ≤ 0.06% according to DIN 54370
292a	97	0.19	4 to 7	60	Cotton linters, low-nitrogen and nitrates, ash content ≤ 0.06% according to DIN 54370
132	80	0.17	5 to 7	55	Cotton linters and refined pulp, low-phosphate and low-potassium, ash content < 0.02% according to DIN 54370
131	80	0.16	3 to 5	100	Cotton linters and refined pulp, low-phosphate and low-potassium, ash content < 0.02% according to DIN 54370

Ordering Information

Filter Discs, 100 Pieces

Dia. in mm	Grade 131	Grade 132	Grade 292	Grade 292a
55		FT-3-329-055	FT-3-205-055	FT-3-215-055
70		FT-3-329-070	FT-3-205-070	FT-3-215-070
90		FT-3-329-090	FT-3-205-090	FT-3-215-090
110		FT-3-329-110	FT-3-205-110	FT-3-215-110
125	FT-3-351-125	FT-3-329-125	FT-3-205-125	FT-3-215-125
150		FT-3-329-150	FT-3-205-150	FT-3-215-150
185		FT-3-329-185	FT-3-205-185	FT-3-215-185
240		FT-3-329-240	FT-3-205-240	FT-3-215-240

Folded Filters, 100 Pieces

Dia. in mm	Grade 131	Grade 132	Grade 292	Grade 292a
110	FT-4-351-110	FT-4-329-110	FT-4-205-110	FT-4-215-110
125	FT-4-351-125	FT-4-329-125	FT-4-205-125	FT-4-215-125
150	FT-4-351-150	FT-4-329-150	FT-4-205-150	FT-4-215-150
185	FT-4-351-185	FT-4-329-185	FT-4-205-185	FT-4-215-185
240		FT-4-329-240	FT-4-205-240	FT-4-215-240

Sheets in 580×580 mm, 100 Pieces Grade 292 Grade 292a

FT-2-205-580580 FT-2-215-580580

Other dimensions are available on request

Filter Papers Basic Filtration

111

Filter Papers

For Qualitative-Technical Analyses

These filter papers are used for routine analyses like clarification, determination of substances, but also as discs with a center hole for technical applications. Grades with a wet burst resistance > 30 kPa are referred to as wet-strengthened and are therefore suitable for pressure or vacuum filtration. They are made of refined pulp and linters with an > 95% α -cellulose content, are very pure with an ash content between <0.1 to 0.15%. Below you will find an overview of the most commonly used grades.

Typical Values

Grade	Surface	Weight (g/m²)	Thickness (mm)	Particle Retention (µm)	Filtration (s)	Wet Burst Resistance (kPa)	Properties
3 hw	Smooth	65	0.14	8 to 12	20	40	Medium fast filtering, filter paper for routine work in the lab
4 b	Smooth	75	0.15	8 to 12	22	> 15	Medium fast filtering, filtration of coarse precipitates, wick paper for seed testing
603/N	Crêped	75	0.25	> 15	8	≥ 50	Fast filtering, filtration of sugar solutions
6	Smooth	80	0.17	10 to 13	15	30	Fast filtering, degassing beer before analysis, clarification of spirits
100/N	Smooth	85	0.18	6 to 8	30	80	Medium fast filtering, ash content < 0.1%, low potassium and sodium content, determination of the sugar content
5 H/N	Crêped	85	0.28	> 40	3	≥ 40	Very fast filtering, wide-pore, filtration of essential oils
3 S/h	Smooth	200	0.36	5 to 7	55	15	Medium fast to slow filtering, narrow-pore, re-wet test for diapers

Ordering Information

Filter Discs

Dia. in mm	Grade 3 hw (100 Pieces)	Grade 4 b (100 Pieces)	Grade 603/N (100 Pieces)	Grade 6 (100 Pieces)	Grade 100/N (100 Pieces)	Grade 5 H/N (100 Pieces)	Grade 3 S/h (50 Pieces)
55	FT-3-303-055	FT-3-309-055		FT-3-312-055	FT-3-328-055		FT-3-307-055
70	FT-3-303-070	FT-3-309-070		FT-3-312-070	FT-3-328-070		
90	FT-3-303-090	FT-3-309-090	FT-3-335-090	FT-3-312-090	FT-3-328-090	FT-3-423-090	FT-3-307-090
110	FT-3-303-110	FT-3-309-110	FT-3-335-110	FT-3-312-110	FT-3-328-110		FT-3-307-110
125	FT-3-303-125	FT-3-309-125	FT-3-335-125	FT-3-312-125	FT-3-328-125	FT-3-423-125	FT-3-307-125
150	FT-3-303-150	FT-3-309-150	FT-3-335-150	FT-3-312-150	FT-3-328-150	FT-3-423-150	FT-3-307-150
185	FT-3-303-185	FT-3-309-185	FT-3-335-185	FT-3-312-185	FT-3-328-185	FT-3-423-185	FT-3-307-185
240	FT-3-303-240	FT-3-309-240	FT-3-335-240	FT-3-312-240	FT-3-328-240	FT-3-423-240	FT-3-307-240

Folded Filters, 100 Pieces

Dia. in mm	Grade 3 hw	Grade 4 b	Grade 603/N	Grade 6	Grade 100/N	Grade 5 H/N
125	FT-4-303-125	FT-4-309-125	FT-4-335-125	FT-4-312-125		FT-4-423-125
150	FT-4-303-150	FT-4-309-150	FT-4-335-150	FT-4-312-150	FT-4-328-150	FT-4-423-150
185	FT-4-303-185	FT-4-309-185	FT-4-335-185	FT-4-312-185		FT-4-423-185
240	FT-4-303-240	FT-4-309-240	FT-4-335-240	FT-4-312-240	FT-4-328-240	FT-4-423-240
270	FT-4-303-270	FT-4-309-270	FT-4-335-270	FT-4-312-270	FT-4-328-270	FT-4-423-270
320	FT-4-303-320	FT-4-309-320	FT-4-335-320	FT-4-312-320	FT-4-328-320	FT-4-423-320




Sheets in 580×580 mm, 100 Pieces

Grade 3 hw	Grade 4 b	Grade 603/N	Grade 6	Grade 100/N	Grade 5 H/N
FT-2-303-580580	FT-2-309-580580	FT-2-335-580580	FT-2-312-580580	FT-2-328-580580	FT-2-423-580580

Glass Microfiber Filters

Without Binder

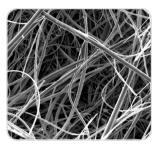

Binder-free glass microfiber filters are recommended for analytical and gravimetric analyses and also as prefilters. These filters combine fast flow rates with high load capacity and the retention of very fine

particles; they are biologically inert, are resistant to most chemicals and withstand temperatures up to 500°C (grade 550-HA up to 550°C).

Typical Values

Grade	Weight (g/m²)	Thickness (mm)	Penetration 0.3 µm (%)*	Particle retention in liquids (µm)	Filtration speed (mL/min)	Fulfills the requirements in EN 872:2005 (weigh loss)
MGA	54	0.23	< 0.001	1.6	510	Yes
MGB	143	0.70	< 0.001	1.0	210	
MGC	54	0.24	< 0.001	1.2	335	Yes
MGD	120	0.47	< 0.1	2.7	920	
MGF	75	0.38	< 0.001	0.7	110	
MGG	64	0.28	≤ 0.001	1.5	600	
13440	88	0.44		0.7	120	Yes
MG 160	75	0.35	< 0.002	1.2	400	
MG 550-HA	65	0.27		1.5	400	

^{*} Measurement according to EN 143 (0.3 µm, 5.3 cm/s, paraffin oil)

Ordering Information

Filter Discs

Dia. in mm	MGA (100 Pieces)	MG 160 (50 Pieces)	MGB (50 Pieces)	MGC (100 Pieces)	MGD (50 Pieces)
21			FT-3-1102-021		
25	FT-3-1101-025		FT-3-1102-025	FT-3-1103-025	FT-3-1104-025
37	FT-3-1101-037	FT-3-01110-037			
47	FT-3-1101-047	FT-3-01110-047	FT-3-1102-047	FT-3-1103-047	FT-3-1104-047
50	FT-3-1101-050	FT-3-01110-050	FT-3-1102-050	FT-3-1103-050	FT-3-1104-050
55	FT-3-1101-055		FT-3-1102-055	FT-3-1103-055	
70	FT-3-1101-070	FT-3-01110-070	FT-3-1102-070	FT-3-1103-070	FT-3-1104-070
80	FT-3-1101-080				
90	FT-3-1101-090	FT-3-01110-090	FT-3-1102-090	FT-3-1103-090	FT-3-1104-090
100	FT-3-1101-100	FT-3-01110-100	FT-3-1102-100	FT-3-1103-100	FT-3-1104-100
110	FT-3-1101-110	FT-3-01110-110	FT-3-1102-110	FT-3-1103-110	FT-3-1104-110
125	FT-3-1101-125		FT-3-1102-125	FT-3-1103-125	FT-3-1104-125
150	FT-3-1101-150		FT-3-1102-150	FT-3-1103-150	FT-3-1104-150
293					FT-3-1104-293

Dia. in mm	MGF (100 Pieces)	MGG (100 Pieces)	MG 550-HA (100 Pieces)	13440*
24			FT-3-01147-024	
25	FT-3-1105-025	FT-3-1106-025		
42				1344042Q
44				1344044Q
47	FT-3-1105-047	FT-3-1106-047	FT-3-01147-047	1344047Q
50	FT-3-1105-050	FT-3-1106-050	FT-3-01147-050	1344050Q
55	FT-3-1105-055	FT-3-1106-055	FT-3-01147-055	
70	FT-3-1105-070	FT-3-1106-070	FT-3-01147-070	
90	FT-3-1105-090	FT-3-1106-090	FT-3-01147-090	
100				13440-100K
110	FT-3-1105-110	FT-3-1106-110	FT-3-01147-110	
125	FT-3-1105-125	FT-3-1106-125	FT-3-01147-125	
130				13440-130K
150	FT-3-1105-150	FT-3-1106-150		13440-150K
293	FT-3-1105-293			13440-293K

* Q = 500 pieces | K = 50 pieces Other dimensions as well as sheets are available on request



Glass Microfiber Filters

With Binder



These filters are mostly used either for monitoring air and gas or as a prefilter. They are manufactured with synthetic binding agents to ensure that the filter has a defined strength.

They are mechanically and chemically stable, have a temperature resistance up to 180°C and – depending on the binding agent used – are either hydrophobic or hydrophilic.

Typical Values

Grade	Weight (g/m²)	Thickness (mm)	Penetration 0.3 μm (%)*	Pressure drop 5.3 cm/s (Pa)	Binding agent
MG 227/1/60	60	0.32	< 0.5	260	Hydrophobic
13430	220	1.25	0.02	360	Hydrophilic
13400	73	0.39	0.015	363	Hydrophilic
MG 227	75	0.40	< 0.01	350	Hydrophobic
MG 161	75	0.40	≤ 0.002	≤ 580	Hydrophilic
MG 400	80	0.38	< 0.001	400	Hydrophilic
MG 1387/1	90	0.37	≤ 0.003	400	Hydrophilic

^{*} Tested and classified according to the Standard EN 143


Ordering Information

Filter Discs

Dia.	MG 221/1/60	13430**	13400**	MG 227	MG 1387/1
in mm	(100 Pieces)			(100 Pieces)	(50 Pieces)
13			1340013S		
20			1340020S		
25			1340025Q		
42			1340042Q		
44			1340044Q		
45			1340045Q		FT-3-01125-045
47		1343047S	1340047Q	FT-3-01120-047	FT-3-01125-047
50			1340050Q		FT-3-01125-050
55				FT-3-01120-055	FT-3-01125-055
80			1340080N		
100		13430-100K	13400-100K		
110				FT-3-01120-110	FT-3-01125-110
120			13400-120K		
124		13430-124K	13400-124K		
125					FT-3-01125-125
127		13430-127K	13400-127K		
130		13430-130K	13400-130K		FT-3-01125-130
142		13430-142K	13400-142K		
150	FT-3-01124-150		13400-150K		
293		13430-293K	13400-293K		

^{**} K = 50 pieces, N = 100 pieces, Q = 500 pieces, S = 200 pieces Other dimensions as well as sheets are available on request

Quartz Microfiber Filters

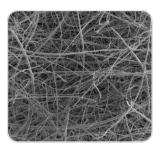

The quartz microfiber material of the Sartorius pre-heated filters, grade Q3400, is made of high-purity quartz microfibers without any addition of glass microfibers or binding agents. In addition, the Q3400 filter grade is tempered to remove all chemically combined water and to give the filters

excellent weight and dimensional stability. Sartorius filters are especially suitable for emissions monitoring at temperatures of up to 900°C and wherever filters of the highest purity are needed.

Typical Values

Grade	Material	Weight (g/m²)		Penetration, 0.3 μm 15 cm/s*	Temperature Resistance
Q3400	100% Quartz microfiber	85	0.43	< 0.002	up to 900°C
	silicium dioxide (Si O_2)				

^{*} Tested and classified according to the Standard EN 143

Ordering Information


Filter Discs

Dia. in mm	Q3400
20	Q340020G
25	Q340025G
30	Q340030G
37	Q340037G
45	Q340045G
47	Q340047G
50	Q340050G
82	Q340082N
90	Ω340090N
142	Q3400-142K
150	Q3400-150K

 * G = 25 pieces, K = 50 pieces, N = 100 pieces Other dimensions as well as sheets are available on request




Membrane Filtration – Quick Selection Guide

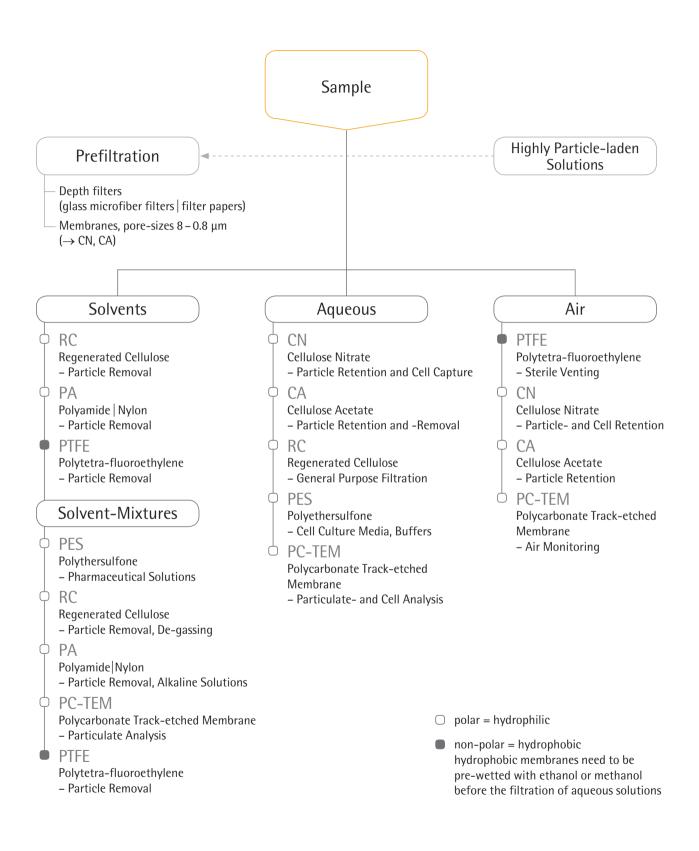

Membrane Filters

Cellulose Nitrate (Mixed Cellulose Ester)

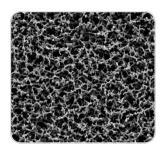

Cellulose nitrate membrane filters are indicated for many general laboratory applications where a membrane with a high non-specific adsorption is suitable. They are hydrophilic, have high flow rates thanks to their symmetric structure and are compatible with aqueous solutions

(pH 4 to 8), hydrocarbons and several other organic solvents. The cellulose nitrate membranes are available in different pore sizes from 0.2 μm to 8 μm .

Typical Values

Туре	Pore Size (μm)	Thickness (μm)	Bubble Point (bar)	Water Flow Rate (mL/min/cm²/bar)	Burst Pressure (bar)
11327	0.2	130	4.2	25	≥ 0.35
11306	0.45	130	2.4	70	≥ 0.3
11305	0.65	130	2	130	≥ 0.25
11304	0.8	130	1.4	200	≥ 0.2
11303	1.2	130	1	200	≥ 0.2
11302	3	130	0.5	430	≥ 0.2
11342	5	130	0.5	570	≥ 0.15
11301	8	130	0.3	750	≥ 0.1


Ordering Information

Filter Discs

- Tiller L	JISCS			
Dia. in mm	11301 (8 μm)*	11302 (3 μm)*	11303 (1.2 μm)*	11304 (0.8 μm)*
13	1130113N	1130213N	1130313N	1130413N
20				1130420N
25	1130125N	1130225N	1130325N	1130425N
37	1130137N			1130437N
47	1130147N	1130247N	1130347N	1130447N
50	1130150N	1130250N	1130350N	1130450N
70	1130170G			
90		1130290G	1130390G	1130490G
100	11301-100N	11302-100G	11303-100G	11304-100G
142	11301-142G	11302-142G	11303-142G	11304-142G
293	11301-293G	11302-293G	11303-293G	11304-293G
Dia. in mm	11305 (0.65 µm)*	11306 (0.45 μm)*	11327 (0.2 μm)*	11342 (5 μm)*
13	1130513N	1130613N	1132713N	1134213N
20		1130620N		
25	1130525N	1130625N	1132725N	1134225N
37		1130637N		
47	1130547N	1130647N	1132747N	1134247N
50	1130550N	1130650N		1134250N
85		1130685N		
90		1130690N		1134290G
100	11305-100N	11306-100N		11342-100G
110		11306-110N		
142		44000 440 N	44007 440 N	11040 140 N
	11305-142N	11306-142N	11327-142N	11342-142N
293	11305-142N 11305-293G	11306-142N 11306-293G	11327-142N	11342-142G

^{*} G = 25 pieces, N = 100 pieces

Cellulose Acetate

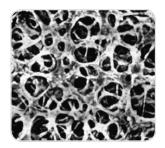

Cellulose acetate membranes combine high flow rates and thermal stability with very low adsorption characteristics, and are therefore excellently suited for use in pressure filtration devices. They are hydrophilic, have high flow rates thanks to their symmetric structure and are compatible with aqueous solutions

(pH 4-8), oils, alcohols and other organic solvents. The 0.2 μm membrane is the filter of choice for sterile filtration of aqueous solutions, such as nutrient media, buffers and sera. The cellulose acetate membranes are available in different pore sizes from 0.2 to 5 μm .

Typical Values

Туре	Pore Size (μm)	Thickness (μm)	Bubble Point (bar)	Water Flow Rate (mL/min/cm²/bar)	Burst Pressure (bar)
11107	0.2	120	2.9	24	0.8
11106	0.45	120	1.9	69	0.7
11105	0.65	120	1.5	115	0.7
11104	0.8	120	1	200	0.5
12303	1.2	140	0.8	320	0.4
12342	5	140	0.4	570	0.25

Ordering Information


Filter Discs

Dia. in mm	11104 (0.8 μm)*	11105 (0.65 μm)*	11106 (0.45 μm)*	11107 (0.2 μm)*	12303 (1.2 μm)*	12342 (5 μm)*
13	1110413N		1110613N	1110713N		
25	1110425N	1110525N	1110625N	1110725N	1230325N	1234225N
30			1110630N	1110730N		
37			1110637N			
45						
47	1110447N	1110547N	1110647N	1110747N	1230347N	1234247N
50	1110450N	1110550N	1110650N	1110750N	1230350N	
70						
85			1110685N			
90	1110490N	1110590G	1110690G	1110790G		
100			11106-100N	11107-100N	12303-100G	
110			11106-110N			
142	11104-142N	11105-142N	11106-142N	11107-142N	12303-142G	12342-142G
293	11104-293G	11105-293G	11106-293G	11107-293G	12303-293G	12342-293G

^{*} G = 25 pieces, N = 100 pieces

Other dimensions and packaging units are available on request

Regenerated Cellulose

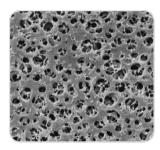

The very low adsorption membranes are hydrophilic, solvent-resistant (pH 3-12) and therefore suited for the particle removal from solvents. The membrane is reinforced with nonwoven cellulose. They are available in two pore sizes: $0.45 \mu m$ and $0.2 \mu m$.

Typical Values

Туре	Pore Size (μm)	Thickness (µm)	Bubble Point (bar)	Water Flow Rate (mL/min/cm²/bar)
18407	0.2	170	4.4	15
18406	0.45	170	2.9	30

Ordering Information

Filter Discs

Dia. in mm	18406 (0.45 μm)*	18407 (0.2 μm)*
13	1840613N	1840713N
25	1840625N	1840725N
47	1840647N	1840747N
50	1840650N	1840750N
90	1840690G	
100	18406-100G	18407-100G
142	18406-142G	18407-142G
293	18406-293G	18407-293G

^{*} G = 25 pieces, N = 100 pieces Other dimensions and packaging units are available on request

Polyethersulfone



Polyethersulfone (PES) membrane filters are hydrophilic, have high flow rates, a low non-specific protein adsorption and are chemically resistant over a pH range of 1 – 14. They are therefore recommended

for the filtration of aqueous solutions as well for protein filtration. Furthermore, the low level of extractables makes them suitable for environmental analysis.

Typical Values

Туре	Pore Size (μm)	Thickness (µm)	Bubble Point (bar)	Water Flow Rate (mL/min/cm²/bar)	
15458	0.1	150	3.8	10	≥ 0.6
15407MI	0.2	150	3.5	25	≥ 0.5
15406	0.45	150	2.6	46	≥ 0.5


Ordering Information

Filter Discs

Dia. in mm	15406 (0.45 μm)*	15407MI (0.2 μm)*	15458 (0.1 μm)*
25	1540625N	1540725MIN	1545825N
47	1540647N	1540747MIN	1545847N
50	1540650N	1540750MIN	1545850N
90		1540790MIK	
142	15406-142G	15407-142MIG	15458-142G
293		15407-293MIG	15458-293G

^{*}G = 25 pieces, K = 50 pieces, N = 100 pieces Other dimensions are available on request

Polyamide

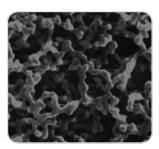

Polyamide membrane filters are hydrophilic and chemically resistant to alkaline solutions and organic solvents. They are therefore recommended for the particle removal from aqueous solutions and solvents for analytical determination such as HPLC, for the sterile filtration of these liquids as well as for applications where a membrane with a relatively high nonspecific adsorption is suitable.

Typical Values

Туре	Pore Size (μm)	Thickness (µm)	Bubble Point (bar)	Water Flow Rate (mL/min/cm²/bar)	
25007	0.2	115	3.2	15	≥ 0.25
25006	0.45	115	2.3	35	≥ 0.23

Ordering Information


Filter Discs

Dia. in mm	25006 (0.45 μm)*	25007 (0.2 μm)*	
13	2500613N	2500713N	
25	2500625N	2500725N	
47	2500647N	2500747N	
50	2500650N	2500750N	
90	2500690G	2500790G	
142	25006-142N	25007-142N	
293	25006-293N	25007-293N	

^{*} G = 25 pieces, N = 100 pieces

Other dimensions and packaging units are available on request

Hydrophobic PTFE

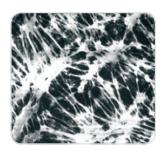

The main application of these membrane filters is the filtration of air, gases or chemicals. They are made of PTFE (polytetra-fluorethylene) only and are therefore permanently hydrophobic. Unlike other (hydrophilic) filter types, they are not wetted by air humidity, allowing unhindered passage of air at low differential pressures as well. PTFE

membrane filters have an excellent chemical compatibility (pH 1 to 14), so that they are also used for the filtration of solvents and acids, to which other filter types are not resistant. Due to their hydrophobic characteristics, they must be pre-wetted with ethanol or methanol before the filtration of aqueous media.

Typical Values

Туре	Pore Size (μm)	Thickness (µm)	Bubble Point (bar)	lsopropanol Flow Rate (mL/min/cm²/bar)	
11807	0.2	65	1.4	11	
11806	0.45	80	0.9	20	
11803	1.2	100	0.45	80	
11842	5	100	0.10	250	

Ordering Information

Filter Discs

Dia. in mm	11803 (1.2 μm)*	11806 (0.45 μm)*	11807 (0.2 μm)*	11842 (5 μm)*
13	1180313N	1180613N	1180713N	
25	1180325N	1180625N	1180725N	1184225N
37	1180337N	1180637N		
42				1184242N
47	1180347N	1180647N	1180747N	1184247N
50	1180350N	1180650N	1180750N	1184250N
90	1180390G	1180690G	1180790G	
100	11803-100G	11806-100G	11807-100G	11842-100G
142	11803-142G	11806-142G	11807-142G	11842-142G
293	11803-293G	11806-293G	11807-293G	11842-293G

^{*} G= 25 pieces, N= 100 pieces

Other dimensions and packaging units are available on request

Polycarbonate Track-Etched



Those white and hydrophilic polycarbonate track-etched membranes are manufactured from high grade polycarbonate film using track-etch technology. Their capillary pore structure is uniform and precise, with a narrow pore size distribution to retain particles on their surface. Track-etched membranes are an excellent choice for accurate fractionation of particulates because of their precise pore size.

Track-etch technology offers the user distinct performance advantages when excellent surface capture and high sample visibility are required. Their main applications are particulate analysis, epifluorescence microscopy, fluid clarification, cytology, cell biology, bioassays, water microbiology and environmental analysis.

Typical Values

Туре	Pore Size (μm)	Thickness (µm)	Bubble Point (bar)	Water Flow Rate (mL/min/cm²/0.7 bar)	Burst Pressure (bar)
23058	0.1	25	7.0	≥ 0.5	≥ 0.7
23007	0.2	25	3.5	≥ 10	≥ 0.7
23006	0.4	25	2.0	≥ 30	≥ 0.7
23004	0.8	25	0.6	≥ 40	≥ 0.7

Ordering Information

Filter Discs, 100 Pieces

Dia. in mm	23004 (0.8 μm)	23006 (0.4 μm)	23007 (0.2 μm)	23058 (0.1 μm)
25	2300425N	2300625N	2300725N	2305825N
47		2300647N	2300747N	2305847N
50			2300750N	

Other dimensions are available on request

Blotting | Chromatography Papers



These papers are made of cotton linters only with α -cellulose content of > 98%. These highly pure papers are not only ideal for blotting and chromatography, but also for a

wide range of absorption applications like those common in the life sciences and diagnostics. Below you will find an overview of the most commonly used grades.

Typical Values

Grade	Weight (g/m²)	Thickness (mm)	Capillary Rise (mm/30 min)	Capillary Rise (mm/10 min)	Properties
FN 4	125	0.24	95		Chromatography paper, ash content < 0.04%
FN 7	150	0.32	145		Chromatography paper, ash content < 0.04%
FN 30	320	0.90	240		Chromatography paper, ash content < 0.04%, paper for antibiotic test strips
FN 100	195	0.35	115	70	The most commonly used chromatography and blotting paper
BF 3	330	0.76	30	130	Blotting paper to increase and maintain the transport of liquids

Ordering Information



Sheets in 580×600 mm

Grade FN 4	Grade FN 7	Grade FN 30	Grade FN 100	Grade BF 3
(100 Sheets)	(50 Sheets)	(25 Sheets)	(50 Sheets)	(50 Sheets)
FT-2-504-580600N	FT-2-507-580600K	FT-2-526-580600G	FT-2-527-580600K	FT-2-520-580600K

Other dimensions and shapes are available on request

Nitrocellulose Membrane for Blotting



Sartorius nitrocellulose membranes are available in two pore sizes, 0.22 μ m and 0.45 μ m. Both versions combine the advantages of high protein binding capacity with low background and high membrane stability, which ensures easy handling. Due to its large surface area, the 0.22 μ m

membrane version is recommended for small proteins. Sartorius blotting membranes are ideal for western blotting, DNA blotting as well as dot or slot blots. They have been optimized for all protein blotting systems, such as electrotransfer, semi-dry or simple capillary blotting.

Typical Values

	0.22 μm	0.45 μm
Material	Cellulose nitrate	Cellulose nitrate
Thickness	130 μm	130 μm
Water flow rate	27 mL/(min.cm ² bar)	70 mL/(min.cm ² bar)
Bubble point	4.4 bar	2.4 bar
Wettability in water	≤ 1 S	≤ 1 S
Extractable content in water	≤ 1%	≤ 1%
Burst pressure	0.8 bar	0.2 bar
Binding capacity for IgG	200 μg/cm²	200 μg/cm²

Ordering Information

	Roll Size	Order No.
NC 0.22 μm	30 cm×3 m	1132741BL
NC 0.45 μm	30 cm×3 m	1130641BL

All indicated data to be understood as typical average values

Re-usable, 13 mm Syringe Filter Holders

For the Ultracleaning of Small Volumes Up to About 10 mL



PTFE Holder for Solvents and Chemicals Made completely of PTFE, this holder is unaffected by chemicals and contains no trace elements which could be released into the liquid being filtered. It is therefore extremely well suited for particle removal from samples and reagents for analytical methods, such as NMR samples. Other

benefits of this application are the low hold-up volume, the easy cleaning and the drying at a temperature of 180°C. The construction of the holder ensures leak proof sealing without a sealing ring, and avoids twisting of the membrane filter when the top is tightened onto the base.



Specifications

Connectors	Female Luer Lock inlet, luer slip outlet
Chemical compatibility	As for PTFE
Filtration area	0.5 cm ²
Materials	PTFE top and bottom parts
Max. operating pressure	5 bar 500 kPa 72.5 psi
Membrane filter diameter	13 mm
Sterilization	By autoclaving (max. 134°C) or by dry heat (max. 180°C)
Hold-up volume	Less than 0.03 mL after overcoming the bubble point (0.3 mL before)

Ordering Information

Description	Order No.	
13 mm PTFE Syringe Filter Holder	16574	



Polycarbonate Holder for Aqueous Solutions

This inexpensive filter holder is made of clear, autoclavable polycarbonate.
The silicone gasket enables a leak-free

filtration at pressures of up to 7 bar by simply screwing it together manually. Filter supports in the top and bottom parts allow filtration in either direction.

Specifications

Connectors	Female Luer Lock inlet, luer slip outlet
Chemical compatibility	As for polycarbonate and silicone
Filtration area	0.5 cm ²
Materials	Polycarbonate top and bottom part, silicone gasket
Max. operating pressure	7 bar 700 kPa 101.5 psi
Membrane filter diameter	13 mm
Sterilization	By autoclaving at 121°C
Hold-up volume	Less than 0.2 mL after overcoming the bubble point (0.3 mL before)

Ordering Information

Description	Order No.
13 mm Polycarbonate Syringe Filter Holder, pack of 12	16514E
Silicon gasket, 10 × 14.9 × 0.5 mm, pack of 10	6980569

Re-usable 25 mm Syringe Filter Holders

For the Ultracleaning and Sterilizing Filtration of Volumes of Up to About 100 mL



Stainless Steel Holder for Solvents and Chemicals

Made of stainless steel, this holder is heatresistant, and the chemical compatibility depends only on the inserted filter type. The top part can easily be mounted on the bottom part using the enclosed tightening tool. Filter supports in the top and bottom parts allow filtration in either direction.

Specifications

Connectors	Female Luer Lock inlet, luer slip outlet
Chemical compatibility	As for stainless steel
Filtration area	3 cm ²
Materials	Stainless steel (1.4305) top and bottom parts
Max. operating pressure	7 bar 700 kPa 101.5 psi
Membrane filter diameter	25 mm
Sterilization	By autoclaving (max. 134°C) or by dry heat (max. 180°C)
Hold-up volume	Less than 0.1 mL after overcoming the bubble point (0.3 mL before)

Ordering Information

Description	Order No.
25 mm Stainless Steel Holder	16214
Tightening tool, Polyman 24/5	6980595



Polycarbonate Holder for Aqueous Solutions

This inexpensive filter holder is made of clear, autoclavable polycarbonate. The silicone gasket enables a leak-free filtration

at pressures of up to 7 bar by simply screwing it together manually. Filter supports in the top and bottom parts allow filtration in either direction.

Specifications

Connectors	Female Luer Lock inlet, luer slip outlet
Chemical compatibility	As for polycarbonate and silicone
Filtration area	3 cm ²
Materials	Polycarbonate top and bottom parts, silicone gasket
Max. operating pressure	7 bar 700 kPa 101.5 psi
Membrane filter diameter	25 mm
Sterilization	By autoclaving at 121°C
Hold-up volume	Less than 0.3 mL after overcoming the bubble point (0.6 mL before)

Ordering Information

Description	Order No.
25 mm Polycarbonate Syringe Filter Holder, pack of 12	16517E
Silicone gasket, 20.5 × 26.5 × 0.5 mm, pack of 10	6980570

25 mm Glass Vacuum Filter Holder

For Hybridization Tests, Particle Testing and Clarification






This filter holder is available in two versions differing from each other only in the type of the filter support. The filter with glass frit ensures uniform distribution of retained particles and is therefore recommended when the residue on the filter surface is of interest. Because it is easy to clean, the device with the PTFE-coated screen support is preferable when the filtrate is required, or when liquids difficult to remove from glass frits must be examined. The PTFE ring, which holds the

glass frit and the screen support, allows for the autoclaving of the devices with a filter in position and protects the edge of the glass frit from breakage and potential leakage. It has a rim around the upper edge to simplify the positioning of the membrane filter when inserted and a silicone O-ring in the underside for a leak-proof seal on the filtrate side. The funnel-shaped top part simplifies filling in the sample.

Specifications

Outlet spout	12 mm dia.
Parts and materials	Borosilicate glass funnel and base PTFE glass filter support (type 16306) or PTFE stainless steel filter support, coated with PTFE (type 16315) Silicone O-ring 25 × 3 mm Anodized Aluminium clamp
Chemical compatibility	As for glass, PTFE and silicone. The silicone O-ring can be replaced by a fluoroelastomer O-ring (order no. 00118)
Funnel capacity	30 mL
Filtration area	3 cm ²
Max. operating pressure	Only for vacuum
Suitable membrane filter diameter	25 mm (or 24 mm)
Sterilization	By autoclaving (max. 134°C) or by dry heat (max. 180°C)

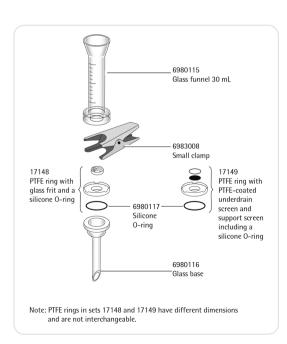

Ordering Information

Description	Order No.
Glass vacuum filtration holder	16306
for 25 mm (or 24 mm) membrane filter,	
with glass frit filter support	
Glass vacuum filtration holder	16315
for 25 mm (or 24 mm) membrane filter,	
with PTFE-coated screen filter support	

50 mm Glass Vacuum Filter Holder

For Particle Testing or Clarification and Sterile Filtration




This filter holder is available in two versions differing from each other only in the type of the filter support. The device with glass frit ensures uniform distribution of retained particles and is therefore recommended, when the residue on the filter surface is of interest. Because it is easy to clean, the device with the PTFE-coated screen support is preferable when the filtrate is required, or when liquids difficult to remove from

glass frits must be examined. The PTFE ring, which holds the glass frit and the screen support, allows the autoclaving of the devices with a filter in position and protects the edge of the glass frit from breakage and potential leakage. It has a rim around the upper edge to simplify the positioning of the membrane filter when inserted and a silicone 0-ring in the underside for a leak-proof seal on the filtrate side.



Specifications

Outlet spouts	15 mm dia.
Parts and materials	Borosilicate glass funnel and base Silicone caoutchouc lid PTFE glass filter support (type 16307) or PTFE stainless steel filter support, coated with PTFE (type 16316) Silicone O-ring 45 × 3 mm Anodized Aluminium clamp
Chemical compatibility	As for glass, PTFE and silicone. The silicone O-ring can be replaced by a fluoroelastomer O-ring (order no. 00124).
Funnel capacity	250 mL
Filtration area	12.5 cm ²
Max. operating pressure	Only for vacuum
Suitable membrane filter diameter	50 mm (or 47 mm)
Sterilization	By autoclaving (max. 134°C) or by dry heat (max. 180°C)

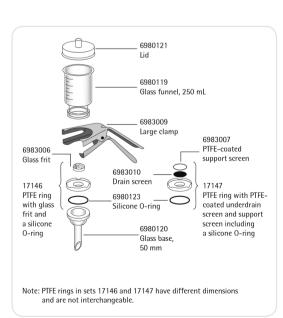

Ordering Information

Description	Order No.
Glass vacuum filtration holder for 50 mm (or 47 mm) membrane filter, with glass frit filter support	16307
Glass vacuum filtration holder for 50 mm (or 47 mm) membrane filter, with PTFE-coated screen filter support	16316

All-Glass Vacuum Filter Holder

For Analytical Determinations, Particle Removal from Solvents




All areas, where liquid and device can come into direct contact, are made of glass or PTFE. The device, in combination with solvent-resistant, hydrophilic RC-membranes, is therefore ideal for ultracleaning and degassing solvents and solvent mixtures for HPLC, GC and AA. Convenience of handling is ensured by several beneficial features. A 6 mm wide

non-ground rim above the ground glass neck of the suction flask prevents the filtrate from contacting grease on the ground glass surface and so avoids its contamination while being poured out of the flask. The hose nipple connector is made of polypropylene for safe connection of the vacuum hose. The filtrate outlet spout ends well below the entrance to this hose nipple.



Specifications

Parts and materials	Borosilicate glass funnel, base and flask, sintered glass frit in a PTFE ring and fluoroelastomer O-ring (45 × 3 mm) underneath, anodized aluminium clamp
Chemical compatibility	As for glass and PTFE
Funnel capacity	250 mL
Capacity of the filtrate flask	1 liter
Filtration area	12.5 cm ²
Max. operating pressure	Only for vacuum
Suitable membrane filter diameter	50 mm (or 47 mm), 40 or 42 mm prefilter
Sterilization (without connector)	By autoclaving (max. 134°C) or by dry heat (max. 180°C)

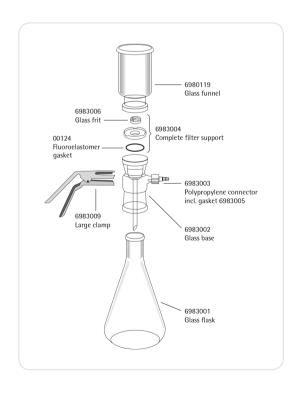

Ordering Information

Description	Order No.	
All-glass vacuum filter holder	16309	
for 50 mm (or 47 mm) membrane filter,		
with vacuum-resistant flask, capacity 1 liter		

Polycarbonate In-Line Filter Holder

For the Filtration of Liter Volumes of Aqueous Solutions




This holder is made of stable, autoclavable polycarbonate. This practical holder is suitable for many simple laboratory filtrations. It can be connected to a peristaltic pump or a pressure container. The bell-shaped base protects the filtrate from repeated contamination while flowing in a receiver. The holder is characterized

by an excellent resistance to pressure and density setting by simple hand-tightening. The transparent top part allows the visual control of the correct fit of the O-ring. The hose nipples can be replaced by luer connectors to use it as a large area syringe filter holder.



Specifications

Chemical compatibility	As for polycarbonate, polypropylene and silicone
Filtration area	12.5 cm ²
Weight	83 g
Threads for connectors	M 12×1 female thread
Materials	Polycarbonate top part, base part and hose nipple, polypropylene filter support, silicone O-ring (40×5 mm)
Max. operating pressure	7 bar 700 kPa 101.5 psi
Suitable membrane filter diameter	50 mm (40 or 42 mm prefilter)
Sterilization	By autoclaving at 121°C The material withstands repeated cycles, provided aggressive cleaning agents are completely washed off and that the boiler water does not contain anti-corrosive or anti-scaling additives.

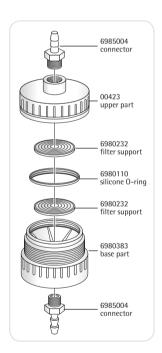

Ordering Information

Description	Order No.
Polycarbonate in-line filter holder for 50 mm membrane filter,	16508B
pack of 5.	

25 mm Stainless Steel Filter Holder

For In-Line Filtration



The G½ connection threads with density barrel, guarantee leak-proof sealing of the hose nipple and the holder without sealing rings. Other connectors, available as accessories, fit the holder onto reducing


valves or pumps with $G\frac{3}{8}$ female thread order no. 01029) or onto pressure tanks with $G\frac{3}{8}$ male thread (order no. 00177).

Specifications

Connectors	Hose nipples DN10
Filtration area	3 cm ²
Weight	ca. 170 g
Materials	Stainless steel, except silicone O-ring (21 × 2 mm) and aluminium closing ring
Max. operating pressure	5 bar 500 kPa 72.5 psi
Suitable membrane filter	25 mm (20 mm prefilter for the filtration of liquids only)
Sterilization	By autoclaving (max. 134°C) or by dry heat (max. 180°C)

Ordering Information

Description	Order No.	
Stainless steel pressure filter holder for 25 mm dia.	16251	
membrane filter.		

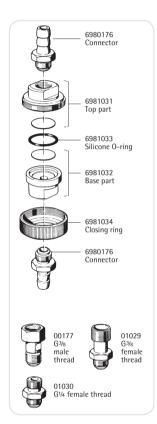

47 mm Stainless Steel Filter Holder

For In-Line Filtration



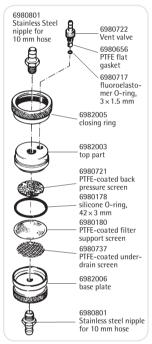


Diagram for 16254

The filter holder is suitable for a pressure of up to 20 bar. The inlet side valve is convenient for the intermittent run-off of waste water. Other connectors, available as accessories, fit the holder onto reducing

valves or pumps with G_{8} female thread (order no. 17089) or onto pressure tanks with G_{8} male thread (order no. 17069) or on taps with G_{4} male thread (order no. 17068).

Specifications

Hose nipples DN10
M12×1
13 cm ²
ca. 490 g
Stainless steel, except silicone O-ring (42×3 mm), PTFE and fluoroelastomer valve seals
20 bar 2,000 kPa 290 psi
47 mm (40 or 42 mm prefilter)
By autoclaving (max. 134°C) or by dry heat (max. 180°C)

Ordering Information

Description	Order No.
Stainless steel filter holder for 47 mm membrane filter (with adapter M12×1 male thread to hose barb DN10, Mat. 316, ref. 6980801) – Replacement parts are shown in the diagram	16254
Stainless steel filter holder for 47 mm membrane filter (with adapter M12×1 male thread to hose barb DN 4 to 5, Mat. 316, ref. 6981132)	16278
Stainless steel back pressure screen	69807211
Stainless steel filter support screen	69801801
Stainless steel underdrain screen	00181
Stainless steel connector M12×1 male thread to hose barb DN 4-5	6981132
Adapter Quick connect nipple length 60 mm male part to male thread M12 x 1, Mat 316	170901

Stainless Steel Pressure Filter Holder

For the Filtration of Up to 5 Liter Volumes




A practical filter holder for many laboratory filtrations. It can be attached to a tripod with the help of a steel rod which can be screwed in. The hose nipple is screwed into the side of the top part, leaving room for a large filling opening. This makes pouring in the sample easier, and the sample can be refilled without removing the tube connection to the pressure source. Leak-proof sealing is achieved by

hand-tightening the closing ring. For the filtration of small volumes (up to about 200 ml of soil samples or viscous liquids, such as oils), the holder is connected directly to a pressure source. For the filtration of up to 5 liter volumes of relatively easily filterable liquids (e.g. buffer solutions, solutions for cell counters and tissue culture solutions), it is used in combination with a pressure tank.



Specifications

Chemical compatibility	As for stainless steel, PTFE and silicone. If required, the silicone O-ring in the filter support can be replaced by a fluoroelastomer O-ring 00179 or a PTFE O-ring 17038 (by reducing the max. operating pressure to 4 bar 58 psi); the silicone O-ring in the top part can be replaced by a fluoroelastomer O-ring 17145.
Filtration area	13 cm ²
Weight	960 g
Threads for connectors	M 12×1 female thread
Materials	Top part, barrel, base part, corrugated iron, closing ring, closure cap, back pressure screen and stainless steel hose nipples 1.4401 (AISI 316), PTFE-coated stainless steel filter support, silicone O-rings, 41×2 mm (top part) and 42×3 mm (filter support), PTFE-sealing (cap).
Max. operating pressure	10 bar 1,000 kPa 145 psi
Suitable membrane filter diameter	47 mm (40 or 42 mm prefilter)
Sterilization	By autoclaving (max 134°C) or by dry heat (180°C)

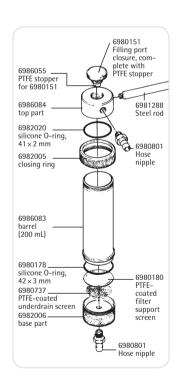

Ordering Information

Description	Order No.
Stainless steel pressure filter holder	16249
Stainless steel pressure filter holder with double jacket	162493

Replacement Parts

Description	Order No.
Fluoroelastomer O-ring, 42 × 3 mm	00179
PTFE O-ring, 42×3 mm	17038
Fluoroelastomer O-ring for upper part, 41 × 2 mm	17145

Other replacement parts are shown in the diagram or on page 2.

Chemical-resistant PTFE Filter Holder

For the Filtration of Aggressive Liquids



The holder hinders the release of trace elements into the filtrate and is resistant to almost all chemicals. The fluoroelastomer O-ring in the top part allows easy hand tightening, and can be replaced by a PTFE

O-ring, order no. 17039. The 6 mm outlet nipple is an integral part of the base, the 10 mm inlet hose nipple can be replaced by a G3/8 connector, order no. 17051.

Specifications


Chemical compatibility	As for PTFE and fluoroelastomer
Filtration area	12.5 cm ²
Thread for inlet connector	M 14×1.5 male thread
Materials	Top part, barrel, base part: corrugated iron, hose nipples and filter support with 40×3.5 mm 0 -ring: PTFE, locking rings: aluminium 39×3.5 mm fluoroelastomer 0 -ring (top part)
Max. operating pressure	5 bar 500 kPa 72.5 psi
Suitable membrane filter diameter	47 mm
Sterilization	By autoclaving (max 134°C) or by dry heat (180°C)

Ordering Information

Description	Order No.
PTFE pressure filter holder, 47 mm, with 200 mL capacity.	16579

Replacement Parts

Description	Order No.
PTFE O-ring, 39 × 3.5 mm	17039

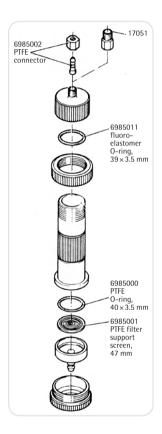

Combisart® Manifolds

1-, 3- and 6-Branch



Made of high-grade stainless steel (B.S. 304S3 | AISI 304); accommodates any type of vacuum funnel. Stainless steel three-way valves (taps) allow the vacuum for each filter holder to be individually

controlled and each holder to be sterilely vented. The low height of the manifold ports is particularly advantageous for working on a clean bench.

16832-CS

16831-CS

Ordering Information

Combisart® Manifolds, without Base Support and Frit	Order No.
Combisart® 1-branch manifold	16844
Combisart® 3-branch manifold	16842
Combisart® 6-branch manifold	16843
Combisart® Sets, Stainless Steel Capacity	Out an Na
Combisart Sets, Stainless Steel Capacity	Order No.
1-branch 1×100 mL	16844-CS
• • •	
1-branch 1×100 mL	16844-CS

In each set stainless steel funnels with lids are preassembled.

Accessories and Replacement Parts

6-branch 6×100 mL

6-branch 6×500 mL

Description	Pack Size	Order No.
Plug, conical, to close the venting hole beside the 3-way valve	10	6980225
Silicone O-ring for manifold female threads	3	6980235
Rubber tubing, 1 m	1	16623

Glass Filter Holders; 30, 250 mL

For Particle Counting



Glass Filter Holders

Two compact vacuum filter holders for easy particulate analysis. Both the top and bottom part of the filter holders are easily and securely fastened together using the

metal clamp. The centering rim on the filter support ensures correct positioning of the membrane filter. The glass frit filter support guarantees uniform distribution of retained particles on the filter surface.



Ordering Information

Description		Order No.
Glass filter holder	30 mL	16306
Filter diameter	25 mm (or 24 mm)	
	Prefilter, 20 mm	
Filtration area	3 cm ²	
Capacity	30 mL	
Outlet	12 mm outer diameter	
Glass filter holder	250 mL	16307
Filter diameter	47 mm (or 50 mm)	
	Prefilter, 40 mm	
Filtration area	12.5 cm ²	
Capacity	250 mL	
Outlet	15 mm outer diameter	



Adapter, 16836 | Adapter, 16837

For use of a glass filter holder, 16306 or 16307, on a Combisart® stainless steel manifold.

Ordering Information

Description	Order No.	
Adapter with 11 mm opening in stopper; for using filter holder 16306 on a Combisart® manifold	16836	
Replacement stopper for 16836	00280	
Adapter with 14 mm opening in stopper; for using filter holder 16307 on a Combisart® manifold	16837	
Replacement stopper for 16837	00281	

Polycarbonate Filter Holders

For Particle Counting



Polycarbonate Filter Holder, 250 mL

This reusable, practical filter holder made of autoclavable plastic is ideal for analytical testing outside the laboratory. For use with 47 mm membrane filters.

Outlet: TR 20×2 mm male thread



Ordering Information

Description	Order No.
Polycarbonate filter holder without receiver flask	16511
Polycarbonate filter holder with receiver flask	16510
Hand vacuum pump with gauge and 60 cm PVC tubing	16673

Ready-to-Use Biosart® 250 Funnels

For Particle Counting



Biosart® 250 Funnel

The Biosart® 250 Funnel has been specially designed for analytical quality assurance. The sterile 250 mL plastic funnel guarantees fast filtration and high sample

throughputs during routine testing. Its large inner diameter allows high flow rates, and the tapered inner walls permit thorough flushing of the funnel, after filtration.



Description	Order No.
Biosart® 250 Funnel, 50 units, sterile-packaged	16407-25-ALK



Single Support, 16840

For adapting a Biosart® 250 Funnel for use on a Combisart® stainless steel manifold.

Ordering Information

Description	Order No.
Stainless steel filter support for stainless steel manifold.	16840

Replacement

Description	Order No.
Stainless steel frit for 50 mm membrane filters	6980102
Stainless steel frit for 47 mm membrane filters	6980103
Silicone flat gasket underneath the frit	6980124
PTFE flat gasket underneath the frit	6980104
Silicone O-ring for 16840 male thread	6980274

Did You know...

... that you can order more lab catalogs?

Get an overview of our high-quality laboratory instruments, high-grade consumables and excellent services.

Our product portfolio focuses on laboratory instruments, such as lab balances, pipettes and laboratory water purification systems. Moreover, we offer the widest range of consumables like laboratory filters and pipette tips.







Products for Microbiological Control Publication No.: SL-1530-e Order No.: 85034-538-81 **Laboratory Weighing Products**Publication No.: WL-0007-e
Order No.: 98649-018-13

Sartorius Liquid Handling ProductsPublication No.: SUL0002-e
Order No.: 85032-542-25



Laboratory Water Purification Products

Turning science into solutions

Laboratory Filtration ProductsPublication No.: SLU0006

Purification ProductsPublication No.: SL-0002-e
Order No.: 85037-548-54

Laboratory Water

Publication No.: SLU0006-e Order No.: 85032-543-83



Download our lab catalogs from our website: www.sartorius.com/lab-catalogs

Europe

Germany Sartorius Lab Instruments GmbH & Co. KG Otto-Brenner-Strasse 20 37079 Goettingen

Phone +49.551.308.0

France & Suisse Romande

Sartorius France 2. rue Antoine Laurent de Lavoisier ZA de la Gaudrée 91410 Dourdan

Phone +33.1.70.62.50.00

Austria

Sartorius Austria GmbH Modecenterstrasse 22 1030 Vienna

Phone +43 1 7965760 0

Belgium

Sartorius Belgium N.V. Rue Colonel Bourg 105 1030 Bruxelles

Phone +32.2.756.06.90

Finland & Baltics

Sartorius Biohit Liquid Handling Oy Lainnatie 1 00880 Helsinki

Phone +358.9.755.951

Hungary

Sartorius Hungária Kft. Kagyló u. 5. 2092 Budakeszi

Phone +3623.457.227

Ireland

Sartorius Ireland Ltd. Unit 41, The Business Centre Stadium Business Park Ballycoolin Road Dublin 11

Phone +353.1.8089050

Italy

Sartorius Italy S.r.l. Via Torino 3/5 20814 Varedo (MB)

Phone +39.0362.5557.11

Netherlands

Sartorius Netherlands B.V.

Phone +31.30.60.53.001

info.netherlands@sartorius.com

Sartorius Poland sp.z o.o. ul. Wrzesinska 70 62-025 Kostrzyn

Phone +48.61.6473830

Russian Federation

LLC "Sartorius RUS" Vasilvevsky Island 5th line 70, Lit. A 199178 St. Petersburg

Phone +7.812.327.53.27

Spain & Portugal

Sartorius Spain, S.A. Avda, de la Industria, 32 Edificio PAYMA 28108 Alcobendas (Madrid)

Phone Spain +34.913.586.095 Phone Portugal +351.800.855.800

Switzerland

Sartorius Mechatronics Switzerland AG Ringstrasse 24a 8317 Tagelswangen (ZH)

Phone +41.44.746.50.00

U.K.

Sartorius UK Ltd. Longmead Business Centre Blenheim Road, Epsom Surrey KT19 9QQ

Phone +44.1372.737159

Ukraine

LLS "Sartorius RUS" Post Box 440 "B" 01001 Kiev, Ukraine

Phone +380.44.411.4918

Americas

IISA

Sartorius Corporation 5 Orville Drive, Suite 200 Bohemia, NY 11716

Phone +1.631.254.4249 Toll-free +1.800.635.2906

Argentina

Sartorius Argentina S.A. Int. A. Ávalos 4251 B1605ECS Munro **Buenos Aires**

Phone +54.11.4721.0505

Brazil

Sartorius do Brasil Ltda Avenida Senador Vergueiro 2962 São Bernardo do Campo CEP 09600-000 - SP- Brasil

Phone +55.11.4362.8900

Canada

Sartorius Canada Inc 1173 North Service Road West, D4 Oakville, ON L6M 2V9

Phone +1.905.569.7977 Toll-Free +1.800.668.4234

Mexico

Sartorius de México, S.A. de C.V. Libramiento Norte de Tepotzotlan s/n, Colonia Barrio Tlacateco, Municipio de Tepotzotlan, Estado de México, C.P. 54605

Phone +52.55.5562.1102 leadsmex@sartorius.com

Peru

Sartorius Peru S.A.C. Avenue Alberto del Campo 411 Floor 12 - The Office 15076 - San Isidro, Lima

Phone +51.1.441 0158

Asia | Pacific

Australia

Sartorius Australia Ptv. Ltd. Unit 5, 7-11 Rodeo Drive Dandenong South Vic 3175

Phone +61.3.8762.1800

China

Sartorius (Shanghai) Trading Co., Ltd. 3rd Floor, North Wing, Tower 1 No. 4560 Jinke Road Zhangjiang Hi-Tech Park **Pudong District** Shanghai 201210, P.R. China

Phone +86.21.6878.2300

Hong Kong

Sartorius Hong Kong Ltd. Unit 1012, Lu Plaza 2 Wing Yip Street Kwun Tong Kowloon, Hong Kong

Phone +852.2774.2678

Sartorius Weighing India Pvt. Ltd. #69/2-69/3, NH 48, Jakkasandra, Nelamangala Tg 562 123 Bangalore, India

Phone +91.80.4350.5250

Japan

Sartorius Japan K.K. 4th Fl., Daiwa Shinagawa North Bldg. 8-11, Kita-Shinagawa 1-chome Shinagawa-ku, Tokyo, 140-0001 Japan

Phone +81.3.3740.5408

Malaysia

Sartorius Malaysia Sdn. Bhd Lot L3-E-3B, Enterprise 4 Technology Park Malaysia Bukit Jalil 57000 Kuala Lumpur, Malaysia

Phone +60.3.8996.0622

Singapore

Sartorius Singapore Pte. Ltd 10 Science Park Rd The Alpha #02-13/14 Singapore Science Park II Singapore 117684

Phone +65.6872.3966

South Korea

Sartorius Korea Ltd. 8th Floor, Solid Space B/D, PanGyoYeok-Ro 220, BunDang-Gu SeongNam-Si, GyeongGi-Do, 463-400

Phone +82.31.622.5700

Thailand

Sartorius (Thailand) Co. Ltd. 129 Rama 9 Road, Huaykwang Bangkok 10310

Phone +66.2643.8361-6

